U.S. Application No. 10/550,797 For IDS filed May 1, 2007

PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 6:

(11) International Publication Number:

WO 96/13610

C12Q 1/68, C07H 21/04, C12P 19/34, A61K 39/395, G01N 33/50, C12N 5/08

(43) International Publication Date:

9 May 1996 (09.05.96)

(21) International Application Number:

PCT/US95/11230

A2

(22) International Filing Date:

24 August 1995 (24.08.95)

(81) Designated States: AU, CA, JP, MX, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).

(30) Priority Data:

08/332,420

31 October 1994 (31.10.94)

Published

US

Without international search report and to be republished upon receipt of that report.

(71) Applicant: GERON CORPORATION [US/US]; 200 Constitution Drive, Menlo Park, CA 94025 (US).

(72) Inventors: LINSKENS, Maarten, H., K.; 1891 Channing Avenue, Palo Alto, CA 94303 (US). HIRSCH, Kenneth, S.; 836 Southampton Drive, Palo Alto, CA 94303 (US). VILLEPONTEAU, Bryant; 1371 Greenbrier Road, San Carlos, CA 94070 (US). FENG, Junli; 1371 Greenbrier Road, San Carlos, CA 94070 (US). FUNK, Walter, 4858 Mendota Street, Union City, CA 94587 (US). WEST, Michael, D.; 15 Wakefield Court, Belmont, CA 94002 (US).

(74) Agents: WARBURG, Richard, J. et al.; Lyon & Lyon, Suite 4700, 633 W. Fifth Street, Los Angeles, CA 90071-2066

(54) Title: METHODS AND REAGENTS FOR THE IDENTIFICATION AND REGULATION OF SENESCENCE-RELATED GENES

(57) Abstract

Identification of senescence-related genes can be accomplished by comparing mRNA expression between young and senescent cells. Probes complementary to such genes can be used to detect senescent cells and distinguish between young and senescent cells as well as in screens to identify compounds that alter expression levels of senescence-related genes.

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
AU	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY.	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgystan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic	SD	Sudan
CG	Congo		of Korea	SE	Sweden
CH	Switzerland	KR	Republic of Korea	SI	Slovenia
CI	Côte d'Ivoire	KZ	Kazakhstan	SK	Slovakia
CM	Cameroon	LI	Liechtenstein	SN	Senegal
CN	China	LK	Sri Lanka	TD	Chad
CS	Czechoslovakia	LU	Luxembourg	TG	Togo
CZ	Czech Republic	LV	Latvia	TJ	Tajikistan
DE	Germany	MC	Monaco	TT	Trinidad and Tobago
DK	Denmark	MD	Republic of Moldova	UA	Ukraine
ES	Spain	MG	Madagascar	US	United States of America
FI	Finland	ML	Mali	UZ.	Uzbekistan
FR	France	MN	Mongolia	VN	Viet Nam
GA	Gabon		-		

Description

Methods and Reagents for the Identification and Regulation of Senescence-Related Genes

Background of the Invention

Field of the Invention

The present invention relates to the fields of molecular biology, gerontology, and medical pharmacology and diagnostics.

Description of the Related Art

There is substantial evidence that somatic cells have a finite replicative capacity (Hayflick, 1965, Exp. Cell Res. 37: 614-636, and Hayflick, 1970, Exp. Geront. 5: 291-303) and that this process is a major etiological factor in aging and age-related disease (Goldstein, 1990, Science 249: 1129-1133; Stanulis-Praeger, 1987, Mech. Ageing Dev. 38: 1-48; and Walton, 1982, Mech. Ageing Dev. 19: 217-244). As cells undergo replicative senescence in vitro and in vivo, they not only lose the ability to divide in response to growth stimuli but also exhibit significant deleterious changes in the pattern of gene expression (West, 1994, Arch. Derm. 130: 87-95).

During replicative senescence, cells exhibit an elongation of the G_i phase of the cell cycle, leading to a longer cell time of cycle transit. As the progression from mitotically active to senescent continues, cells fail to respond to mitotic signals and remain instead in G_i . The inability of senescent cells to enter the cell cycle represents a major difference between young and old cells in that young cells become quiescent entering G_0 until such time when they are induced to reenter the cell cycle and divide. Senescent cells exhibit changes in morphology, increasing in size and volume. However, senescent cells remain viable and are metabolically

active. Another characteristic of replicative senescence is that changes in the pattern of gene expression become more dramatic as the cell reaches the end of replicative life. These changes are reflected in a decrease in the expression of "young-specific" genes with an increase in the expression of "old-specific" "senescent-specific" genes. For purposes of the present invention, any gene whose product is differentially expressed between young quiescent cells and senescent cells is a "senescent-related" gene. Not only to do these changes affect the structure and function of the senescent cell, but such changes can also influence the physiology of surrounding cells and tissue matrix by altering the extracellular environment or in a paracrine fashion through the release of different proteins or through changes in cell-cell interactions.

Several senescent specific genes have been described in the scientific literature. Dermal aging is illustrative example. Dermal aging is characterized by changes in the structure and function of extracellular matrix (ECM) proteins. Many of these same changes have been observed in experiments conducted with fibroblasts either grown to senescence or in cells derived from older individuals. Cells derived from older individuals exhibited up to a 4.4-fold greater level of fibronectin mRNA when compared to levels expressed in fetal cells. Similarly, the synthesis of fibronectin is increased in cells grown to senescence. In late passage cells, the fibronectin that is synthesized is structurally different from that observed in younger cells. These changes may reflect age-related changes in the processing Functionally, these changes translate fibronectin mRNA. in a decreased capacity to mediate cell adhesion, cell spreading, and contact formation. When comparing younger older cells, the fibronectin lattice different with the lattice in the older cells tending to

be less well organized. See Eleftheriou et al., 1991, Cellular ageing related proteins secreted by human fibroblasts, Mutat. Res. 256: 127-38; Kumazaki et al., 1993, Enhanced expression of fibronectin during in vivo cellular aging of human vascular endothelial cells and skin fibroblasts, Exp. Cell Res. 205: 396-402; Hara et al., 1993, DNA-DNA subtractive cDNA cloning oligo(dT)30-Latex and PCR: identification of cellular genes which are overexpressed in senescent human diploid fibroblasts, Analyt. Biochem. 214: 58-64; and Martin et al., 1990, Fibronectin and collagen gene expression during in vitro ageing of pig skin fibroblasts, Exp. Cell Res. 191: 8-13.

Another characteristic of dermal aging is that the expression of interstitial collagenase, also known as fibroblast collagenase, has been reported to increase in senescent cells as well as in cells derived from older donors. Not only is there an increase in the collagenase mRNA, but the activity of the enzyme is also increased. These effects appear to be at the transcriptional level and may in part be mediated by interleukin-1 (IL-1), which itself appears to be upregulated during senescence. See Sottile et al., 1989, Regulation of collagenase and collagenase mRNA production in early- and late-passage human diploid fibroblasts, J. Cell. Physiol. 138: 281-290; West et al., 1989, Replicative senescence of human skin fibroblasts correlates with a loss of regulation and overexpression of collagenase activity, Exp. Cell Res. 184: 138-147; Burke et al., 1994, Altered transcriptional regulation of human interstitial collagenase in cultured skin fibroblasts from older donors, Exp. Gerontology 29: 37-53; and Lafyatis <u>et</u> al., 1990, Interleukin-1 stimulates and all-trans-retinoic acid inhibits

collagenase gene expression through its 5' activator protein-1 binding site, Mol. Endo. 4: 973-980.

In addition, PAI-1 expression appears to be regulated as a function of dermal aging at both the mRNA and protein levels, although the mechanism is as yet unclear. See Shay et al., 1992, Re-expression of senescent markers deinduced reversibly immortalized cells. Exp. <u>Gerontology 27</u>: 477-492. Stromelysin mRNA and protein are over-expressed in senescent cells (see Millis et al., 1992, Metalloproteinases and TIMP-1gene expression during replicative senescence, Exp. Gerontology 27: 425-428; and Millis et al., 1992, Differential expression of metalloproteinase and tissue inhibitor of metalloproteinase genes in aged human fibroblast, Cell Res. 201: 373-379), as is tPA (see West, 1994, The cellular and molecular biology of skin aging, Arch. Dermatol. 130: 87-95). The levels of TIMP-2 protein and mRNA were studied in early and late passage human fibroblasts and found to be upregulated during senescence and Millis, 1994, Expression 72-kDa gelatinase and TIMP-2 in early and late passage human fibroblasts , Exp. Cell Res. 213: 148-155).

Several other genes have been found to be over-expressed in senescent cells. Some of these genes appear to play a role in cell growth and signaling, and alteration in such genes may contribute significantly to an alteration in tissue physiology. IL-1 is upregulated during senescence, which can affect the transcription of several ECM genes, including stromelysin, PAI-2, and collagenase. See Kumar et al., 1993, Expression of interleukin-1-alpha and & in early passage fibroblasts from aging individuals, Exp. Gerontology 28: 505-513; and Kumar et al., 1992., Expression of interleukin 1-inducible genes and production of interleukin 1 by aging

human fibroblasts, Proc. Natl. Acad. Sci. USA 89: 4683-7. IFN gamma can act to decrease the expression of several genes that are down-regulated during senescence. See Eleftheriou et al., 1991, supra; and Eleftheriou et al., 1993, A group of three fibroblast secreted polypeptides suppressed by cellular ageing and interferon-gamma, Biochim. Biophys. Acta 1180: 304-12. For example, IFN gamma functions to decrease the expression of collagen and fibronectin.

Fibroblasts are responsible for elastogenesis (see 1989, Braverman, Elastic fiber and microvascular abnormalities in aging skin, Clin. Geriat. Med. 5: 69-90), and an examination of the elastin produced by cultured dermal fibroblasts from individuals increasing age reveals that, in the sixth decade, there is a marked reduction in the synthesis and repair of elastin fibers (see Fazio et al., 1988, Isolation and characterization of human elastin CDNAs, and associated variation in elastin gene expression cultured skin fibroblasts, Lab. Invest. 58: 270-7; and Dalziel, 1991, Aspects of cutaneous ageing, Clin. Exp. Dermatol. 16: 315-23). Cultured dermal fibroblasts also exhibit an age-related decrease in collagen synthesis and increase in degradation either when senescence or when derived from people of various ages. See Mays et al., 1990, Similar age-related alterations in collagen metabolism inrat tissues in vivo and fibroblasts in vitro, Biochem. Soc. Trans. 18: 957: Furth, 1991, The steady-state levels of type I collagen mRNA are reduced in senescent fibroblasts, J. Gerontol. 46: B122-4; and Takeda et al., 1992, Similar, but not identical, modulation of expression of extracellular

matrix components during in vitro and in vivo aging of human skin fibroblasts, J. Cell. Physiol. 153: 450-9. Specifically, collagen type 1, pro alpha 1 and 3 chains, as well as type 3 pro alpha 1, are all down regulated during senescence at both the mRNA (see Hara et al., 1993, supra) and the protein level. See Dumas et al., 1994, In vitro biosynthesis of Type 1 and III collagens by human dermal fibroblasts from donors of increasing age, <u>Mech</u>. <u>Age</u>. <u>Develop</u>. <u>73</u>: 179-187. Similar results were observed in experiments using cultured fibroblasts from pigs (see Martin et al., 1990, supra). In these studies, there was an increase in type III collagen, while type I collagen was synthesized but Thus, the ratio of type I:type III collagen degraded. These results demonstrate that age-related was altered. changes may be species-specific.

At least two proteins involved in the maintenance of the ECM are down-regulated in dermal aging: tissue inhibitor of metallo proteinase 1 (TIMP-1) is downregulated (see West et al., 1989; and Millis et al., 1992, supra) as is osteonectin, a structural ECM protein induced during proliferation. See Reed et al., 1994, TGF-beta 1 induces the expression of type I collagen and SPARC, and enhances contraction of collagen gels, by fibroblasts from young and aged donors, J. Cell. Physiol. 158: 169-79. Early Passage Clone 1 (EPC1, see Pignolo et al., 1993, Senescent WI-38 cells fail to express EPC-1, a gene induced in young cells upon entry into the G_0 state, J. <u>Biol</u>. <u>Chem</u>. <u>268</u>: 8949-8957), which is identical to Pigment Epithelium Derived Factor (PEDF, see Steele et 1993, Pigment epithelium-derived Neurotrophic activity and identification as a member of the serine protease inhibitor gene family, Proc. Natl. Acad. Sci. USA 90: 1526-1530), is an example of a growth

down-regulated in senescent factor that is cells. Another gene that is down-regulated during senescence is a ribosomal protein, L7 (see Seshadri et al., Identification of a transcript that is down-regulated in senescent human fibroblasts. Cloning, sequence analysis, and regulation of the human L7 ribosomal protein gene, J. Biol. Chem. 268: 18474-80), while c-fos induction is also repressed during senescence. See Seshadri and Campisi, 1990, Repression of c-fos transcription and an altered genetic program in senescent human fibroblasts, Science 247: 205-209.

Thus, as an individual grows older, senescent cells make up an increasing percentage of the cells present in the tissues of the aging individual. The altered pattern of gene expression senescent cells exhibit contributes significantly to age-related pathologies. As the number of aged individuals is expected to increase dramatically in the near future, the cost of health-care for the aged will likewise increase dramatically. Reversal, partial reversal, or modulation of senescent gene expression can provide effective therapies for diseases, disease conditions, and pathologies in which replicative cell senescence is an etiological factor.

Consequently, there is profound а need therapeutic agents and treatment regimes based upon the underlying biology of aging and age-related diseases, particularly the biology relating to the fundamental changes in gene expression that contribute to cell senescence and the development of age-related disease. The present invention helps meet that need by providing new methods for culturing senescent cells for use in cell-based assays and screens; for discriminating between genes expressed by young proliferative cells, quiescent cells, and non-proliferating senescent cells; for separating young cells from senescent cells and using

those cells to conduct high-throughput screens based on cell senescence that can be used to identify compounds that, by reversing the senescent phenotype, treat or diminish age-related disease or pathologies; and for treating age-related human disease, as well as providing compounds and reagents useful in those methods.

Summary of the Invention

In a first aspect, the present invention provides a method for identifying and isolating senescence-related genes and gene products, which method comprises: isolating mRNA from senescent cells and young quiescent cells; (b) amplifying said mRNA in a polymerase chain reaction to produce amplified gene sequences; separating said amplified gene sequences by gel electrophoresis; and (d) analyzing said amplified gene sequences separated in step (c) to detect an amplified gene sequence from young quiescent and young dividing cells that is present at a level different from that observed in amplified gene sequences from senescent cells. With this method, one can readily identify and isolate senescence-related genes. For instance, senescence-related genes or gene tags can be physically removed from the gel and sequenced, either directly or after cloning into a suitable recombinant DNA vector.

Thus, in a second aspect, the present invention provides useful nucleic acids in isolated form, which nucleic acids include portions of senescence-related genes and are useful as nucleic acid probes in diagnostic methods, as nucleic acid primers, and as components of recombinant DNA cloning and/or expression vectors. The present invention also encompasses the gene products of senescence-related genes.

In a third aspect, the present invention provides diagnostic methods for detecting senescent cells in culture and <u>in vivo</u> and for distinguishing senescent

cells from non-senescent cells. These methods comprise the steps of: (a) contacting the mRNA present in a cell or tissue with a labelled nucleic acid probe that comprises a sequence of a senescence-related gene under conditions such that complementary nucleic' hybridize to one another; (b) determining whether specific hybridization has occurred; and (c) correlating the presence of senescent and non-senescent cells with the occurrence of hybridization.

These methods are especially useful in conjunction with therapeutic regimens and strategies. For example, one can use a probe (e.g., a nucleic acid or modified nucleic acid complementary to a senescence-related gene or RNA transcript or an antibody against a senescencerelated gene product) to identify the senescent (or young) cells in a tissue sample. By appropriately labeling the probe for detection (e.g., fluorescent molecule or a molecule that serves as a binding partner for another molecule), one can label the senescent (or young) cells in a tissue and then separate the young from the senescent cells. In addition. detectable reagent (e.g., a contrast molecule for magnetic resonance imaging) can be attached an antibody or other substance which is specific for a senescence-related gene product, and can specifically label senescent cells. With such a preparation of cells enriched for either the young or senescent phenotype, one can then perform many useful procedures, including reintroduction of young cells into a host or treatment of senescent cells for reintroduction to the host.

In a related aspect, the identification of senescence-related genes can provide new therapeutic opportunities. Toxic substances (or "toxins") can be provided which are specific for senescent-specific (or senescence-related) gene products. Such toxins can kill cells expressing these senescent-specific or senescence-

related gene products. Toxins can be provided in the form of toxic substrates, wherein the toxic substrates can be activated (become toxic to cells) by senescencerelated gene products. For example, beta-galactosidase is known to be a senescent-specific gene product, so one can prepare a toxic substrate (or "pre-toxin") that requires activation by beta-galactosidase to be activated and become toxic. Such a toxic substrate could be used to eliminate senescent cells in vivo or in vitro. further aspect, antibodies to senescence-related gene products can be linked with toxins or toxic substrates and delivered to cells by methods known to those in the art to destroy cells expressing senescence-related gene In this aspect, it is not necessary that the products. toxic substrates linked to the antibodies are activated by the senescence-related genes, as any type activation known to those in the art (e.g., radiation, etc.) can be used. In another aspect, antisense oligonucleotides, can be targeted to senescence-related genes. Such antisense oligonucleotides can be comprised of ribonucleic acids. deoxyribonucleic acids, modified nucleic mixtures. Those of skill in the art will recognize that the senescence-related genes and gene products of the invention provide a wide array of such agents that can be used to target or direct therapeutic or reagents to young or senescent cells.

In a fourth aspect, the present invention provides a method for screening compounds to identify compounds that can alter gene expression in senescent cells, which method comprises: (a) contacting senescent cells with a compound; (b) determining mRNA expression patterns in said senescent cells; and (c) correlating an alteration in mRNA expression of a senescence-related gene with a compound that can alter gene expression in senescent cells. The present invention also encompasses the

compounds identified by this method and the use of those compound to alter gene expression in senescent cells.

These and other aspects of the invention are described in more detail below.

Brief Description of the Drawings

The figures will first be briefly described.

shows three examples of enhanced differential display autoradiographs. The primer combinations are indicated above the gel. The lanes of the gel are as follows: BJ cells old, Doubling Level (PDL) 90.3, grown in 0.5% serum; (2) BJ cells young, PDL 40, grown in 0.5% serum; (3) IMR90 old, PDL 54, grown in 0.5% serum: (4) IMR90 young, PDL 21.4, grown in 0.5% serum; (5) IMR90 old, PDL 53, grown in 10% serum; and (6) IMR90 young, PDL 27.4 grown in 10% serum.

Figure 2 shows an analysis of the frequency of mismatches in the 3' most 8 bases of the 5' primers used in EDD: (a) frequency of mismatches detected in 34 genes; and, (b) distribution of the position of the mismatch in those sequences that have a single mismatch. Figure 3 shows differentially displayed bands in EDD analysis of young (Y) and senescent (O) fibroblasts. loading order is identical to that of the gel shown in Figure 1. The arrows on the left indicate the differentially displayed band, with the numbers indicating the lane numbers with elevated expression. The sequence of the nucleic acids corresponding from bands comprising fragments of known genes are indicated on the right. The primer sets that were used for detection of the bands were 02 and C for Col3a1, 16 and C for Laminin A, 18 and C for ALDH-1 and 01 and D for IFN gamma.

Figure 4 shows Northern blots of RNA from young and senescent cells. 4A shows the loading order (for 4B, 4C

and 4D) of the RNA samples; 20 µg of total RNA were loaded in each lane, and loading was verified by ethidium bromide staining of an identical gel. The probes shown in 4B target genes previously reported (prior to September, 1995) to be expressed differentially between young and senescent cells. 4C shows probes for known genes that were identified by enhanced differential display to be expressed differentially; and 4D shows probes for novel genes that were identified by enhanced differential display to be expressed differentially.

Description of the Preferred Embodiments

In a first aspect, the present invention provides a method for identifying and isolating senescence-related genes and gene products, which method comprises: (a) isolating mRNA from senescent cells and young quiescent cells; (b) amplifying said mRNA in a polymerase chain reaction to produce amplified gene sequences; separating said amplified gene sequences qel electrophoresis; and (d) analyzing said amplified gene sequences separated in step (c) to detect an amplified gene sequence from young quiescent and young dividing cells that is present at a level different from that observed in amplified gene sequences from senescent This method differs from prior art methods in cells. that the method allows one to identify and isolate senescence-related genes rapidly and efficiently. "senescence-related gene" refers to a gene that expressed at a different level in a senescent cell than in a non-senescent cell of the same cell type. cases, a senescence-related gene will be expressed in a senescent cell and will either be expressed at a lower level or not be expressed at all in a non-senescent cell, in which case the gene is referred to as an In other cases, a senescence-related specific" gene. gene will be expressed in a non-senescent cell and will

either be expressed at a lower level or not be expressed at all in a senescent cell, in which case the gene is referred to as a "young-specific" gene.

The advantages of this method in part result from a comparison of the mRNA populations of different populations, i.e., young quiescent cells and senescent interest and relevance Of greatest comparison between young quiescent and senescent cells, because in vivo, most cells are generally in a quiescent state, unless there is need for proliferative activity, such as during wound-healing or tissue regeneration. comparison between young quiescent cells and senescent cells therefore reflects the in vivo situation most accurately, and genes that are identified in this way have a high likelihood of being differentially expressed Those of skill in the art recognize that in the tissue. growth conditions can be modified to select for population of cells that are mitotically active (typically, high serum concentrations and frequent passaging to keep the cells in a non-confluent state are used to keep the cells dividing, e.g., for fibroblasts, about 10% serum is adequate for this purpose) as opposed to quiescent (typically, low serum concentrations and contact inhibition, i.e., confluency, are used to keep the cells in a quiescent state, e.g., for fibroblasts, about 0.5% serum adequate for is this purpose). Typically, fetal or embryonic cells are not used unless one wants to detect developmentally-regulated genes.

The comparison of the mRNA population produced by each of these cell populations allows one to identify senescence-related genes. In other embodiments of the method, additional cell populations are used to provide additional information. Thus, while the method typically comprises comparing the mRNA population of young cells cultured in low serum (0.5%) with that of old cells cultured in low serum, one can enhance the method by

including in the comparison young cells cultured in high serum (10%) and old cells cultured in high serum. cells cultured in high serum with frequent passaging should be dividing and mitotically active, allowing one to detect mRNAs of growth-specific genes, such as those that encode proteins that regulate the cell cycle. cells cultured in high serum can produce mRNA species not seen in the mRNA populations of young cells cultured in high or low serum or old cells cultured in low serum. Furthermore, mere comparison of the mRNA population of young dividing cells with that of senescent cells might result in the identification of a gene product involved in cell cycling and cell division as a senescence-related gene product, because senescent cells do not divide. using the mRNA population expressed by a young quiescent (non-dividing) cell for the comparison, one can avoid mistaking cell cycling and cell division gene products as senescence-related gene products.

The method is especially advantageous when used in conjunction with an mRNA preparation methodology known as "Enhanced Differential Display" or "EDD" and described fully in U.S. patent application Serial 08/235,180, filed 29 April 1994, incorporated herein by reference. As the name implies, EDD is an improvement of methodology known as "Differential Display" or "DD". (see Liang, & Pardee, 1992, <u>Science 257:</u> involves the use of PCR amplification of DNA fragments that represent the mRNA of a given cell population. of the two primers used in the PCR is complementary to the poly-A tail of the mRNA, while the other primer (the 5' primer) has a randomly selected sequence intended to be complementary to an internal sequence within an mRNA. The annealing conditions for the 5'-primer, which is ten bases long, are chosen to be degenerate, so that only the last six to eight bases determine the sequence homology. Twelve different poly-T primers anchored by

additional nucleotides are run in separate reactions in combination with defined but randomly selected 5'-primers.

Under these conditions, the assay will generate a display of about 30 to 50 to 100 bands that typically range in size from 100-400 base pairs (bp) per PCR reaction when resolved on a standard sequencing gel. application of a sufficient number of primer combinations for each mRNA sample thus generates a catalog of bands (each band is called a "genetag") that represents the 3'end of mRNA from expressed genes (or internal fragments of mRNA molecules that comprise an internal poly A The displays of different mRNA populations can tract). then be compared and differentially-displayed bands These bands can then be cut out of the gel. sequenced and/or cloned, and the DNA fragment can then be used as a probe to isolate cDNA clones for a particular gene, as illustrated in the examples below.

Although DD has been used by many laboratories in its original form, the method has a tendency to exhibit extensive variability from experiment to experiment (see Liang & Pardee, 1993, Nuc. Acids Res. 21: 3269-3275; and Bauer et al., 1993, Nuc. Acids Res. 21: 4272-4280). reproducibility gives rise to many positives, presumably due to the degeneracy of the PCR protocol (see Sun et al., 1994, Cancer Res. 54, 1139-As a result, DD does not consistently allow the generation of a catalog of one mRNA population that can be compared to a catalog generated separately for another mRNA population. The improved process of the present invention, EDD, addresses the basis for the generally poor reproducibility of the DD technique.

The reliability of the DD technique appears to be in part dependent on primer-length and the temperatures at which primer extension is conducted. EDD essentially differs from DD in the use of longer primers and higher

temperatures for primer extension after a limited number of low temperature primer extension steps. primers" it is meant that the primers are greater than 18 nucleotides in length, and preferably are between 18 and 30 nucleotides in length. In DD, all of the primer extension steps are conducted at low temperatures. the first few (2 to 6) primer extension steps in EDD, the temperature of primer annealing and extending is low, allowing primers to bind readily and not necessarily with complete complementarity to mRNA sequences. temperature" or "low stringency" or "low fidelity" conditions, it is meant that the primer annealing/extending temperature is about 37°C to 50°C. Under these low stringency conditions, only the 3'-most 8) nucleotides of the primer dictate specificity of hybridization. However, the primers used are much longer (20 to 30 or more nucleotides) than the hybridizing region active in the first few temperature cycles, and the later cycles of primer annealing and extending are conducted at temperatures that favor higher specificity and result in the replication of only primer extension products formed in the first few cycles, leading to more reproducible results. By "high temperature" or "high stringency" or "high fidelity" it is meant that the annealing/extending temperature is about 55°C to 75°C.

After a number of experiments with EDD, sufficient data was obtained to permit detailed analysis of primer annealing. One major question was the number of bases required at the 3' end of the arbitrary primer in determining specificity. From the described catalog (see below) and other experiments, 34 genes with known sequences were obtained. The site for the arbitrary primer annealing was determined by matching the gene sequences with primer sequences and band size. The complementarity for the eight bases at the 3' end were

determined. Matches with the gene sequences for bases 9-21 (3'-5') in the primer were highly variable and were not considered to contribute to the initial annealing of the primer. Figure 2 summarizes the mismatches in the 3' most 8 bases for 34 genes. Mismatches in the last eight bases were permitted and in some cases more than one mismatch was found. Mismatches occur more frequently toward the 5' end of the primer, as expected. Overall, the data in Figure 2 suggests that a 7 out of 8 match is the typical behavior of a 5' primer in EDD.

Using EDD, one can readily clone and sequence a 100 to 400 bp DNA fragment derived from the 3'-end of the mRNA from an expressed gene. Typically, one employs a total of 100 to 300 different primer combinations and generates from 1,000 to 10,000 (typically 4000-8000) different genetags, although one can get multiple genetags from the same mRNA. The genetags can be readily separated and visualized using gel electrophoresis or other techniques, and the genetags that are differentially expressed can be readily isolated, i.e., cut from the gel, and sequenced either directly or after Typically, about 50% of the sequences can be cloning. determined directly (from a gel isolated band), while the remaining bands must first be cloned to determine the sequence.

Each differentially expressed genetag (sometimes referred to as "band" due to the appearance of the genetag in a gel) isolated in this manner will typically first be used as a probe for Northern analysis to confirm the senescence-related nature of the genetag. RNA samples are prepared from young and old and mitotic and quiescent cells, and then probed with a labelled genetag probe to verify that the genetag is from an mRNA that is differentially expressed between young and senescent Genetags that are indeed specific for young or old cells can then be used as probes for in situ RNA

analysis in tissue or organ sections of young and old donors, including both diseased and normal tissues or organs, to discriminate between young and old cells. Such analysis can be conveniently carried out using the reagents and in situ hybridization protocol described in the SureSite™ II System Manual, commercially available The genetags can also be used as probes from Novagen. for Northern analysis of RNA from cells of a type other than the cell-type from which the gene-tag was derived. Novel genes whose expression changes in association with aging can be cloned and further characterized using methodology well known in the art (see, e.g., Sambrook et Molecular Cloning, A Laboratory Chapters 8 <u>et seg</u>. (2d ed., CSH Press, Cold For instance, the DNA sequence of novel genes can be analyzed by comparison to existing genes using standard molecular biology techniques and further establish the function of the molecule analyzed to encoded by the mRNA from which the genetag is derived. Antibodies can be raised against novel gene products to facilitate this analysis and to provide an antibody-based method for distinguishing young from old cells, or for other purposes, e.g., drug delivery, as described above. The gene products can also be incorporated in a senescent cell-based drug screen, as described more fully below.

The EDD methodology allows one to detect genes that exhibit a difference in the steady-state level of mRNA produced from those genes. Steady-state mRNA levels can regulated at the transcriptional and posttranscriptional level. Old cells can also differ from young cells by altered steady-state mRNA levels and by altered levels of a protein or the activity of a protein, which can be due to alterations in mRNA translation or protein structure. Thus, regulation of gene expression occur by a variety of mechanisms. the

transcriptional level, the production of mRNAs can either increase or decrease. The level of translation or changes in post-translational modification can lead to an increase or decrease in the abundance of proteins. activity of a protein can be modulated or the turnover rate of the protein can change. Each of these mechanisms can in turn be regulated. The methods of the present invention have resulted in the identification of a number of known and previously unreported genes and products, the expression or abundance of which is controlled at least in part by the aging process.

These gene products were identified in a study initiated to examine cell senescence in fibroblasts. Cultured dermal fibroblasts were selected for analysis in study in part due the extensive research to conducted to date on cell senescence in this cell type (see, e.g., Harley et al., 1990, Nature 345: 458-460). Although cultured dermal fibroblasts were used, those of skill in the art will recognize that any cell type could be used in the methods of the invention. exhaustive list of potentially useful cells are; (a) cells of the central nervous system, including astrocytes, endothelial cells, and fibroblasts, which are involved in such age-related diseases as Alzheimer's disease, Parkinson's disease, Huntington's disease, and stroke: (b) cells of the integument, including fibroblasts, sebaceous gland cells, melanocytes, keratinocytes, Langerhan's cells, and hair follicle cells, which are involved in age-related diseases of the integument, such as dermal atrophy, elastolysis and skin wrinkling, sebaceous gland hyperplasia, senile lentigo, graying of hair and hair loss, chronic skin ulcers, and age-related impairment of wound healing; (c) cells of the articular cartilage, such as chondrocytes, the senescence of which leads to the overexpression of the destructive proteins collagenase and stromelysin, which destroy

articular cartilage in osteoarthritis, and lacunal. synovial, connective tissue fibroblasts, which involved in degenerative joint disease; (d) cells of the bone, such as osteoblasts, osteoclast progenitor cells, bone marrow stromal fibroblasts, and osteoprogenitor cells, which are involved in osteoporosis; (e) cells of the immune system such as B and T lymphocytes, monocytes, neutrophils, eosinophils, basophils, NK cells and their respective progenitors, which are involved in age-related immune system impairment; (f) cells of the vascular system, including endothelial cells, smooth muscle cells, and adventitial fibroblasts, which are involved in agerelated diseases of the vascular system including atherosclerosis, calcification, thrombosis, and aneurysms; (g) cells of the eye, such as retinal pigmented epithelium, lens epithelial cells, iris muscle cells (myoblasts), and vascular endothelial cells, which are involved in loss of vision, i.e., age-related macular degeneration; (h) muscle satellite cells involved muscular dystrophy; and (i) cells of the gut, such as intestinal epithelial cells, which are involved malabsorption syndromes

The methods and reagents of the present invention in part arise out of the recognition that the structural and functional changes in organs and tissues intrinsic to the aging process can be attributed to an alteration in the pattern of gene expression accompanies cell senescence. Research into cellular aging has provided insight into the mechanisms through which the lifespan of cells is regulated. The <u>in</u> <u>vitro</u> culture of normal diploid fibroblasts has served as a model studying cellular system for senescence immortalization. Hayflick and Moorhead reported in 1961 that, with continuous passage, human diploid fibroblasts reach replicative senescence at a characteristic number of population doublings. Somatic cells derived from the

tissue of a young individual and grown in culture can divide a maximum of 50-100 times before reaching senescence. Furthermore, the upper limit in the number of cell divisions is inversely related to the age of the donor. Replicative senescence thus appears to be a genetically-programmed series of changes exhibited by normal cells that culminates in exit from the cell cycle and expression of a senescent phenotype.

As the body ages, the proportion of senescent cells within the skin increases. The accumulation of such cells is likely to have both direct and indirect effects that contribute to age-related changes and pathologies. As a cell becomes senescent, changes in the pattern of gene expression lead to functional changes. changes can then influence the physiology of surrounding cells by altering the extracellular environment or in a paracrine fashion through the release of proteins. For instance, the consequence accumulation of senescent cells within the skin is a progressive decrease in skin structure and function.

The present methods can be used in a number of applications to ameliorate the problems associated with the accumulation of senescent cells. For instance, the method can be used in the treatment of age-related lipofuscin accumulation in the retinas of humans and other animals. Lipofuscin accumulation correlates with Recently, Eldred and Lasky (1993) macular degeneration. characterized at least one of the components from the retinas of donors aged 52 to 98 years of age as Nretinylidene-N-retinylethanolamene, a Schiff base between retinaldehyde and ethanolamine. The source of ethanolamine was proposed to be the rod outer segments that are rich in phosphatidylethanolamine phagocytosed by the retinal pigmented epithelium (RPE). One theory for age related accumulation of lipofuscin is the decreased conversion of retinaldehyde to retinoic

acid in the aging RPE. The gene aldehyde dehydrogenase-1 (ALDH-1) is a senescence-related gene, as it is upregulated in young quiescent cells. Because senescent cells show a phenotype of M1 activation, this gene cannot be induced in senescent cells. One of the activities of ALDH-1 is the conversion of retinaldehyde to retinoic During cell senescence, the down-regulation of ALDH-1 can lead to increased retinaldehyde accumulation and therefore increased Schiff base formation with ethanolamine to form N-retinylidene-Nretinylethanolamine. The problem of decreased ALDH-1 in the formation of lipofuscin in aged RPE can be solved by therapy, i.e., by transfecting a recombinant constitutively expressed ALDH-1 into RPE cells to restore activity in aging cells and tissues, thereby decreasing the levels of lipofuscin. ALDH-1 activity could also be targeted in a screen for therapeutic agents that inhibit lipofuscin accumulation by increasing ALDH-1 activity. In addition, the gene for ALDH-1 could be added to aging cells by gene therapy .

A further use for the instant methods relates to treatment for AIDS. From recent studies, it appears that immune system dysfunction can result from the accumulation of senescent cells, as observed in AIDS patients and in the very old. There is evidence that senescent cells are refractory to apoptosis, and it has been reported that there is an increased number lymphocytes in aged mice that are refractory This correlates with decreased expression of When FAS is added transgenically, FAS in the aged cells. there is a marked restoration of normal immune function, even in old animals. It may be that the senescent cells that accumulate in aging humans and in AIDS patients are likewise refractory to apoptosis via a down-regulation of FAS levels or activity. If these long-lived, senescent lymphocytes have an inappropriate recognition of self, or

inhibitory signals arowth to prevent the proliferation of progenitor cells, then their presence could be deleterious. The instant method could be used to detect such senescent cells and treat the problems associated with their accumulation. For senescent cells can be eliminated via the administration of the FAS ligand or a similar molecule that is capable of activating the FAS receptor. This could stimulate apotosis in these cells, clearing them where appropriate, allowing younger cells to replenish the circulating pool. The administration of such agents could decrease the percentage of senescent cells and decrease the phenomenon autoimmunity or other immune system dysfunction. Molecules capable of activating the pathway via receptor or other means could be discovered by placing senescent T lymphocytes in a cell-based screen for agents that induce apoptosis.

The methods of the instant invention also include the sorting of cells to distinguish senescent cells from young cells. Young cells (cells which have undergone fewer divisions) can be isolated from the tissue of donors, using a method of the instant invention, and prepared for use, i.e., reintroduction to the same donor. The cells with the greatest replicative capacity can be isolated and then grown in a culture medium which slows the replicative senescence of these cells. One method of slowing senescence is by the lengthening of telomeres, as in "Methods and Reagents for Lengthening discussed Telomeres," PCT/US94/13130. WO 95/13383, hereby incorporated by reference herein.

One cell type of particular relevance to the above described method is the melanocyte, a dendritic cell that inhabits the epidermis and produces melanosomes for inclusion in keratinocytes. The melanocyte is thought to replicate along with the basal keratinocytes and to possess a finite replicative capacity comparable to that

of other somatic cells. Senescent cells produce more melanin than the comparable young cell. The loss of melanocytes through replicative senescence may account for the hypopigmentation observed in aged skin and the poor tanning response, and the highly pigmented patches known as solar lentigines, or liver spots, may represent patches of senescent cells that overexpress melanin. Senescent melanocytes from a donor can be separated from young cells using the instant methods, cultured with an agent which slows or reverses replicative senescence, and reintroduced to the donor. Alternatively, one could administer a pre-toxin that is activated only senescent cells.

These and many other therapeutic benefits can be realized with the present invention. The invention is illustrated by example below with respect to fibroblasts. To identify the senescence-related genes in fibroblasts, which may be responsible for or contribute to this decrease structure and function of aged skin, the cells used to illustrate the present method were BJ (foreskin) and IMR90 (lung) fibroblasts cultured in media containing either 0.5% or 10% serum. One could also fibroblasts derived from skin. such as fibroblasts isolated from fetal dorsal hand tissue. A series of 20 different 5'-primers and 12 different 3'-primers were used to amplify the mRNA, so 240 different primer sets were employed. The primers used in the amplification, together with their sequences and alpha-numeric designations, are shown in Table 1, below.

Table 1
Primers used in EDD for Human Fibroblasts

__3'-(T-rich)-primers:

```
5'-GCG CAA GCT TTT TTT TTT TTC T-3' (SEQ ID NO. 1)
    5'-GCG CAA GCT TTT TTT TTT TTC C-3'
B:
                                         (SEQ ID NO. 2)
    5'-GCG CAA GCT TTT TTT TTT TTC G-3'
                                         (SEO ID NO. 3)
    5'-GCG CAA GCT TTT TTT TTT TTG T-3'
                                         (SEQ ID NO. 4)
    5'-GCG CAA GCT TTT TTT TTT TTG G-3'
                                         (SEQ ID NO. 5)
    5'-GCG CAA GCT TTT TTT TTT TTG A-3'
                                         (SEQ ID NO. 6)
    5'-GCG CAA GCT TTT TTT TTT TTA T-3'
                                         (SEQ ID NO. 7)
    5'-GCG CAA GCT TTT TTT TTT TTA C-3'
H:
                                         (SEQ ID NO. 8)
    5'-GCG CAA GCT TTT TTT TTT TTA G-3'
J:
                                         (SEQ ID NO. 9)
    5'-GCG CAA GCT TTT TTT TTT TTA A-3'
K:
                                         (SEQ ID NO. 10)
    5'-GCG CAA GCT TTT TTT TTT TTC A-3'
                                         (SEQ ID NO. 11)
    5'-GCG CAA GCT TTT TTT TTT TTG C-3' (SEQ ID NO. 12)
   5'-(randomly-selected)-primers:
00: 5'-CGG GAA GCT TAT CGA CTC CAA G-3'
                                         (SEQ ID NO. 13)
01: 5'-CGG GAA GCT TTA GCT AGC ATG G-3'
                                         (SEQ ID NO. 14)
02: 5'-CGG GAA GCT TGC TAA GAC TAG C-3'
                                         (SEQ ID NO. 15)
03: 5'-CGG GAA GCT TTG CAG TGT GTG A-3'
                                         (SEQ ID NO. 16)
04: 5'-CGG GAA GCT TGT GAC CAT TGC A-3'
                                         (SEQ ID NO. 17)
05: 5'-CGG GAA GCT TGT CTG CTA GGT A-3'
                                         (SEQ ID NO. 18)
06: 5'-CGG GAA GCT TGC ATG GTA GTC T-3' (SEQ ID NO. 19)
07: 5'-CGG GAA GCT TGT GTT GCA CCA T-3'
                                         (SEQ ID NO. 20)
08: 5'-CGG GAA GCT TAG ACG CTA GTG T-3'
                                         (SEQ ID NO. 21)
09: 5'-CGG GAA GCT TTA GCT AGC AGA C-3'
                                         (SEQ ID NO. 22)
10: 5'-CGG GAA GCT TCA TGA TGC TAC C-3'
                                         (SEQ ID NO. 23)
11: 5'-CGG GAA GCT TAC TCC ATG ACT C-3'
                                         (SEQ ID NO. 24)
12: 5'-CGG GAA GCT TAT TAC AAC GAG G-3'
                                         (SEQ ID NO. 25)
13: 5'-CGG GAA GCT TAT TGG ATT GGT C-3'
                                         (SEQ ID NO. 26)
14: 5'-CGG GAA GCT TAT CTT TCT ACC C-3'
                                         (SEQ ID NO. 27)
15: 5'-CGG GAA GCT TAT TTT TGG CTC C-3'
                                         (SEQ ID NO. 28)
16: 5'-CGG GAA GCT TTA TCG ATA CAG G-3'
                                         (SEQ ID NO. 29)
17: 5'-CGG GAA GCT TTA TGG TAA AGG G-3'
                                         (SEQ ID NO. 30)
18: 5'-CGG GAA GCT TTA TCG GTC ATA G-3'
                                         (SEQ ID NO. 31)
19: 5'-CGG GAA GCT TTA GGT ACT AAG G-3' (SEQ ID NO. 32)
```

The amplification involved 4 cycles of degenerate (low temperature and fidelity) amplification, each cycle comprising 94°C for 45 sec.; 41°C for 60 sec.; and 72°C for 60 sec.; and 18 cycles of high temperature and fidelity amplification, each cycle comprising 94°C for 45 sec.; 60°C for 45 sec.; and 72°C for 120 sec. The amplified products were separated and visualized by polyacrylamide gel electrophoresis, and the differentially-displayed bands

were assigned a band number and then excised from the gel and either sequenced or cloned or both. An example of an EDD gel is shown in Figure 1 and examples of differentially displayed bands are shown in Figure 3. About 150 young and old-specific genetags of senescence-related genes were identified by this process. These genetags are summarized in Table 2, below.

The band number in Table 2 is a 4 digit number, the first two digits identify the 5'-primer used to generate the band, the third digit is a letter identifying the 3'primer used to generate the band, and the fourth digit was assigned according to the number of differentiallydisplayed bands in a particular lane on the gel. number in Table 2 reflects the cells and media conditions in which the mRNA corresponding to the band is observed, according to the formula: 01 is BJ senescent cells, 0.5% serum; O2 is IMR90 senescent cells, 0.5% serum; O3 is IMR90 senescent cells, 10% serum; Y1 is BJ young cells, 0.5% serum; Y2 is IMR90 young cells, 0.5% serum; and Y3 is IMR90 young cells, 10% serum. If available, Genbank locus designation is provided, and if the Genbank locus designation is not known, the term "Novel" is used to indicate that sequence information from the genetag is available but does not match any Genbank locus, and the term "Unknown" is used to indicate that sequence information is not yet available.

Table 2
Genetags Identified in EDD of Fibroblasts

Band No.	Age No.	Genbank Locus	Size (bp)
00C2	010203	HUMTIMPR	225
00D3	0102	Unknown	150
00E1	Y1	Unknown	135
00H1	03	Unknown	185
00K1	Y1Y2Y3	HUMC1A2	450
00L1	01	Unknown	450
00M2	Y1Y2	Novel	370
01C1	010203	HUMTPA	170
01C2	Y1Y2	Unknown	135

Table 2 (continued)

Genetags Identified in EDD of Fibroblasts						
5						
Band	Age	Genbank	Size			
No.	No.	Locus	<u>(qd)</u>			
01C3	01	Unknown	235			
01D1	01	HUMSECP3	230			
01E1	Y1	Novel	450			
01E2	Y1Y2	Novel	320			
01E4	Y1Y2	HSCOL3A1	137			
01M4	01 (0203)	HUMINFGAMM	235			
01M5	Y1	HUMBGR1A	145			
01M6	Y1Y2	Unknown	130			
02A1	Y1Y2(Y3?)	Unknown	>>500			
02 A 2	0102 (03?)	MIT1HS	285			
02B1	Y2	Novel	355			
02B2	Y1Y2	Unknown	175			
02C1	Y1Y2	HUMCG1PA1	200			
02E3	03	Unknown	148			
02M4	Y1	Unknown	140			
02M6	0203	Unknown	225			
03C1	Y1Y2	Novel	>450			
03C2	0102	Novel	380			
03E1	0102	Unknown	>400			
03F1	0203	Novel	>400			
03F2	02	Novel	200			
03J1	01	Unknown	250			
03J1r	01	Unknown	205			
03J3	Y1Y2	Novel	240			
03J3r	Y2	Unknown	190			
03J4	0102	Unknown	330			
03M1	0102 (03)	Unknown	245			
03M2	Y1Y2	Unknown	190			
03M3	010203	Novel	180			
04D3	Y1Y2	Novel	200			
04E2	Y2Y3	Unknown	180			
04E2	Y1	unknown	175			
04F2	0203	Unknown				
04F3	0203	Unknown	160			
	V1	Novel	170			
04L3 04M1	Y1Y2	Unknown	135 240			
05B1	Y2	Unknown	270			
05C1	Y1Y2	Unknown				
05C1	01	Unknown	450 350			
05C2 05C3	0102		350			
		Novel MnSOD	280			
05C4	01		255			
05D1	Y2	Novel	300			
05D2	Y1Y2	Novel	260			
05J1	0102 (03)	Unknown	160			
05J2	Y1Y2	T08744	>500			
05K1	Y1	Unknown	140			
06D1	010203	Novel	215			
06E1	010203	HUMTFPA	180			
06E2	01	Novel	150			

Table 2 (cont.)

Genetags Identified in EDD of Fibroblasts					
Band	Age	Genbank		Size	
No.	No	Locus		(qd)_	
06J1	01	Novel		190	
06L1	010203?	Unknown		>500	
06L2	01	Unknown		240	
07C1	03	Unknown		230	
07C2	Y1Y2	Novel		190	
07E1	010203Y3		(666)	140	
07J1	01020313	HUMIGFBP5	(000)		
07J2	Y1Y2Y3	Novel		>500	
0762 07L1	03	Novel		180	
07L2	V1	Unknown		220	
07B2	(0102)03			220	
		Novel		215	
08B1	02	Unknown		175	
08D1	01	Unknown		180	
08D2	01	Unknown		165	
08D3	Y1Y2	Unknown		130	
08D4	03	Novel		115	
08D5	01	M78570		90	
08E1	01	Unknown		285	
08E2	01	Novel		230	
08E3	0203	HUMSGP3		200	
08F1	01	M78570		295	
08L2	010203	Unknown		150	
08M1	Y1	Unknown		210	
08M2	02	Novel		130	
09B1	Y1Y2	Unknown		230	
09D1	010203	Novel		350	
09D2	010203	Novel		120	
09E1	Y1	T06399		155	
09E2	Y2Y3	Novel		145	
09J1	Y1Y2	Unknown		180	
10C1	010203	Unknown		220	
10D1	01	HUMPAI2B		>500	
10F1	(Y1)Y2	Novel		100	
10J1	01	Novel		190	
10M1	(01)02	Unknown		240	
10M2 10M3	Y1 Y1Y2	Unknown		155	
10M3 10M4	01	Unknown		140	
11B1	02	Novel Novel		115	
11E1	03	HUMCILA		285	
11E1	03			320	
11E2 11E3	01	Unknown Novel		225	
11K1	Y1Y2	Unknown		150	
11K1 11K2	010203?	Unknown		190	
11K2 11M1	02	Unknown		170	
11M1 11M2	010203	Unknown		240	
12F2	0203	HSCDN7		170	
13C1	Y1	Unknown		100	
1701	1.1	OTIVITOMII		125	

Table 2 (cont.)
enetags Identified in EDD of Fibrobles

Genetags Identified in EDD of Fibroblasts						
Band	Age	Genbank	Size			
No		Locus	(qd)_			
13D1	03	Unknown	185			
13F1	01	Unknown	130			
13M1	Y1Y2	Novel	450			
13M2	Y1Y2	Unknown				
14F1	Y1	Unknown	145			
14F2	0102	Unknown	130			
14M1	01	Novel	110			
15M1	0203	Unknown	175			
15M2	Y1	Unknown	100			
16B1	010203	Unknown	90			
16C1	Y1Y2		230			
16C2	Y1Y2	HSLAMA3	230			
		T09243	90			
16E1	Y2	Unknown	220			
16F1	Y1	Novel	200			
16F2	0102	HUMHERGC	160			
16F3	03	HUMCD44B	120			
16F4	Y1Y2	Unknown	110			
16H2	Y1Y2	Unknown	100			
16J1	Y1Y2	Unknown	125			
16K1	010203?	Unknown	280			
16K2	01	Unknown	190			
16L1	Y2	Unknown	170			
16M1	Y1	Unknown	250			
17B1	01 (02) 03	Unknown	330			
17F1	01	Unknown	320			
17F2	01 (0203)	Unknown	200			
17H1	01	Unknown	260			
17K1	Y1	Unknown	250			
17M1	Y1	HUMSPARC	180			
18C1	Y1	HUMALDHA1	330			
18D1	01	Unknown	290			
18H1	01	Novel	320			
18H2	Y2	Unknown				
18H3	¥2	Unknown	220			
18J1	01	Unknown	190			
18M1	Y1	Unknown	270			
18M2	¥2	Unknown	330			
18M3	Y1	HUMKCS	200			
18M4	01		140			
19E1	01	Unknown	110			
19E1	0203	Unknown	330			
19F1 19F1		Unknown	280			
	0203	Unknown	280			
19M1	Y1	Unknown	250			
19M2	Y2Y3	Unknown	140			

As demonstrated by Table 2, EDD of fibroblast cells resulted in the identification of many different

senescence-related genetags. Many of the genetags were from known genes, including both those known to be regulated with age and those not previously known to be regulated with age. For instance, IFN gamma has not previously been shown to be regulated with age. 00C2 corresponds to an mRNA with sequence homology to Genbank locus HUMTIMPR, which encodes TIMP (tissue metallo-proteinases) inhibitor of also known erythroid-potentiating activity glycoprotein and collagenase inhibitor, and may encode TIMP 2, which has been reported to be present at higher levels in senescent cells (see Zeng & Millis, 1994, Exp. Cell Res. 213: 148). In addition, genetags specific for the mRNA of PAI2 and tPA, which have been previously reported to be oldspecific gene products, and genetags for procollagen chains for type 1 and type 3 collagen, which have been previously reported to be young-specific gene products, were also identified by EDD.

A number of known genes were detected by several different primer sets, and it is expected that the same gene can be detected by different primer sets for the genes corresponding to the "Novel" and designations in Table 2. For instance, Genbank locus HUMTPA (encodes tissue plasminogen activator, also known as tPA) was identified using the primers sets defined by the following band numbers: 01C1, 01E3, 01F1, 02D1, and 03C3; Genbank locus HSCOL3A1 (encodes human pro-alpha 1 type 3 collagen) was identified using the primers sets defined by the following band numbers: 01E4, 01F2/3, 02C2/3/4, 02D2, 02E4, 02F3/4, 02H2, 02J1, 02K1, 02K2, 04C2, 04D2, 04L1, and 18K1 (and 01C2 may also be derived from this gene); and Genbank locus HUMCG1PA1 (encodes human pro-alpha I chain of type I procollagen) identified using the primers sets defined by the following band numbers: 02C1 and 02E2.

Other known genes for which genetags were identified using EDD of fibroblasts include HUMC1A2 (encodes human pro-alpha 2 chain collagen type 1, band no. HUMSECP3 (encodes human JE gene, which encodes a monocyte secretory protein, band no. 01D1, which comprises a sequence also in human interferon gamma, band no. 01M4); HUMBGR1A (encodes human glutamate receptor, band no. 01M5); MIT1HS (encodes mitochodrial RNA, band no. 02A2); HUMTFPA (encodes human tissue factor, band no. HUMIGFBP5 (encodes human insulin-like growth factor binding protein 5, band no. 07J1; band no. 11H1 corresponds to Genbank locus HUMIGFBP5X); HUMSGP3 (encodes human secretory granule core proteoglycan, also (HSHPCP) hematopoietic proteoglycan (HUMSERG) a serglycin gene, band no. 08E3); HUMPAI2B (encodes human PAI-2, band no. 10D1); HUMCILA (encodes human lipoprotein-associated coagulation inhibitor; this genetag also comprises sequences of Genbank HUMOS14E01, from a human HepG2 3'- directed Mb01 cDNA. clone s14e01, band no. 11E1); HUMHERGC (encodes human heregulin beta 2 gene, a specific activator of p185-erbB2 (see Science 256: 1205 (1992), band no. 16F2); HUMCD44B (encodes human cell adhesion molecule CD44, band no. 16F3); HUMSPARC (encodes human osteonectin, identified as aortic endothelial RNA, band no. 17M1); HUMALDHA1 (encodes human aldehyde dehydrogenase 1, band no. 18C1); HUMKCS (encodes human 80K-L protein, which is also known as calmodulin binding protein, protein kinase C substrate, band no. 18M3).

Other genetags could be correlated with sequences in Genbank that have not yet been associated with a known gene or function. Thus, the present invention provides for the first time a utility for synthetic oligonucleotides comprising these gene For instance, band no. 05J2 corresponds to sequences. Genbank locus T08744 from Expressed Sequence Tag (EST)

06636, which shares sequence homology with EST02797 and EST00675, as well as to the GOS2 gene and/or alcohol dehydrogenase and human suilisol mRNA (HUMSUIISO); band no. 07E1 corresponds to Genbank locus T03598 IB568 (666) from a human cDNA clone known as IB568; band nos. 08D5 08F1 correspond to Genbank locus M78570 from EST00718, a cDNA clone homologous to tubulin alpha; band corresponds to Genbank locus from EST04288, a human cDNA clone HFBDS91; band no. 12F2 corresponds to Genbank locus HSCDN7, a CDNA clone isolated using differential display as differentially expressed between androgen dependent and independent prostate carcinoma cell lines; and band nos. 16C2 and 16L2 correspond to Genbank locus T09243 from human sequence tag EST07136, a 3'-end clone HIBBR16. **EST** T09243 has now been identified in Genbank as human microfibril associated glycoprotein 4, HUMMFAPA, number L38486. The genetag 05C4, shown above in Table 2, has been identified as known gene Manganese Super Oxide Dismutase (MnSOD, see Kumar et al., 1993, Exp. Gerontol, <u>28</u>: 505-13).

Other genetags identified using this method share homology with known genes. For instance, band no. 00D3 has sequence homology with the CD44 gene; band no. 00M2 has some homology with the human aromatase cytochrome P-450 gene; band no. 03M3 has a rich G(1-2)A(1-3) stretch and so shares sequence homology with the human pepsinogen gene, the PAI1 gene, and some human sequence tags (the GA region seems to be a repeating motif in some genes); and band no. 16C1 share sequence homology with the human laminin A gene.

With respect to the genetags identified in Table 2 as "Unknown", these genetags can be cloned and then used as a probe for Northern analysis of RNA samples of young and old, mitotic and quiescent fibroblast cells, to verify that the genetags are differentially expressed. Once the

verification is made, the genetags can be sequenced and then identified as either a known or previously unknown Once identified, senescence-related genes, probes specific for those genes, the gene products of those genes, and antibodies to those gene products can be used as markers for detecting senescent cells, for distinguishing between and/or separating young and old cells, and for screening and/or therapeutic purposes. For screening purposes, the gene product of a senescent gene will be a useful target for therapeutic intervention if that gene product is involved in disease pathology or if a change in its expression parallels that of gene products involved in disease pathology. quantitate changes in the level of gene expression caused by a compound using high-throughput screening techniques. Using active compounds, one can determine whether and at what level coordinate modulation of gene occurs (i.e., does the compound affect senescence-related genes globally, in groups, or individually), and if by group, to which group an individual gene belongs. set of genes initially chosen for use in such screening can be modified as screening results accumulate.

Table 3 below summarizes the data collected from the EDD performed on fibroblast cells.

Table 3
Summary of EDD Results on Fibroblasts

Summary of EDD Results on Fibroblasts						
	K	nown	Novel EDD			
<u>Cells</u>	G	enes	Genetag	Bands	Total	
Young IMR90	0		3	9	12	
Young IMR90/BJ		collagen 1, pro a 2 collagen 3, pro a 1 (14x) collagen 1, pro a 1 (2x) laminin A EST06636 EST07136 (2x)	10	16	32	
Young BJ	5	aldehyde dehydrogenase glutamate receptor 80K-L protein osteonectin EST04288	2	13	20	
Old IMR90	5	lipoprotein-associated coagulation inhib. hematopoetic proteoglycan CD44 HSCDN7 tPA (3x)	6	13	24	
Old IMR90/BJ		tPA (2x) TIMP 2 human tissue factor (EST) IB568 heregulin mitochondrial RNA	7	15	28	
Old BJ	6	IGF binding protein 5 PAI-2 (uPAI) MnSOD interferon gamma (2x) EST00718 (2x)	9	18	32	
Total	27		37	84	148	

Most of the genes and gene products noted in Tables 2 and 3, above, have not previously been identified as products of senescence-related genes. Many of the gene products are secreted proteins, which is consistent with the alteration in extracellular matrix observed in aging tissues. Thus, the present invention provides novel methods and reagents for identifying senescent cells in tissue or culture, which methods generally comprise determining whether a cell expresses a senescence-related gene product, which can include an mRNA or other RNA or a

protein, and correlating the presence of that gene product with the state of senescence of the cell or tissue.

Typically, such methods will be practiced using oligonucleotide probe hybridization to the mRNA of the cell, either in situ or in a cell extract. In one such method, probes specific for the mRNA corresponding to a senescence-related gene are immobilized on a membrane or Then, the cells of interest are cultured under conditions conducive to gene expression and flash-frozen. The cells are then thawed in the presence of a labelled mRNA precursor, so that the label is incorporated into transcripts that were being transcribed when the cells were frozen. The labelled mRNA is then harvested from the cell and hybridized to the immobilized probes on the filter. The pattern of hybridization will identify whether senescence-related genes are being expressed by Those of skill in the art readily understand how to make probes specific for a particular gene product provided the sequence of the gene or mRNA produced by the gene is known.

Consequently, although the sequence of known genes is repeated herein, but the sequence of corresponding to novel senescence-related provided in Table 4, below. The sequences are identical to the RNA identified by EDD but for the substitution of deoxyribonucleotides for ribonucleotides and are identified by band no. and in the shown 5'-to-3' direction (N is any base and indicates that the identity of the nucleotide at that position is not known).

Table 4 Novel Sequences from EDD of Fibroblasts

00M2

CATTTATTCA TTCATTGAGA CACTCAA (SEO ID NO. 33)

01E1

ACAGAAAGGC CACTCAGGAT GTCCTTTGTG TCCATTGATG TCATTCAGCA GTCAGTCCCC CAATAATCCT TAAACTAGCT AAAACCAAAG GTAGTCNTTA GAAGATCTGC TT (SEQ ID NO. 34)

01E2

TTGAGTAGTT ACTGGAACCT TGACATTGCC TTTTAATGAG GTACTTCCAA AAAAAGGACC CCTAACAATG GCATAATAGT GAGGTCTCTC TGTGCGTGTA CATAATATA (SEQ ID NO. 35)

02B1

CAAAGATAAG AAACCAAGGA AGAAAGCAA (SEQ ID NO. 36)

03C1

CTGACGCCAN CCGCATACNC CGCANCCACA (SEQ ID NO. 37)

03C2

AGATAAAGCA ATTAGAAGAT GCATTAAAAG ATGTGCAGAA GAGGATGTAT GAGTCAGAAG (SEQ ID NO. 38)

03F1

ATAATAAAC TCTTCATTTT GCGAATTATA GAAGCTACTT TTTATAAAGC CATATTTTT TAGGGAAACT AAGGAGTGAC ATAGAA (SEQ ID NO. 39)

03F2

AACTGCATTT TGATGTTATC GCTTATGTTT AATAGTTAAT TCC (SEQ ID NO. 40)

03J3

CTATTGCCTC TCCTCCTGCA GAGACCATG (SEQ ID NO. 41)

03M3

GAGAAGAAAG GAAAGAAAGG NCACAGAGAT GGAAGGCCA (SEQ ID NO. 42)

04D3

GTTTCTGAAT TACATGAATT GTTGCAGAGC AAAGAAACTT ATGGAAATCT TTCCATTTAT (SEQ ID NO. 43)

04L3

GTAGGCTTCT ATATTGCATT TAACTTG (SEQ ID NO. 44)

05C3

AATGAGGTAG AAGTAGAAAG GAAGAAAAC TCAAAGAATT CTAAAAGGAT TCATAGCAAC ATAATGTGTC CC (SEQ ID NO. 45)

Table 4 (cont.) Novel Sequences from EDD of Fibroblasts

05C4

TCTCACATTC AGTCATACCC TAATGATCCC AGAAAGATAA TCAT (SEQ ID NO. 46)

05D1

AGAAGCCCCA GCAAGATTTA TTCCTTTTTG CTTCTTCATA ACCATGAAGC CATTGAAC (SEQ ID NO. 47)

05D2

CTACCTCCCA CATTAATTTT CATATGT (SEQ ID NO. 48)

06D1

AGGGCACAGC ACCAGATGAA TTGTTGTATA T (SEQ ID NO. 49)

06E2

AAATTAGCTT TCATCACAGA TTTAGGAAACT TGTCT (SEQ ID NO. 50)

06J1

AAACTACTGA ACNGTTACCT AGGTTAACAAC CCTGGTTGAG TATTTGC (SEQ ID NO. 51)

07C2

TTGNATATTG NATTTGTAGT AATATTCCAAA AGAATGTAAA TAGG (SEQ ID NO. 52)

07J2

AAATTGTATA TTGTATTTGT AGTAATATTCC AAAAGAATGT (SEQ ID NO. 53)

07L1

TATGAATNTC ACATTTGAAT TCTTCGATCTC TAA (SEQ ID NO. 54)

07M1

TATGTATAAA AGCATATGTG CTACTCATCTT TGCTCAC (SEQ ID NO. 55)

08D4

AATGTCTAAT TTTCTTTCCG ACACATTTACC AAA (SEQ ID NO. 56)

08E2

ACAACAGCAA ACAAAAAGGT GAAGTCTAAAT GAAGTGCACA (SEQ ID NO. 57)

08M2

AAAAGAATTG GCAGTTACAT TCATACTTT (SEQ ID NO. 58)

Table 4 (cont.)

Novel Sequences from EDD of Fibroblasts

09D1

AAGAATGTGC ATTCCAGTGC CATAGATAGT ATATTGAA (SEO ID NO. 59)

09D2

TTGCTACGGA CTTACGAAAG GACAAAGCGA AGAGCTG (SEQ ID NO. 60)

09E2

AAATAATTTA TTCATTGCAG ATACTTTTTA GGTTGCATTT TATTCATTTC C (SEQ ID NO. 61)

10F1

AGATGATGAT GTTAACCCAT TCCAGTACAG TATTCTTTT (SEQ ID NO. 62)

10J1

AGTATAGTGA ATGANTATGC CTTCCTACTG CAG (SEQ ID NO. 63)

10M4

AGAAATATAA AGATTTTNAT ACCTGCCACA TGG (SEQ ID NO. 64)

11R1

GAAGANATTA TGTTGTGANC NGGAGTNACA CAAA (SEQ ID NO. 65)

11E3

AGGGGCACAA GAGTTTGCGG TTATTGAATC CTGAGANAA (SEQ ID NO. 66)

13M1

GTTGAAGAGA CAGAGACAAG TAATTTGC (SEQ ID NO. 67)

14M1

CCGTGAATAC CCNTTTCTCG ACCAAAGA (SEQ ID NO. 68)

16F1

ATGGAGTTGT GGATGAAAGC CATGTTAGNT G (SEQ ID NO. 69)

18H1

GATCATATAA ACANNNCCGA GTTCTACCTC AGAGTCG (SEQ ID NO. 70)

Those of skill in the art will recognize that the complete coding sequence of a gene corresponding to a genetag of the invention, as well as the endogenous promoter and other regulatory elements of the gene, can be readily isolated once the sequence of the genetag or

the genetag itself (which can be generated using the primers indicated above) is known or available, provided by the present invention. Such genes can be either directly or after suitable modification using standard techniques of molecular biology, not only to express the mRNA or protein encoded by the gene but also to express antisense oligonucleotides or ribozymes that can be used to prevent deleterious expression of senescence-related genes. Those of skill in the art recognize that a wide variety of expression plasmids can be used to produce useful nucleic acids of the invention and that the term "plasmid", as used herein, refers to type of nucleic acid (from a phage, chromosome, etc.) that can be used to carry specific genetic information into a host cell.

To verify that the genetags are differentially expressed, probes for known and novel genes were prepared and tested on Northern blots. Northern and Southern blot analysis was performed according to standard procedures (Sambrook, J., Fritsch, E.F., & Maniatis, T. Molecular Cloning; a laboratory manual, 2d ed., Spring Harbor Laboratory Press, Cold Spring Harbor). Nucleic acids were transferred to charged membranes, then cross-linked using a UV STRATALINKER (Stratagene). Probes for novel genetags were prepared by restriction digestion of the appropriate plasmid, which contains the novel genetag cloned into it, with HindIII followed by isolation of the insert band on low-melting agarose and then radioactive labeling of the DNA using the random hexamer-primed method (Feinberg & Vogelstein, 1983, Anal. Biochem. 132:6-13). Following hybridization to the probe, the filter blot was washed twice in 1x SSC and 0.5% SDS at 65°C for 30 minutes. The blots were then analyzed using a PhosphorImager exposed and 425E (Molecular Dynamics). Probes for known genes prepared by kinase end labeling of approximately 40-mer

oligos that were designed to be complementary to the 5' coding region of the message. The following probes for known genes were used:

80K-L protein, Human, GenBank Accesion Number (GBAN): D10522, bases: 497-458, 5'-CCG TTT ACC TTC ACG TGG CCA TTC TCC TGT CCG TTC GCT T-3' (SEQ ID NO.71);

Aldehyde Dehydrogenase 1, Human, GBAN: K03000, bases: 186-147, 5'-AGG AAC AAT ATT CAC TAC TCC AGG AGG AAA CCC TGC CTC T-3' (SEQ ID NO. 72);

Cell Adhesion Molecule (CD44), Human, GBAN: M59040, bases: 166-127, 5'-CCG AGA GAT GCT GTA GCG ACC ATT TTT CTC CAC GTG GTA T-3' (SEQ ID NO. 73);

Collagenase, Human, GBAN: X05231, bases: 181-142, 5'-CCA GGT ATT TCT GGA CTA AGT CCA CAT CTT GCT CTT GTG T-3 (SEQ ID NO. 74)';

Collagen 1 alpha 1, Human, GBAN: K01228, bases: 990-951, 5'-ATC AGC ACC TTT GGG ACC AGC ATC ACC TCT GTC ACC CTT A-3' (SEQ ID NO. 75);

Collagen 1 alpha 2, Human, GBAN: J03464, bases:519-481, 5'-AAG GTT ACT GCA AGC AGC AAC AAA GTC CGC GTA TCC ACA A-3' (SEQ ID NO. 76);

Collagen 3 alpha 1, Human, GBAN: X14420, bases: 142-103, 5'-CGA GAA GTA GCC AGC TCC CCT TTT GCA CAA AGC TCA TCA T-3' (SEQ ID NO. 77);

Elastin, Human, GBAN: M36860, bases: 183-144, 5'-CCT GGA TAA AAG ACT CCT CCA GGA ACT CCA CCA GGA ATG G-3' (SEQ ID NO. 78);

EPC-1/PEDF, GBAN: M76979, bases: 256-217, 5'-AAG AAA GGA TCC TCC TCC ACC AGC GCC CCT GTG CTG T-3'
(SEQ ID NO. 79);

Heregulin-beta 2, Human, GBAN: M94167, bases: 654-615, 5'-GAG GAG TAT TCA GAA CTG GTT TCA CAC CGA AGG ACT AGT T-3' (SEQ ID NO. 80);

Human Tissue Factor, GBAN: M16553, bases: 357-318, 5'-CTC GTC GGT GAG GTC ACA CTC TGT GTC TGT TGT GTA AAA C-3' (SEQ ID NO. 81);

IGF Binding Protein 5, Human, GBAN: M62403, bases: 538-499, 5'-CTG GTG CTC CGG TCT CGA ATT TTG GCG AAG TGC TTC TGC A-3' (SEQ ID NO. 82);

Laminin A, Human, GBAN: X70904, bases: 193-154, 5'CTC CAT ATT GAT AGG CGT GCT CTA TTG CTC TAG GGC TGT T-3'
(SEQ ID NO. 83);

PAI-1, Human, GBAN: M16006, bases: 340-301, 5'-TCT TGA ATC CCA TAG CTG CTT GAA TCT GCT GCT GGG TTT C-3' (SEQ ID NO. 84);

PAI-2, Human, GBAN: M18082, bases: 74-34, 5'-ATT GAG GGC AAA GAG TGT GTT TGC CAC ACA AAG ATC CTC C-3' (SEQ ID NO. 85);

TPA, Human, GBAN: X13097, bases: 299-259, 5'-TTG CTT CTG AGC ACA GGG CGC AGC CAT GAC TGA TGT TGC TG-3' (SEQ ID NO. 86); and,

UPA, Human, GBAN: K03226, bases: 405-366, 5'-ATC TGT GGG CAT GGT ACG TTT GCT GAA GGA CAG TGG CAG A-3' (SEQ ID NO. 87).

Genetags that could be readily sequenced were used for further analysis of the data generated by Oligonucleotide probes for Northern analysis were designed that were complementary to many of the known genes listed in Table 3 (see above). Five probes were prepared for genes that were previously characterized as being differentially regulated in young and senescent cells (PAI-1 (Goldstein, et al., 1994, <u>J. Cell Physiol.</u> 161: 571-9), elastin (Fazio, M.J., et al., 1988, Lab. <u>Invest.</u> <u>58</u>: 270-7.), EPC-1 (Pignolo, et al., 1993, <u>J.</u> Biol. Chem. 268: 8949-8957), collagenase (West, et al., 1989, Exp. Cell Res. 184: 138-147.) and urokinase type plasminogen activator (Shay, et al., 1992, EXD. Gerontology 27: 477-492.); see Figure 4B). RNA was prepared from both young and senescent IMR90 and BJ fibroblasts grown in either 0.5% or 10% FBS. These 8 RNA samples allow a comprehensive analysis of the expression

of genetags identified by EDD. Northern analysis of previously characterized genes with altered expression in young or senescent cells was performed for all 5 genes and the reported differential expression was confirmed (see Figure 4B). For reasons of convenience, the 8 RNA samples were loaded for Northern analysis in the order shown in Figure 4A. Of the 21 known genes identified by EDD (not including ESTs), Northern analysis confirmed differential expression for 12 genetags (see Figure 4C), including the four differentially displayed bands from Figure 3, while 7 genetags failed to show the predicted differences (data not shown). The hybridization pattern of one probe was non-specific, and one gene (mitochodrial RNA) was not analyzed in this analyis.

During the cloning of novel genetags, observed that a band that generated an identifiable sequence sometimes gave rise to plasmids with different In such cases, several clones for each band were isolated and a clone that contained the initial sequence was used to probe Northern blots. Of the 37 novel genetags (Table 3), 31 were recovered after cloning. Two genetags were identified twice in EDD, leaving 29 unique genetags to be tested. The conditions for specific hybridization of probes for novel genetags were first tested in Southern blot analysis. are most often generated from the untranslated 3' end of the message, a region rich in repetitive-type sequence. Of the 28 probes that detected specific bands in Southern analysis, 5 probes did not give a signal in Northern analysis (20 µg of total RNA was loaded per lane), 12 probes did not confirm the EDD observations, probes (see Figure 4D) were in agreement with the initial EDD observation. This analysis again demonstrates that about half of the genetags tested confirmed the initial EDD observation.

Novel genes which were shown to be differentially expressed using the instant methods are shown in Table 5 below, which also shows the approximate size of mRNA detected in Northern blots as well as the size of insert cloned into the genetag cloning vectors. The reference numbers for each genetag are also shown.

Table 5

Geron	Approximate	CDNA	Note
Reference	mRNA Size	Insert	
		Cloned	
09D2	2.9 kb	2.9 kb	identical to 08E2
05C4	2.9 kb	2.4 kb	identical to MnSOD2
08E2	2.9 kb	2.9 kb	identical to 09D2
10J1	5.3;3.6 kb	3.9 kb	
10M4	3.9 kb	NA	
11E3	3.3;2.7 kb	2.5 kb	
18H1	6.0;3.8 kb	4.3 kb	
03J3	14;11;3.8 kb	3.8 kb	
04L3	3.7 kb	3.5 kb	
16F1	3.7 kb	2.0 kb	
07L1	3.5 kb	2.0 kb	identical to 08D4
10F1	5.4;3.6 kb	5.4 kb	
05D1		6.5 kb	

Those in the art will recognize that the given mRNA sizes in the above table are only approximations. Figure 4D shows examples of Northern blots of mRNA similar to blots from which the above listed mRNA sizes were obtained.

Sequence data available for these senescence-related novel genes is given below in Table 6; the sequences are also identified by Geron reference number. The term "contig" identifies a sequence of a genetag for which another noncontiguous sequence is also available. The

sequences are shown in the 5'-3' direction. N is any base and indicates that the identity of the nucleotide at that position is not known. M indicates that the nucleotide at that position is either A or C. R indicates that the nucleotide at that position is either A or G. W indicates that the nucleotide at that position is A or T/U. Y indicates that the nucleotide at that position is C or T/U. K indicates that the nucleotide at that position is G or T/U.

Table 6

Senescence-Related Novel Gene Sequences

08E2, AT7

GCGGCGGCCA TGGCGGGACA GGAGGATCCG GTGCAGCGGG AGATTCACCA GGACTGGGCT AACCGGGAGT ACATTGAGAT AATCACCAGC AGCATCAAGA AAATCGCAGA CTTTCTCAAC TCGTTCGATA TGTCTTGTCG TTCAAGACTT GCAACACTAA ACGAGAAATT GACAGCCCTT GAACGGAGAA TAGAGTACAT TGAAGCTCGG GTGAC (SEQ ID NO. 88)

08E2, contig 1

CCAGCAATCT ATCATGGATC CTAATCAGAA CGTGAAATGC AAGATAGTTG TGGTGGGAGA CAGTCAGTGT GGAAAAACTG CGCTGCTCCA TGTCTTCGCC AAGGACTGCT TCCCCGAGAA TTACGTTCCT ACAGTGTTTG AGAATTACAC GGCCAGTTTN GAAATCGACA CACAAAGAAT AGAGTTGAGC CTGTGGGACA CTTCGGGTTC TCCTTACTAT GACAATGTCC GCCCCCTCTC TTACCCTGAT TCGGATGCTG TGCTGATTTG CTTTGACATC AGTAGACCAG AGACCCTGGA CAGTGTCCTC AAAAAGTGGA AAGGTGAAAT CCAGGAATTT TGTCCCAAAT ACCAAAATGC TCTTGGTCGG CTGCAAGTCT GATCTGCGGA CAGATGTTAG TACATTAGTA GAGCTCTCCA ATCACAGGCA GACGCCAGTG TCCTATGACC AGGGGGCAAA TATGGCCAAA CAGATTGGAG CAGCTACTTA TATCGAATGC TCAGCTTTAC AGTCGGAAAA TAGCGTCAGA GACATTTTTC ACGTTGCCAC CTTGGCATGT GTAAATAAGA CAAATAAAAA CGTTAAGCGG AACAAATCAC AGAGAGCCAC AAAGCGGATT TCACACATGC CTAGCAGACC AGAACTCTCG GCRGTTGCTA CGGACTTACG AAAGGACAAA GCGAAGAGCT GCACTGTGAT GTGAATCTTT CATTATCTTT AATGAAGACA AAGGAATCTA GTGTAAAAAA CAACAGCAAA CAAAAAGGTG AAGTCTAAAT GAAGTGCACA GCCAAAGTCA TGTATACCAG AGGCTTAGGA GGCGTTTGAG AGGRTACTCA TCTTTTTGGG AATCCTGACC TTAGGTTCGG CATGTAGACC AAGTGATGAG AAGTGAATAC ATGGAAGAGT TTTTAAGTGT GACTTGAAAA ATATGCCAAA AAATGAGAGA TACAAATGAG CTAGAGGAAG ATGAGGGGGG ATGCGAGTAC CTCCAAGAAG

AAAAATCACA CTCTGAATGG TGCTTGCATT TTGAGGTTTN
NNNNCNNNNG GGGTATAATC TATCATGGAT CTCCACTTGG
ATTAATTTTA AATGTTTAAT CTCCTTTACA AAAAGTATAC
GTTAATATAC CGTCCTCAAG GGGGACTGGC ACTGTGACCT
AGCATTAGTT TCTAGAGGAT GTGATCTAAT TCTTCTAGCT
CATCATAAAA AGGAATTGTA TCAGGACCCA TGGGATATAT
CCAGAGGCAA CTTATGAGGC TTGGAATCTG GCTTCCTGAA
GATAGCTGAG TAGGATGGTC TAAGGAAAGC CTTGGAATCT
TGCAAGATTG GTGGACCAGC ACTACAAAGA TCGCATAGAT
CAAATAGGAA AAAAATGTCG ATTTTTATTC AGTCTGATGG
TTCTGTTCTC ATGGTGGATG GTCATAAAAA GTGG
(SEQ ID NO. 89)

08E2, contig 2

GGCACGAGCG GAATTGGACT TGGGAGGCGC GGTGAGGAGT CAGGCTTAAA ACTTGTTGGA GGGGAGTAAC CAGCCTGCTC CTCTCGCTC TGCGCCGCGT TTCAGAGGTT GCCCATCAGC CTTGTGATTT ATTTTATAT CTGCTTTTTA TAAAGAGAGA AATATATATA TATATATATA TATTTTTTT TTCTTCTTAA GAGAAAATTC CTGTTCCAAG AGAAAATAAG GCAACATCAA TGAAGGAGAG AAGAG (SEO ID NO. 90)

08E2, contig 3

08E2, DT7

TCCTCTCGCT CTCCTCCTCG TCTGCGCCGC TTTCAGAGAG AAAATTCCTG TTCCAAGAGA AAATA (SEQ ID NO 92)

10J1

ACCATGGAAA GTATAGTGAA TGAATATGCC TTCCTACTGC
AGCAAAACTC AAAAAAGCCC ATGACAAATG AGAAACAAAA
TTCCATTTTG GCCAACATTA TTCTGAGTTG TCTAAAGCCC
AACTCCAAGT TAATTCAACC ACTTACCACG (SEQ ID NO. 93)

10J1, 3T3

GTTTTTTTG AGTTTAACAC AGATTTTATT GCCCTATAGA CAGTTATGAT GTGACCAGTG GATATCAATG AAACTTCTTA ATTATTTGAG TCTGAAAATG CATATTTA (SEQ ID NO. 94)

10J1, 3T7

ACAAACCACA GTATTTCCAC TTTAAATATA GAACTGGTAA ACAGCACTAT CCTTAAACTA AAATCGGTGA CTCGGTCATC AAGAAGGTTT TTGCCCGCCC GTGGATC (SEQ ID NO. 95)

10J1, 4T3

TTTTTTTTT TTTTTTTTT TTTTTTTT GAAAGAATAG GTTTAATTTA TTAGTTGCTC TTTAGCAAAG GCTATATAGA ACA (SEQ ID NO. 96)

10J1, 4T7
GCTAAACCAA ACCAACTCCT CTGCTTTGTC CCTTAGGTGG
AAAAGAGAGG TAGTTTAGAA CTCTCTGCAT AGGGGTGGGA ATTAAT
(SEQ ID NO. 97)

10M4

TGTTATTATT AGCTTTCCAT GTGGCAGGTA TTAAAATCTT TATATTTCTG AAAAATACTC TTCTCTT (SEQ ID NO. 98)

PCT/US95/11230

11E3

GATGGGGACG TCCTGATTTA CCAGATCAAA GTATGGTAAG GCTGTAGATA GCACACTAGT TTTCTCAGGA TTCAATAACC GCAAACTCTT GGTGCCCCTA TTG (SEQ ID NO. 99)

11E3, 4T3

TTTTTTTTT TTTTTTTT TTAGGAAAGA AAGAATTCTT TTATTCCCAC ATGACAGCCC AATTTTTTAA AATGGTTATC TTAAGTCAGG CCAGTTTTAT TTTATTGACC ATGTATATAT AACATCAGAT ATTTCTAAGA AAGAGAAGAG AACCTGATTG ATGTCTCTCA TGT (SEQ ID NO. 100)

11E3, contig 1

CCCACCGGGG CCATGGCGTC CGCCGACATC GTCGTGGGCG
GGGTGGCCCA CGGGCGGCCC TACCTCCAGG ATTATTTTAC
AAATGCAAAT AGAGAGTTGA AAAAAGATGC TCAGCAAGAT
TACCATCTAG AATATGCCAT GGAAAATAGC ACACACAA
TAATTGAATT TACCAGAGAG CTGCATACAT GTGACATAAA
TGACAAGAGT ATAACGGATA GCACTGTGAG AGTGATCTGG
GCCTACCACC ATG (SEQ ID NO. 101)

11E3, contig 2

18H1, 1T3

TTTTTTTTT TTTTAAACTT AAAAATCAGT TTATCTTAAA ATTTTTTTGA ATTAGAAGAA ATTGATTTCA CATGAAAATA TAACATTCTG AAGGTAAGTA TTTTAAACAT TACAACTGTT (SEQ ID NO. 103)

18H1, 2T3

TGAAAATGCC AAAAATAGTG AACTTAACCT TGTTGGCTAT
GATCTGCTAC ACCAGCCAAC ATTATTGTGC CTTTGCATTT
CACACCAAAA GTTCCAAATG GGTATTTTTT GATGATGCAA ATGTGAAAG
(SEQ ID NO. 104)

18H1, 3T3

18H1, 5T7

CATGTTTGTC AGGGAGATGC CTGTATATTT TGTGCATTGA AGACGATATT TGCACAGTTC CAACACAGTC GAGAAAAAGC ACTTCCCTCA GATAACATAA GGCATGCTCT TGCAGAAAGT TTCAAAGATG AGCAGCGATT TCAACTTGGC CTTATGGATG ATGCTGCGGA GTGCTTTGAA AATATGTTGG AGAGGATTCA TTTTCACATA GTGC (SEQ ID NO. 106)

18H1, contig 1

TAAAAATACT ATCAACATTT AATGTTTCTC TCATAGGTTT CCAACCTCG TTCCGGCTTT TACTGCTGCT ATCACAGCTG TTTCCTCTAT CCCTAGAATC TTGGCTGCTG TTCTGTGTCA TATCCAGTGC CATTGTCACT CTTAACTTTG CCAGTTATTT TCTCTCCTGG AGCAAGAATC TGGGATTTAC TTGAACTTAT TATTTGTGCA GAAGCTMGAC TYTGAGGTAC AACTCGGTCR TGTTTATATG ATCCTTTTCC TTGA (SEQ ID NO. 107)

03J3, 1T3

GCGCGCCGC AGCTGCAGGC TCAGGCCATC TCCAGTCTTC CTGCCGCGCC CGCCTCTCCT CCTGGGATTC TCCTCCTCCT CCTGGACTTC CCCGCAGCTG CCCACGCCTC GGCGGCCGCC AGTGCTCCTG GAGTGCAGA (SEQ ID NO. 108)

03J3. 2T3

CTTTTGCGTC GCCAGCCGAG CCACATCGCT CAGACACCAT GGGGAAGGTG AAGGTCGGAG TCAACGGATT TGGTCGTATT GGGCGCCTGG TCACCAGGCT GCTTTTAACT CTGGTAAAGT GGATATTGTT GCCATCAATG ACCCCTTCAT TGACCTCAAC TACATN (SEQ ID NO. 109)

03J3, 2T7

03J3

CTGCCAATTC TTGCTCAATT CTAAGGCAGA TAGACTGTGT GAGTTCAAAA GATATCCTCT GAAAAGCATC AAAATCTTCC ACTGTGAACA CATGGGTCTC TGCAGGAGGA GAGGCAATAG CTTCCAATTC TGAGCGAACG GCATCCTC (SEQ ID NO. 111)

04L3, 11T3

GCAGGGGCTG AAATAACCCA GATGCCCCCA CCCTGCCACA TACTAGATGC AGCCCATAGT TGGCCCCCCT AGCTTCAGCA GTCACTAT (SEO ID NO. 112)

04L3, 17T3

CTGAAATGCA GCTCCCTGTC CAAGTGCCTT GGAGAACTCA
CAGCAGCACG CCTTAATCAA AGGTTTTACC AGCCCTTGGA
CACTATGGGA GGAGGGCCAA GAGTACACCA ATTTGTTAAA
AGCAAGAAAC CACAGTGTCT CTTCACTAGT CATTTAGAAC
ATGTTATCAT CCAAGACTAC TCTACCCTG (SEQ ID NO. 113)

04L3, 24T3

CCTCCGCTTA CAGCTCGCTG CCGCCGTCCT GCCCCGCGCC CCCAGGAGAC CTGGACCAGA CCACGATGTG GAAACGCTGG CTCGCGCTCG CGCTCGCGCT GGTGGCGGTC GCCTGGTCCC GCCGAGGAAG ACTAAGACAA ATCAAGATCT T (SEQ ID NO. 114)

04L3, contig 1

TTTTTTTTT TTWWTTTYTT TATTCATCAA TAGTATCCGA AAAGGAAGAA TCAGGAGTTA CAAAAACAAG TTAAATGCAA TATAGAAGCC TACTAAATAC AAATACAAGT CACAAACACA TATGCAAGAG AAACTTGTTT AGATTG (SEQ ID NO. 115)

16F1, 1T3

CCGGGCCGGG GAGGCGCGCT CGCTCCGCGC TCCCTTCGCT CGCTCGTTTC CTCCTCCTC GGCAGCCGCG GCGCAGCAG GAGAAGGCGG CGGCGGCGC TAGGGATCAG ACATGGCGGC GGATCTGAAC CTGGAGTGAT CTCCT (SEQ ID NO. 116)

16F1, 1T7

TTTTTTTTT TTTTTTTAA GAACATCAAC ATTTATTTAA CATGATAAAA AAAGAAATGA GATATGAACA TTTGCATTTA AACAATAGTA AGTAGCCTTT AATACATTAC ATGTGCTCAT TGTATAATAT ATACACAATG AACATAATTA CATTTGTACA CAAACTAAGT ACCGGATTTG GAAACCTGCT TATTGCTGTA CACATGTATT CCAATG (SEQ ID NO. 117)

16F1, 3T3

CACCACGATC AAAAGGGACA AGCATCAAGC ACGCAGCAAT
GCAGCTCAAA ACGCTTAGCC TAGCCACACC CCCACGGGAA
ACAGCAGTGA TTAACCTTTA GCAATAAACG AAAGTTTAAC
TAAGCTATAC TAACCCCAGG TTGGTCAATT TCGTGCCAGC
CACCGCGTCA CACGATTACC AAGTCATAGA G (SEO ID NO. 118)

16F1, contig 1

TTTTTTTTT TTTTTTTK YTYWTTTYYT TGTTTTTTG
ATTGTTTTGG ATCTCTGGTT TAATTAGCAC TCTATGGTTG
GGAATGTTAT TGGTTTCTTT AGTTGGTGCA TTTTCAGATG
TAATCTTGTC CACTCTTYTC ACAGGTTCTG TCTGTACTAG
GGCAGCATCT AACATGGCTT TCATCCACAA CTCCATTTCC
TTTCCTGTAT CAGTGCAGAA AAAAGG (SEO ID NO. 119)

16F1, contig 2

CGGATCTGAA CCTGGAGTGG ATCTCCCTGC CCCGGTCCTG
GACTTACGGG ATCACCAGGG GCGGCCGAGT CTTCTTCATC
AACGAGGAGG CCAAGAGCAC CACCTGGCTG CACCCCGTCA
CCGGCGAGGC GGTGGTCACC GGACACCGCG GCAGAGCACA
GATTTGCCTA CTGGCTGGGA AGAAGCATAT ACTTTT
(SEQ ID NO. 120)

07L1, 2T3

AGATGGTCCA GCTGCCAGGA CTACTTTGGC AGGCAGCGTG
CTACAGGACG AAAATGTAAG AGAAGTCTAT TAAGGCTGGA
CAGCCCAGGG TTATTTATAC TCTCTCAGCC CCAAGTCCCC
CGGACTAAAG ACCTAAAGGC TGATTGACTC ATTCCTGATT
GATTTAATGG AAAGTCTCCC ACCCCATCAT CATTTGCCAG AGTAC
(SEQ ID NO. 121)

07L1, 2T7

CAAAACTCAG CAGTGCTTCT GGTGCTGGTG ATCAGTGCTT CTGCAACCCA TGAGGCGGAG CAGAATGACT CTGTGAGCCC CAGGAAATCC CGAGTGGCGG CTCAAAACTC AGCTGAAGTG GTTCGTTGCC TCAACAGTGC TCTACAGGTC GGCTGCGGGG CTTTTGCATG CCTGGAAAAC TCCACCTGTG ACACAGAT (SEO ID NO. 122)

07L1

TTAGATGCAA GAAGATGCAG GCTCAAAGTC TGGTTGGACA GCCAGGCTCA AGCAATTTGG TAAATGTGTC GGAAAGAAA TTAGACATTG GAGGATCAAG ACCATAAGAC ACTAGCTCAT TAGAGATCAA GAATTCAAAT GTGACATTCA TATTCGTCC (SEQ ID NO. 123)

PCT/US95/11230

10F1

ATTTTAAAAG AATACTGTAC TGGAATGGGT TAACATCATC ATCTTTGGCA TCCT (SEQ ID NO. 124)

10F1, 2T3

GACATTCGCC CTGATATAAA AGATGATATA TATGACCCCA CCTACAAGGA TAAGGAAGGC CCAAGCCCCA AGGTTGAATA TGTCTGGAGA AACATCATCC TTATGTCTCT GCTACACTTG GGAGCCCTGT ATGGGATCAC TTTGATTCCT ACCTGCAAGT TCTACACCTG GCTTTGGGGG GTATTCTACT ATTTTG (SEQ ID NO. 125)

10F1, 5T3

TCGCACTTTG CCCCTGCTTG GCAGCGGATA AAAGGGGGGCT GAGGAAATAC CGGACACGGT CACCCGTTGC CAGCTCTAGC CTTTAAATTC CCGGCTCGGG GACCTCCACG CACCGCGGCT AGCGCCGACA ACCAGCTAGC GTGCAAGGCG CCGCGGCTCA GCGCGTACCG GCGGCTTCG AAAC (SEQ ID NO. 126)

10F1, 5T7

TTTTTTAAAT TACAACACTT TATTGCAGCA TCGGCAAAGG TCAGATTTCT GAAGCTGGTG AAGATTGGGC AGCATTTCCA TGTGAAATGT TACAACTTTA CAAGTTTTGT TTTTTATTTA AATCTACATG CAGAAACTGA AACATGGTAA AAGAAAAAAT GCAAAATAGC TAGAAAAAAA GATGTAATCA AGTTGTCGCA TACAGATGTG CTCTCCG (SEQ ID NO. 127)

10F1, 6T3

CTTGTCTGAC CTAGAAGCTG AGAAACTGGT GATGTTCCAG AGGAGGTACT ACAAACCTGG CTTGCTGATG ATGTGCTTCA TCCTGCCCAC GCTTGTGCCC TGGTATTTCT GGGGTGAAAC TTTTCAAAAC AGTGTGTTCG TTGCCACTTT CTTGCGATAT GCTGTGGTGC TTAA (SEQ ID NO. 128)

10F1, contig 1

TTTTTTCAGA TTCACTTCAC TTTTATTATG AACAAACACA ATCTCAGATT AGTACAATTA GCTTCAGAGT TGATATTAAT AGAAATTATT CCAAAATTAT TCTTGTCACA AGTAACTACT ATATCCCACA TAAAAAGGGA AAAAATCCCA CCCAATCACA GAAAAGGCAT CCTCTGTATG TTTCCGTGGC AATGCGTTGT TTATGTATTC TCAAATTTTG TCTGGCTAGT TATC (SEQ ID NO. 129)

05D1

TTGGTGAGGG TGGGGAGGAG TGCTTCGAAG GGAGAAGCCC CAGCAAGATT TATTCCTTTT TGCTTCTTCT TCTCCCTGTC CCTGCCATAA CCATGAAGCC TTGAACAAAC CACCCAAATC TCAGGATCTT AGTGTTTTCT CTGTAAATTG TAATATGAAC TTATAAAGAT CCTCCATTGC TGATAGTCTC AGGTTCTGTG AGTAACAGCA AAAAAACTTT GTATCTAACT TCAACCAGAG CAGGCTGTAC CCTTAAGCTC T (SEQ ID NO. 130)

The sequences shown above can be used in gene therapy to produce RNA directly in cells or tissues, as nucleic acid probes in diagnostic methods, as nucleic acid primers, and as components of recombinant DNA cloning and/or expression vectors. Further, those of skill in the art will recognize that the complete coding sequence of a gene corresponding to a genetag of the invention can be isolated once the sequence of the genetag is known. Thus, those in the art could use the above sequence information to isolate the coding sequence of the gene and to express the mRNA or protein encoded by the gene. Diagnostic and therapeutic applications for an expressed mRNA and corresponding protein are known in the art and are applicable to the practice of the present invention. For example, antisense oligonucleotides or ribozymes can be designed to target the mRNA, and antibodies can be generated to target the protein. Such antisense oligonucleotides, ribozymes, and antibodies can be used for the detection of senescence-related genes or to prevent the deleterious expression of such genes, described above.

As noted above, probes and/or primers comprising the sequences shown in Table 4, or those shown in Table 6, or sequences from other senescence-related genes identified according to the methods of the present invention can be used in diagnostic methods to detect the presence of young or old cells in a tissue or other sample. Several probes that were identified to be young or senescent

specific by Northern analysis have been used in in situ analysis of human skin tissue sections. These probes can be used to identify senescent cells or young cells in tissue (in vivo; in situ). Probes for Collagen 1α1. Collagen $1\alpha3$, and for MnSOD were used in a commercially available in situ hybridization protocol (SureSite™ II System Manual, Novagen). Preliminary results show that the probes for Collagen $1\alpha1$ and Collagen $1\alpha3$ detect RNA in dermal fibroblasts in human skin sections of two young donors (<15 years age), while probes for MnSOD detects RNA in dermal fibroblasts in human skin sections of two old donors (>69 years age). These results confirm the expectation that young specific genes detect cells in tissue from young, but not normal old, donors and that senescent-specific genes detect cells in tissue from old, but not normal young, donors. These findings support the belief that replicative senescent dermal fibroblast cells accumulate with age in human skin.

This in situ data illustrates the value of the nucleic acid primers and probes of the invention. Primers and probes are oligonucleotides that complementary, and so will bind, to a target nucleic acid. Although primers and probes can differ in sequence and length, the primary differentiating factor is one of function: primers serve to initiate DNA synthesis, as in PCR amplification, while probes are typically used only to bind to a target nucleic acid. Typical lengths for a primer or probe can range from 8 to 20 to 30 or more nucleotides. Modified, synthetic, and/or non-naturally occurring nucleotides can also be used in whole or in part in the oligonucleotides of the invention. or probe can also be labelled to facilitate detection (i.e., radioactive or fluorescent molecules are typically

used for this purpose) or purification/separation (i.e., biotin or avidin is often used for this purpose).

Depending in part on the length and/or intended function of the primer, probe, or other nucleic acid comprising sequences from a senescence-related gene, expression plasmids of the invention may be useful. For instance, recombinant production of RNA corresponding to a genetag or senescence-related gene of the invention can be carried out using a recombinant DNA expression plasmid of the invention that comprises a nucleic acid comprising the nucleotide sequence of the genetag positioned for transcription under the control of a suitable promoter. Host cells for such plasmids can be either prokaryotic or eukaryotic, and the promoter, as well as the other regulatory elements and selectable markers chosen for incorporation into the expression plasmid will depend upon the host cell used for production.

One important use of probes derived from the genetags and corresponding genes of the present invention relates to a method for screening compounds to identify compounds that can alter gene expression in senescent cells, which method comprises: (a) contacting senescent cells with a compound; (b) determining mRNA expression patterns in said senescent cells; and (c) correlating an alteration in mRNA expression of a senescence-related gene with identification of a compound that can alter gene expression in senescent cells. Preferably, determination of mRNA expression pattern involves determination of mRNA expression levels of two or more Thus, this screening method senescence-related genes. identifies compounds with the capacity to partially reverse, or modulate the pattern of gene expression that is altered as a consequence senescence. The present invention also encompasses the compounds identified by this method and the use of those compounds to alter gene expression in senescent cells.

Such screening can also identify compounds that activate young-specific genes or prevent cells from entering a this senescent state. In method, the novel oligonucleotide probes of the invention serve as indicators of whether a test compound can alter the expression levels of a senescence-related gene.

Compounds ideally suited for testing in this method include compounds identified in primary screens based on the expression of a specific senescence-related gene In general, the basic format of the screen is product. as follows. Senescent cells are cultured in 96-well microtiter plates. After an incubation period, i.e., three days in culture, the medium will be removed and assayed for one or more senescence markers, providing a "before treatment" baseline. The medium will be replaced with fresh medium containing a test agent or its vehicle. The cells will be cultured for an additional two to four days or more in culture, presence of the test agent. The cells and/or medium will then be assayed for the senescent markers treatment" measurement). Samples that normalize the senescence markers can be presented to cultures of young cells and their effects measured in a similar fashion. Compounds that act selectively on the senescent cells will proceed to additional screening.

As noted in Table 2 and in the scientific literature, a number of known genes are senescence-related genes. For instance, the activity of ß-galactosidase is elevated in senescent fibroblasts. (See, Cristofalo, et al., 1975, Mech. Ageing Dev. 4:19-28; and Dimri et al., 1995, Natl. Acad. Sci. USA, 92 (in press)). Consequently, one can first conduct a primary screen of compounds to determine whether that inhibits ß-galactosidase activity in senescent cells. one embodiment of this screen, fibroblasts are grown to senescence, plated in 96 well plates, and incubated with

At the end of the incubation period, a test compound. analyzed enzyme for activity colorimetric assay based on the ability of the enzyme to cleave a colorless substrate into a colored reaction Compounds identified in this screen product. Example 14, below) as active compounds will then be tested in a secondary assay to determine that the active compounds are inhibiting the senescence-specific increase of activity of the enzyme and not merely inhibiting the enzyme itself. Other primary screens can be conducted using the senescence-related genes identified in Table 2, identified according to the methods invention, or known from the scientific literature. For instance, one could conduct a primary screen to identify compounds that have the capacity to induce the downregulation of collagenase activity, an enzyme that is known to be elevated in senescent fibroblasts.

Compounds active in these screens can then be tested according to the screening method of the invention to determine whether the compound inhibits the expression of other senescence-related, specifically old-related, genes or activates the expression of young-related genes, or both. The method can employ Northern analysis to examine the effects of the lead compounds on panels of genes that show altered expression or abundance in senescence as indicated by EDD. Based on the results of this screen. one can determine which compounds normalize expression of those genes that are altered in senescence and believed to contribute to age-related pathologies. Furthermore, it will be possible to determine the level at which the compound acts to reverse the pattern of altered expression. Complete reversal to a young pattern of gene expression would suggest that a single common mechanism is involved. Reversal of defined groups of genes would indicate that several mechanisms operating and that each is affecting a different bank of

genes. A compound may also act to modulate the activity of individual genes, suggesting the absence of a common mechanism. This information will in turn influence primary screening strategy. If, for example, all active compounds seem to reverse the altered expression of batteries of genes, or of only individual genes, then the screen can be expanded so that many more markers, including members from each of the putative batteries, if appropriate, are included.

Screens such as those described above can be based on many different cell types. For example, several types of skin models are useful for such screening methods. Examples of useful skin models include the wounded dermal model (see, Genever, et al., 1993, Exp. Dermatol. 2: 266-273), keratinocytes proliferated and differentiated on collagen sponge containing fibroblasts (see, Maruguchi et 1994, Plast. Reconstr. Surg., 93: 537-544), cultured skin substitute composed of fibroblasts and keratinocytes with a collagen matrix (see, Kuroyanagi, et al., 1993, Ann. Plast. Surg., 31: 340-351), reconstituted non-crosslinked collagen, or reconstituted collagen chemically crosslinked with glutaraldehyde, aluminum alginate or acetate (see, Middelkoop et al., 1995, Cell Tissue Res., 280: 447-453), and human keratinocytes grown on polyglactin mesh-cultured fibroblast substitute (see, Hansbrough, et al., 1993, J. Burn Care Rehabil. 14: 485-494).

In addition, the gene aldehyde dehydrogenase 1 (ALDH-1), found by EDD to be a quiescent- upregulated gene, can also be used in a screen, ALDH-1 is believed by those in the art to have a 17 beta dehydrogenase activity. The down-regulation of this gene in senescent cells could therefore lead to altered steroid metabolism (i.e., more active androgen). Since the cells in the prostate, hair follicle, and hair sebaceous gland respond to androgen,

this increased androgen activity could lead to benign prostatic hyperplasia, male pattern hair loss, and sebaceous gland hyperplasia (conditions known to be agerelated). The gene ALDH-1 can be used in a screen to find agents that inhibit the oxidase, but not the reductase activity of the enzyme, thereby allowing the enzyme to more efficiently metabolize steroids present in cells.

Using cultured cells for screening requires that a number of technical challenges be met. First, the cells must be kept viable during the screening. Second, the metabolism of the cells must not be perturbed by the assay conditions. Third, particularly in a multiple-day it is vital that sterility be maintained. screen, Culture conditions that will produce old senescent and young quiescent and mitotic cells must be carefully selected. taking into consideration the following criteria for senescence: (1) the cells typically will exhibit a change in the morphology that is characterized by the enlargement and flattening of the cell as it reaches senescence; (2) the cell will irreversibly leave the cell cycle and will be incapable of proliferation for a minimum of three weeks as measured by the population doubling time of the cells in culture; and (3) the cells typically will exhibit a nuclear labeling index below 1%, as measured by the incorporation of a labeled DNA precursor in the nuclei over a period of 24 hours (see Example 1, below). As the cells approach senescence, their generation time increases. This means that the interval between passaging cells is constantly varying and must be determined by continuously monitoring cell density. The time at which cell replication ceases (replicative senescence) must be determined.

Cell-based screens have traditionally been labor intensive and so have not often been used for high-throughput screening. However, the present method is

amenable to high-throughput screening. Liquid handling operations can be performed by a Microlab 2000™ pipetting station (Hamilton Instruments). equipment needed for the screen (e.g. incubators, plate washers, plate readers) can either be adapted automated functioning or purchased as automated modules. Movement of samples through the assay will be performed by an XP^{TM} robot mounted on a 3m-long track (Zymark).

Through these screens, libraries of synthetic organic compounds, natural products, peptides. and oligonucleotides can be evaluated for their capacity to modulate the expression of genes that reflect contribute to the disease process. Specifically, compounds can be identified that will down-regulate genes that are up-regulated during senescence or, conversely, will increase the expression of genes that are downregulated during senescence. Active compounds can be optimized, if desired, via medicinal chemistry. Initially, one can define a pharmacophore(s), modern computational chemistry tools, representative of the structures found to be active in the high throughput Once a consensus pharmacophore is identified, screens. design focused combinatorial one can libraries compounds to probe structure-activity relationships. Finally, one can improve biopharmaceutical the properties, such as potency and efficacy, of a set of lead structures to identify suitable compounds clinical testing.

Thus, the present invention provides novel methods for identifying senescence-related genes, methods and reagents for identifying senescent and young cells and for distinguishing senescent from young cells in tissue, and compounds and therapeutic methods for treating diseases and conditions resulting from cell senescence targeting senescent cells. The following examples describe specific aspects of the invention to illustrate

the invention and provide a description of the present methods for those of skill in the art. The examples should not be construed as limiting the invention, as the examples merely provide specific methodology useful in understanding and practice of the invention.

Example 1

Determining Mitotic Index of Cells with Immunohistochemical Staining

____ This assay allows one to quantitate the fraction of cells in S phase of the cell cycle via detection of incorporated 5-bromo 2-deoxyuridine into DNA and can be used to determine when cells are senescent.

__A. <u>Preparation and Treatment</u>

- 1. Grow putative senescent cells on sterile coverslips (Corning No. 1, 18 mm sq.) to 60% confluence. Allow cells to "recover" for one day before treating them with 5-bromo-2'-deoxyuridine (BrdU). Do not wash cells just prior to incubation with BrdU; this will slow the growth of the cells during the incorporation phase of the procedure.
- 2. Add BrdU (Sigma #B-5002) at a final concentration of 10 mM to cells in growth media. Keep the BrdU shielded from light during the addition step. Incubate the cells for a defined time period (typically 2 to 24 hours) in the dark at 37°C. One should have 1 mM stocks of BrdU dissolved in PBS and aliquoted into light protective tubes (or wrapped in foil protected from light) already made prior to performing the assay.
- 3. After the BrdU treatment, aspirate media into a waste bottle containing bleach, and wash the cells 3 times with PBS. After the third wash, add 3 ml of PBS to each well.

B. Fixation

1. The plates should contain 3 ml of PBS in each well and be kept on ice. Add 3 ml of fix solution, which is ice cold methanol:glacial acetic acid (3:1), to each well (on ice).

- 2. Remove 3 ml from each well and discard. Add 3 more ml of fix solution to each well.
 - 3. Repeat step 2 two more times.
- 4. Remove 5 ml from each well and add 6 ml of fix solution to each well. Leave the plates on ice for 15 minutes; then, repeat this step once again.
- 5. Remove coverslips from plate and allow to air dry overnight or for several hours. The solvents should be completely evaporated. Store coverslips in a covered box in the dark for staining the next day. If coverslips are to be kept longer than one day before staining, store the coverslips frozen in a dessicator under an atmosphere of N_2 .

C. Immunohistochemical Staining

- 1. Treat fixed cells with 0.01 N NaOH for 3 minutes to denature DNA and expose the antigen. Do not overtreat cells with base; the cells will fall off the coverslips.
- 2. Neutralize the base with PBS at pH=8.5. Wash the cells two times with PBS to remove all of the base.
- 3. Block the cells with 1.5% Horse Serum (Vector Labs #S-2000) in PBS (use serum from species in which secondary antibody was generated) for 15 minutes.
- 4. Carefully aspirate or tilt coverslips to drain horse serum from cells. Do not wash the coverslips at this point.
- 5. Add 400 ml of anti-bromouridine monoclonal antibody (Sigma #B-2531, made in mouse, IgG) at a dilution of 1/500 in 1% BSA, 0.05% Tween 20, in PBS, to cover cells completely. Incubate for 2 hrs. at room temperature in a humid chamber away from light.

- 6. Rinse cells in PBS three times.
- 7. Add 400 ml of secondary biotinylated horse antimouse IgG (Vector Labs #BA-2000) for 30 min. at room temperature in a humid chamber free from light. The antibody should be diluted in 1.5% horse serum in PBS at 10 mg/ml.
 - 8. Wash cells three times with PBS.
- 9. Incubate cells with 400 ml of 30 mg/ml Fluorescein Avidin D (Vector Labs A-2001) in 10 mM HEPES, 0.15 M NaCl, pH=8.5, for 20 min. at room temperature in a humid chamber protected from light.
 - 10. Wash the cells with PBS 3 times.
- 11. Add 4',6-diamidino-2-phenylindole, DAPI, (Sigma #D-9542) at a concentration of 1 mg/ml in 2X Standard Saline Citrate (SSC) buffer for 5 minutes to stain nuclei. Standard Saline Citrate is prepared by first making 20X SSC as follows: dissolve 175.3 g NaCl and 88.2 g sodium citrate in 800 ml of water; adjust pH to 7.0 with 10 N NaOH; and adjust volume to 1 liter with water; then, dilute 20X SSC to 2X SSC by adding 10 ml 20X SSC to 90 ml of distilled water. It is convenient to have stock solutions of DAPI prepared and stored in light protected tubes at 1 mg/ml dissolved in distilled water.
 - 12. Aspirate DAPI, and rinse cells 3 times with PBS.
- 13. Mount coverslips on slides using Vectashield™ (Vector Labs, #H-1000) mounting medium to reduce quenching of FITC. Do not use too much of the mounting medium, as it will not dry completely. Seal coverslips to slides using clear nail polish, and allow to dry for 5 minutes.
 - 14. View slides under the fluorescent microscope.

D. <u>Calculation of Mitotic Index</u>

The mitotic index is equal to the ratio of the number FITC labelled cells to the number of DAPI labelled cells in a given field.

E. Staining with Acridine-Orange to Identify Senescent Human Diploid Fibroblasts

In another method for staining cells to detect senescent cells the stain acridine orange is used. Acridine orange (AO) is a fluorescent dye which has unique spectral properties. In this method, growing cells are stained by adding AO in aqueous solution directly to culture medium and incubating the cells at 37°C (5% CO_2) for approximately 30 minutes. The cells are then viewed by epifluorescence microscopy using a 488nm excitation filter, 520 long-pass barrier filter ("FITC Filter Sets"). The fluorescence emitted from senescent cells differs from that of replicatively young cells in that (1) there is an increase in total fluorescence in the cytoplasm and a shift in emission spectrum towards longer wavelengths with senescent cells, and (2) there is a shift in emission spectrum towards longer wavelengths in the nuclei of senescent cells. Further, there is a change in the pattern of fluorescence in the nuclei of senescent cells which is consistent with the structural change in chromatin organization with senescence. pattern can be diagnostic of senescent cells. staining method can be applied to the analysis of tissues in situ.

Example 2

Cell Culture and RNA Preparation

Human fibroblast cells are split at the appropriate density according to standard tissue culture techniques. Cells are grown in DMEM medium plus 10% bovine calf serum (BCS). The last split before RNA isolation is 1 to 8 for young cells and 1 to 2 for senescent cells. After the cells are split and seeded in DMEM medium plus 10% BCS, two protocols are followed.

(1) If mitotic cells are required, cells are grown at 37°C in DMEM medium plus 10% BCS for two days. Then, RNA is isolated (see below).

(2) If quiescent cells are required, the DMEM medium plus 10% BCS is aspirated 4-8 hours after the seeding, when the cells have attached. The medium is replaced with DMEM medium plus 0.5% BCS, and the cells are grown for 3 days at 37°C. The medium is changed with fresh DMEM medium plus 0.5% BCS, and the cells are grown for 2 more days at 37°C.

For RNA isolation, cells are quickly washed once with PBS; then, the wash is aspirated completely, and about 1.5 to 2 ml of GITC solution (200 ml of GITC solution are prepared by adding 94.53 g of guanidine isothiocyanate to 1.67 ml of 3 M sodium acetate (pH=6), adding DEPC-water to 200 mL, sterile-filtering (0.22 μ m filter), and adding, in a fume hood, 1.67 ml of beta-mercaptoethanol) is added to a 15 cm plate. The cells lyse in this solution. After a few minutes of rocking the plate back and forward, to cover all areas, the (slimy) lysate is collected and prepared for CsCl₂ centrifugation.

Example 3

Spinning Guanidine Isothiocyanate (GITC) RNA Isolation

- ____ Before beginning this procedure, make sure the RNA+GITC and CsCl are at room temperature. The CsCl solution is prepared by adding 95.97 g of CsCl to 0.83 ml of 3 M sodium acetate, pH=6, and adding water to a final volume of 100 ml.
- 1. Use Ultraclear 14 x 89 mm polycarbonate tubes (Beckman 34059).
- 2. Add 4 ml of 5.7 M CsCl to the bottom of the tubes.
- 3. Gently, without disturbing the interface, add 7 ml of the Guanidine Isothiocyanate + RNA to these tubes.

4. Balance the tubes, and then place the tubes in the swinging buckets that are part of the SW41 rotor.

- Spin at 32K for 20 hours at 20°C.
- 6. After the spin, aspirate the supernatant (RNA is pelleted at the bottom of the tube).
- 7. After inverting tubes and letting them drain for about 30 minutes, cut off the bottoms of the tubes.
- 8. Resuspend the pellets in 180 μ l of diethylpyrocarbamate-treated deionized H,0 (DEPC-water).
- 9. Then do a second wash of the tube bottoms with another 180 μ l.
- 10. Pool these washes (total of 360 μ l) and add to them 40 μ l of 3 M sodium acetate, pH=6 (filter-sterilized, made with DEPC-water).
- 11. Vortex and add 1 ml of cold 100% ethanol, mix, and place at -80°C for at least 30′. To precipitate RNA, use 2.5 volumes of ethanol.
- 12. Spin for 30' at 14K at 4°C. Carefully aspirate off the ethanol. Don't let the tubes evaporate, because this will leave salt deposits behind.
- 13. Resuspend in DEPC-water. Decide on volume based on the pellet size.
- 14. Take the O.D. (260), i.e., 1 μ l in 100 μ l and use a 100 μ l quartz cuvette.

To check the RNA for quality control, one can run a gel, using the following procedure.

- 1. Heat 1-2 μg of RNA in 10 μl of DEPC-water at 70°C for 2 minutes. Incubate at room temperature, and add a standard loading dye to the sample.
- 2. Run a 1.2% agarose gel made with 1xTAE buffer made with DEPC-water.
- 3. Run the gel hot and fast, i.e., 150 volts for a $5" \times 7"$ gel or 75 volts for a mini-gel.

Example 4

Enhanced Differential Display

A. Synthesis of first-strand cDNA

The annealing reaction is conducted by mixing 1 mg total RNA with 2.5 ml of 20 μ M 3'-primer (dT₁₂ mer); and 9.5 ml of DEPC-H₂O. The resulting solution is heated for 10 minutes at 75°C, then cooled on ice for 7 minutes, and then spun to collect the mixture at the bottom of the tube.

The elongation reaction is conducted by adding to the tube 5 ml of 5X first strand synthesis buffer; 1 ml RNAsin (Promega or Pharmacia); 2.5 ml of 0.1 M DTT; 2.5 ml of 0.25 mM dNTP; and 1 μ l of reverse transcriptase (SuperScriptTM II RT, BRL). The resulting solution is incubated for 70 min. at 37°C. Then, the solution is heated to inactivate the enzyme by incubating the mixture for 10 min. at 95°C. The reaction mixture can be stored at -20°C for later use.

B. PCR amplification of cDNA

The reaction mixture is prepared using 1 ml of cDNA (3' primer carried over from cDNA); 2 ml of 10X PCR buffer (500 mM KCl, 100 mM Tris pH 8.3, 20 mM MgCl); 1.5 ml of 0.1 mM dNTP; 1.25 ml of 20 mM 5' primer; 1 ml of a 1 to 5 dilution of a-32P dATP; 0.5 ml of Taq polymerase; and 12.75 ml H₂O. About 70 ml of mineral oil are layered on top of the reaction mixture, which is briefly centrifuged to collect the reaction mixture at the bottom of the tube.

The PCR machine (Perkin-Elmer) is programmed to conduct 4 cycles of 94°C for 45 sec.; 41°C for 60 sec.; and 72°C for 60 sec., and then 18 cycles of 94°C for 45 sec.; 60°C for 45 sec.; and 72°C for 120 sec. The tubes are centrifuged briefly to collect the reaction mixture at the bottom of the tube and can be stored at 4°C.

C. Differential Display Gel Analysis

___ 1. Mix 3 ml of PCR product with 2 ml of running dye (formamide dye).

- 2. Heat samples for 3 min. at 90°C, pulse spin, and load on a 6% sequencing gel (see part D) in 0.6X TBE and run gel at 2000 V (current of ~50 mA).
- 3. Run the gel until the second dye reaches the bottom (this can be varied depending on what size range of bands one wants to compare).
- 4. Dry gel, and when the gel is dried, tape the gel and the film together, punch holes at the three corners of the gel, and expose to film overnight.

D. Sequencing gel

To prepare the sequencing gel, mix 36 grams of urea with 11.25 ml of 40% acrylamide/bis solution (19:1) and 4.5 ml of 10X TBE, and add DEPC-H₂O to 75 ml, and allow the components to go into solution. Then, filter the mixture through Nalgene 100 ml disposable filterware (CN), and add 330 ml of 10% ammonium persulfate and 33 ml of Temed, and pour the gel immediately.

Example 5

Removing Differentially Displayed Bands

- 1. Line film up with gel.
- 2. Poke holes through the film and gel on each side of the band with an 18 gauge needle.
- 3. Draw a line between the dots on the gel representing the band to be removed (this makes it easier to see the saran wrap in later steps).
- 4. Cut band out with a clean razor blade and manipulate band with clean forceps (rinse both items with water in between retrieving each band).
 - 5. Use a new razor blade every six bands.

6. Place band into a 1.5 ml Eppendorf™ tube and add 1 ml of TE buffer. Let soak for 15 minutes.

- 7. Aspirate the buffer and separate the saran wrap and paper from the gel slice.
- 8. Then add back 1 ml of TE buffer and aspirate again immediately (to dilute out urea in the gel slice).
 - 9. Add 40 μ l of TE buffer + 100 mM NaCl
- 10. Heat for 10 minutes at 95°C (boiling water bath).
 - 11. Let tubes cool overnight at room temperature.
- 12. Pulse spin at 14K and remove 5 $\,\mu l$ for PCR amplification.
- 13. Add 5 μ l 10X PCR buffer, 2.5 μ l 1 mM dNTP, 3 μ L 20 μ M 5'-primer, 3 μ L 20 μ M 3'-primer, 1 μ l Taq polymerase, DEPC-water to 50 μ l.
- 14. Run PCR (Perkin-Elmer machine) for 25 cycles, each cycle consisting of 94°C for 45"; 60°C for 1'; and 72°C for 1'; and then incubate at 70°C for 15'.

Example 6

PCR Sequencing Differentially Displayed Bands

- 1. Heat gel slice at 95°C to liquefy.
- 2. Remove 3.5 μ l and place into a 1.5 ml tube containing 1.5 μ l of the appropriate 20 μ M 3'-primer.
- 3. Add 5 μl of the dideoxynucleotide termination mix to the wells in the microtiter dish.
- 4. Make up a cocktail containing 10x sequencing buffer, ^{32}P -alpha-ATP, Taq polymerase, and water.
- 5. Add 18 μ l of this reaction cocktail to the tube containing the PCR-amplified band and primer.
- 6. Heat at 95°C for approximately 20" immediately before adding 5 μl of this cocktail to the appropriate termination mixes in the microtiter dish.
 - 7. Overlay with 20 μ l of mineral oil.
- 8. Add a drop of mineral oil to each of the wells in the PCR machine before inserting the microtiter dish.

9. Program the machine to conduct a 95°C soak for 5' and 30 cycles, each consisting of 95°C for 30"; 60°C for 30"; and 72°C for 1 min. Check to be sure that the sample probe heats up quickly enough.

- 10. When finished add 5 μ l of stop mix.
- 11. Denature samples in PCR machine by soaking 5 minutes at 95°C. Immediately load samples on the gel.

Example 7

Cloning Differentially-Displayed Bands in Bluescript SK+

A. Digestion of bands with HindIII

Choose 4 old- or young-specific bands to clone. Digest 5 μ l of each solution of band DNA, previously PCR-amplified from an acrylamide gel slice of a DD gel. The reaction mixture comprises: 5 μ l of band DNA, 5 μ l of 10X restriction buffer B (Boehringer Mannheim); 39 μ l of deionized H₂O; and 1 μ l of restriction enzyme HindIII (10 U/ μ l, Boehringer Mannheim) in 50 μ l total reaction volume at a temperature of 37°C for 2-3 hrs. After digestion, heat-inactivate the enzyme by incubation at 70°C for 20 minutes.

B. Preparation of Bluescript vector for cloning

It is useful to prepare a stock solution of $\sim 50~\mu g$ of pBluescript HindIII-digested (Stratagene) concentration of $\sim 0.25 \, \mu \text{g/ml}$, which can be stored in the freezer and used as a stock to aliquot from when preparing fresh calf-intestinal alkaline phosphatasetreated (CIPed) vector every 2-3 weeks. If more HindIIIdigested pBluescript has to be prepared, digest 30-50 µg of pBluescript with HindIII. Do not place all of this DNA tube, but rather use several tubes. in one containing 5 µg in a reaction volume of 20 µl. Digest at 37°C for at least 3 hours to ensure that digestion is

complete. Digestion must be complete, or too many blue colonies will appear after transformation. After digesting, run 1 ml of each digest on a 1% agarose gel to verify that digestion is complete. If digestions are not complete (i.e., if you see an additional band indicative of supercoiled DNA), combine all of the digests in one tube, add 5-10 μ l of HindIII (do not add any more buffer) and incubate at 37°C for 2-3 hrs. Check 2-4 μ l on a gel to ensure again that digestion is complete.

Prepare fresh CIPed HindIII-digested pBluescript SK+ every 2-3 weeks by reacting 20 μ l, i.e., ~5 μ g, pBluescript ($\underline{\text{Hin}}$ dIII-digested) with 1 μ l CIP (1 U/ml, Promega), 3 μ l of 10X CIP buffer (Promega), 6 μ l of deionized H_2O in a total volume of 30 μ l at 37°C for 1 hr. Add EDTA to a final concentration of 5 mM (i.e. add 0.5 μl of 0.31 M EDTA) and incubate at 70°C for 30 min. and then phenol extract the solution once. Then, add 1/10 volume of 3 M NaAcetate (pH=7.0) and 2 volumes of 100 % EtOH. Place on dry ice for 20 min. or overnight at -20°C. Spin down the DNA in a microfuge by 10 centrifugation, and wash with ~200 μl of 70% EtOH (ice cold). Do another 5 min. spin, and dry the DNA pellet in a dessicator or speed-vacuum centrifuge. Resuspend the DNA (~5 μ g) in 25 μ l of deionized water (dH,0) to a final concentration of ~0.2 $\mu g/ml$. Check 1 μl of this solution on a 1% gel to make sure the DNA is recovered in good yield. Use ~1 µl per ligation reaction.

C. Gel Purification of HindIII-digested band and ligation into HindIII-digested. CIPed vector

After digesting bands with <u>HindIII</u> and heat-inactivating the enzyme, load the entire digestion reactions on a 2% low melting point agarose gel. To prepare the gel, follow the procedure below.

Note that an $8" \times 10"$ gel box that will house a 30 well comb typically holds 300 ml of gel mix.

- 1. For a 2% gel, add 6 g of seaplaque low melt agarose to 300 ml of 1x TBE buffer. Note that 1% is 1 g/100 ml of gel.
- 2. Heat the solution in the microwave at level 7 for about 4 min. When finished, place the solution in the 65°C water bath to cool to a temperature suitable for handling the flask.
- 3. Add 3 μl of 5 mg/ml EtBr to the gel, mix by shaking, and add 4 μl to the 1x TBE running buffer.
- 4. Load .5 μl of the ØX-174 RF DNA digested with <u>Hae</u>III (Pharmacia) on both sides of your samples.
- 5. Run gel between 100 and 150 volts for about 1.5 to 2 hours.

Note: when adding loading dye to samples, add 1 μl of dye to 5 μl of sample.

Load 2 μ l of the 123 bp DNA ladder (BRL) as a marker. Take a picture of the ethidium bromide-stained, and UV-irradiated gel, and confirm that the size of the bands on the gel matches the sizes indicated on the data sheets. Cut out bands using coverslips and place the cut out bands into Eppendorf^m tubes. To remove the bands from the gel, follow the procedure below.

- 1. Clean a 365 nm UV light box by wiping the surface with ethanol.
- 2. Place the gel on the UV light box, and take a photograph using a hand-held camera with the aperture set on 8 and the time (which the shutter stays open) on B (hold trigger down for approximately two seconds).
- 3. Examine the photograph to make sure the bands are the proper size.
- 4. Cut bands out using glass coverslips, and remove with an ethanol-sterilized spatula.

5. Place the band in a labelled 1.5 ml Eppendorf™ tube, and discard the glass coverslip, in addition to resterilizing the spatula.

One should wear a face shield, lab coat, and gloves when cutting out the bands to protect against UV exposure.

Set up ligation reactions by mixing 2 μ l of Band DNA (melt down agarose at 65°C, 10 min. before adding) with 1 μ l of CIPed, <u>HindIII</u>-digested Bluescript (~0.2 μ g/ml), 2 μ l of 10 mM ATP, 2 μ l of 10X One Phor-All^m buffer (Pharmacia), 12 μ l of dH2O, and 1 μ l of T4 DNA ligase (Pharmacia) in a total volume of 20 μ l, and incubate at 37°C for 2-3 hrs.

D. Transformation of subcloning efficiency DH5alpha competent cells

- ___ To prepare the competent cells (DH5alpha cells are available from BRL), follow the procedure below.
 - Grow an overnight culture in 3 ml of media.
- 2. Place all of the overnight culture into 500 ml of LB media.
- 3. Allow the culture to grow to an optical density of ~ 0.4 .
- 4. Place all of the culture into 200 ml disposable centrifuge tubes (conical).
- 5. Centrifuge the cells at 4°C for 5 min. at 2000 rpm.
- 6. Combine and resuspend the cell pellets in 100 ml of cold 50 mM $\rm CaCl_2$.
- 7. Aliquot the resuspended cells into six 50 ml conical tubes that have been pre-chilled on ice.
 - 8. Spin at 2000 rpm at 4°C for 5 min.
- 9. Gently resuspend the cells in 20 ml of cold 50 mM $CaCl_2$ plus 15% glycerol.
- 10. Aliquot 100 μl of the resuspended cells per tube (enough for 1 transformation).

- 11. Flash freeze the tubes on dry ice.
- 12. Store the flash-frozen tubes at -80°C.

To conduct the transformation, follow the procedure below.

- 1. Remove competent cells (DH5alpha) from a -70°C freezer and place the cells on ice to thaw.
- 2. Place Eppendorf^m tubes on ice to chill. When thawed, add 200 μl of competent cells to each of the chilled tubes. The number of tubes corresponds to the number of transformations. Refreeze any unused cells in a dry ice/ethanol bath for 5 min., and return to the -70°C freezer.
- 3. Add all of each ligation reaction (20 μ l) to each of the tubes containing 200 μ l of cells. Mix gently.
- 4. As a control to test transformation efficiency, also add 5 μ l of 0.1 ng/ μ l (i.e. 0.5 ng) of control pUC19 DNA provided by BRL to 50 μ l of competent cells.
- 5. Incubate cells on ice for 30 minutes. Gently mix cells after 15 min. of incubation.
- 6. Heat shock cells at 37°C for 1 minute. Do not shake.
 - 7. Place cells on ice for 2 minutes.
 - 8. Add 800 μ l of L.B. medium.
- 9. Shake at 225 rpm for 1 hour at 37°C to express ampicillin resistance.
- 10. In the interim, remove the appropriate number of LB-Amp plates out of cold storage, and spread 20 μl of 50 mg/ml X-GAL and 100 μl of 100 mM IPTG onto each plate.
- 11. After 1 hr. of incubation, pellet cells (except for the tube containing the pUC19 control DNA) in a microcentrifuge by centrifugation for 1 minute. Decant supernatant.
 - 12. Resuspend cell pellets in 100 μ l of L.B.
- 13. Plate 50 μ l of each suspension on a L.B.-Amp plate containing X-GAL and IPTG. Plate 10 μ l of the 1 ml of pUC19 control transformation.

14. Incubate the plates overnight at 37°C.

E. PCR characterization of inserts

- ____ 1. With a sterile toothpick, pick a white colony and lightly touch a plate of LB-Amp with it (do not streak) to transfer some cells onto the plate.
- 2. Immediately dip the same toothpick into an EppendorfTM tube containing 25 μl of 5 mM Tris-HCl, 0.1 mM EDTA (pH=8.0). Shake the toothpick to suspend cells and then discard the toothpick.
- 3. With a fresh toothpick, streak out cells that were patched onto the LB-Amp plate. Incubate the plates overnight at 37° C.
- 4. Boil buffer-suspended cells for 5 min. by placing them in a beaker of boiling water over a Bunsen burner flame.
- 5. Spin down cellular debris for 1 min. in a microcentrifuge.
- 6. Transfer 3 μl of supernatant to wells of a microtiter dish for the PCR reaction.
- 7. For the PCR reaction, prepare a master mix by combining: 20 μ l 10 mM dNTP mix; 2.74 μ l of Universal primer (1 mg/ml); 2.33 μ l of Reverse primer (1.1 mg/ml); 100 μ l of 10X Taq polymerase buffer (Boehringer Mannheim); 10 μ l of Taq DNA polymerase (5 U/ml, Boehringer Mannheim), and dH₂0 to 1 ml (i.e. 864.93 μ l).
- 8. Add 47 μl of the PCR master mix to each well of the microtiter dish containing 3 μl of boiled cell lysate.
 - 9. Mix and cover with mineral oil.
- 10. Place microtiter dish in a PCR machine, and perform the PCR reaction by first incubating at 95°C for 5' and then performing 30 cycles of 94°C for 10"; 54°C for 30"; and 72°C for 30".
- 11. Add loading dye to the PCR reactions, and load onto a 2% agarose gel with the 123 bp marker DNA (2 μ 1).

12. Examine gel to confirm that the bands are of the expected size.

Note that about 220 bp of vector sequence is being amplified in addition to insert sequence, so the size expected for the band on the gel is equal to the DD band size plus 220 bp. Those clones that appear to contain inserts of the correct size can then be sequenced.

F. PCR Sequencing of bacterial colonies

___ For this protocol, one should use fresh colonies (i.e., one or two day old plates).

- 1. As in the procedure for characterization of inserts of transformants, pick colonies with sterile toothpicks, and resuspend the cells in 25 μ l of 5 mM Tris-HCl, 0.1 mM EDTA (pH=8.0). This time, however, cells are picked from plates containing streaked-out colonies from transformation plates, not from the original transformation plates.
- 2. Boil cells for 5 min., place on ice, and spin down cellular debris for 1 minute in a microcentrifuge.
- 3. Transfer 10 μl out of supernatant to a fresh Eppendorf^m tube to sequence.
- 4. To 10 μ l of the supernatant, add the following: 1 μ l of 1 mM Universal primer (1 pmol); 4 μ l of 10X sequencing buffer; 1 μ l (10 mCi) of [alpha-32P]-dATP; 1 μ l (2 U) of Taq polymerase; and dH₂O to 30 μ l (i.e. 13 μ l).
- 5. Add 3 μ l of each of the four ddNTP (G, A, T, C) mixes to the wells of a microtiter dish, according to the number of sequencing reactions being performed.
- 6. Aliquot 7 μ l of the sequencing reaction mixture from step 4 into each of the 4 termination tubes containing 3 μ l of ddNTP. Mix and overlay with one drop of mineral oil.
- 7. Run reactions in the PCR machine under the following conditions: an incubation at 95° C for 5', followed by 30 cycles of 95° C for 30"; 60° C for 30"; and 72° C for 1'.
- 8. When reactions are complete, add 5 μl of stop mix to each well.
- 9. Heat-denature the samples for 5' at 95°C in a PCR machine, and load 5-6 μl on a 6% acrylamide sequencing gel.
- 10. Run the gel until the xylene cyanol dye front is approximately 5 cm from the bottom of the gel, so that

the sequence corresponding to the <u>HindIII</u> site will be near the bottom of the gel.

G. Analysis of DNA sequence from clones of a band

- ____ 1. Locate the sequence of the <u>HindIII</u> site in the clones, thus locating the beginning of the insert sequence. Confirm that vector sequences are present.
- 2. Scan the different clones to see if any match another. Also compare sequences to sequences previously determined directly from acrylamide slices for matches. Pick clones that match previous sequences for mini-prep analysis.

At this stage, there is no need to read and record the sequences in detail, as this can be done when the plasmids from small-scale plasmid preparations ("minipreps") are sequenced. Instead, scan sequences for comparison purposes to determine the clones that should be differentially expressed.

Sequences determined directly from acrylamide slices were obtained using the 5'-Differential Display primers as sequencing primers. Therefore, these sequences are from the 5' end of the bands. The bands cloned in this procedure were inserted into pBluescript directionally. Therefore, the sequences obtained from the Universal primer can start from the 5' or 3' One can distinguish which end is being sequenced by presence or absence of a long stretch of T residues (complementary to the poly A tail of the 3'-end) immediately after the <u>HindIII</u> site. The presence of the stretch of T's indicates the 3'-end of the gene is being read, whereas its absence means the 5'-end of the gene is being read. If the 5' end of the gene is being read, then this sequence can be compared directly to the previously derived sequences obtained from the gels. If, however. the 3'-end of the gene is being read, the sequence has to

be "inverted" to make a comparison with the gel-derived sequences. The clone may have to be sequenced from the "other end" with the Reverse primer to obtain the 5'-end sequence so that a direct comparison can be made.

Where sequence information was not directly obtainable from the gels (i.e., the "Unknown" genetags in Table 2, above), steps 1 and 2, above, apply, except that one does not initially know what sequence corresponds to the differentially displayed band. One should look for which sequences match, and pick the sequence that appears the most frequently for mini-prep analysis. mean picking more than one band to mini-prep and test by Northern analysis to determine which the differentially expressed band.

H. Alkaline lysis miniprep procedure

___ This procedure is adapted from a procedure described in Book 1 of <u>Molecular Cloning</u> by Sambrooke, Fritsch and Maniatis (pp. 1.25-1.28).

- 1. Prepare overnight cultures of cells containing the plasmid of interest in 5 ml of LB-Amp (5 ml of LB-Amp should contain 5 μ l of 100 mg/ml Ampicillin for a final concentration of 100 μ g/ml).
- 2. Transfer 1.5 ml of culture to a fresh Eppendorf™ tube and pellet cells in a microcentrifuge for 1 minute. Transfer another 1.5 ml of the culture to the same tube, and pellet these cells, so that the tube contains the combined pellets from 3 ml of culture.
- 3. Resuspend pellets in 100 μl of Solution 1. Vortex briefly.
- 4. Add 200 μ l of 1% SDS, 0.2 M NaOH. This solution should be freshly made (just before using) from stocks of 20% SDS and 10 M NaOH. Mix gently by inversion to resuspend cells. Leave on ice for 5 minutes.
- 5. Add 150 μ l of Solution 3. Mix gently by inversion and incubate 5' on ice.

6. Centrifuge the solution in a microcentrifuge for 2-4 min.

- 7. Transfer the supernatant to a fresh tube and discard the pellet. Phenol/chloroform extract (400 μ l) the supernatant, shaking the tube by hand for 20" and then centrifuging the tube for 2 min. Transfer the top (aqueous) layer to a new tube.
- 8. Add 1 ml of 100% ethanol to the solution and mix.
- 9. Pellet the DNA by centrifugation in a microcentrifuge for 10'.
- 10. Remove the supernatant. Add 200 μl of 70% Ethanol to wash the pellet, and centrifuge for 2' in a microcentrifuge.
- 11. Remove the supernatant by vacuum suction, making certain that all of the EtOH is removed.
 - 12. Dissolve the pellet in 50 μ l of TE buffer.
- 13. Digest 8 μ l with <u>Hin</u>dIII in the presence of 1 μ l of DNase-free RNase (concentrations of 500 μ g/ml 3 mg/ml are satisfactory).
- 14. Run entire digest on a 2% low melting point agarose gel, and confirm the insert is of the right size (using the 123 bp marker). Cut out the band from the gel, and place the band in an Eppendorf^m tube. This band can then be used to prepare a probe to screen Northern Blots.

Example 8

Making Oligonucleotide Probes from Known Genes

— When one has identified a putative senescence-related gene or genetag, one can readily verify that the gene or genetag is senescence-related by Northern analyses or in situ hybridization, as discussed above. For both of these procedures, however, one requires a probe that is relatively specific, for if the probe hybridizes to a

sequence that is rather abundant in the RNA population of a cell, then the results of any procedure involving probe hybridization to an RNA in that population could be ambiguous. Consequently, there are advantages to first checking the specificity a probe sequence prior to performing a Northern blot or in situ hybridization experiment.

Where the putative senescence-related gene is a known gene, one can readily prepare a variety of synthetic oligonucleotide probes from the known sequence. Typically, such probes are 20 to 60 nucleotides in length, with longer probes preferred for specificity. The specificity of such probes can be conveniently analyzed by Southern hybridization against genomic DNA. If the probe only hybridizes to a few bands (less than 5, and preferably less than 3 bands), then the probe is specific enough for use in verifying that a gene is senescence-related and for use in screens to determine whether a compound affects expression levels senescence-related genes in senescent cells.

To prepare an oligonucleotide probe from a known gene, one can conveniently label the probe with a radioactive label by a kinase reaction procedure. To 20 µl kinase reaction, prepare a use 600 oligonucleotide (6 μ l of a 40 mer oligonucleotide solution 100 $ng/\mu l)$, at and add the following 7 μl of water; 2 μl of 10X One-phor-all™ ingredients: buffer (Pharmacia); 2 μ l of Pharmacia T4 polynucleotide kinase; and 5 μ l of $gamma-^{32}P-ATP \qquad (3000)$ Incubate the reaction mixture at 37°C for 30', and purify the labelled probe over a Pharmacia S-200 Sephacryl™ spin column.

The following Example describes how to prepare probes from genetags of novel genes, and Examples 10 and 11 describe how to perform the Southern analysis to check for the specificity of probe hybridization.

Example 9

Making Probes from EDD-Cloned Bands

— When one identifies a genetag of a novel gene as a putative senescence-related gene, then only the genetag sequence is initially available for use in probe design. To perform a rapid analysis of whether the genetag sequence can be used as a probe with a high degree of specificity, one can conveniently prepare probes from the genetag clone (see Example 7) by digesting the plasmid with HindIII (for plasmids prepared as per Example 7), separating the resulting fragments on a low-melting agarose gel (as per the procedure in Example 7), and separating and removing the band from the gel as follows.

- 1. Melt 10-15 μ l of the gel slice at 65°C for 2'.
- 2. Add 10 μl of 20 μM 5'-primer and 10 μl of 20 μm 3'-primer to the gel slice solution.
- 3. Heat this solution at 100°C for 10° and then quick-cool on ice.
- 4. Then, add: 2 μ l of BSA (10 mg/ml); 1.5 μ l of Klenow enzyme; and 10 μ l of a 5X oligonucleotide-labelling buffer containing neither primers nor dCTP. The 5X Labeling Buffer contains 250 mM Tris (pH=8); 25 mM MgCl₂; 5 mM beta-mercaptoethanol; 2 mM dATP, dTTP, dGTP; 1 M HEPES (pH=6.6); and 5 μ l of alpha-32P-dCTP (3000 Ci/mmol).
 - 5. Incubate this solution at 37°C for 30'.
- 6. Heat this solution briefly to melt the gel slice if the gel slice has solidified, and then purify the labelled probe over a Pharmacia S-200 spin column.

Example 10

Genomic Digests, Gel Electrophoresis, and Transfer for Southern Analysis

___ To check the specificity of a probe for use in Northern analyses or in situ hybridization, one can hybridize the probe to restriction enzyme-digested genomic DNA. To prepare the genomic DNA, one first isolates (or purchases) genomic human DNA by any of a variety of standard methods and then digests the DNA with restriction enzymes. For best results, the restriction enzyme digestion is conducted with a variety of different restriction enzymes that have 6-base recognition sequences, i.e., HindIII, EcoRI, and BamHI.

Restriction digest conditions are 500 μ l of total reaction volume containing 10 μ g of human genomic DNA, 50 μ l of 10X restriction enzyme buffer, and 5 μ l of restriction enzyme, and the reactions are incubated at 37°C for 30 min., at which point, one can add an additional 5 μ l of enzyme and let the reaction continue for 6 hrs. The digested DNA is precipitated by adding 50 μ l of 3 M NaAcetate and 1 ml of cold 100% ethanol and centrifuging for 10 min. at room temperature at 14K. Aspirate the ethanol, and let the DNA pellet air dry for 20 min., or dry the DNA in a speed-vac centrifuge for 10 min. Resuspend the DNA in 20 μ l of water by incubation at 4°C overnight (pellet may require more time to resuspend).

For an 0.8% agarose gel, add 2.4 g of agarose to 300 ml of 0.5% TBE buffer. Add 3 μ l of loading dye to each sample, and load the sample into the 0.8% agarose gel. Run the gel at 35 V for 6 hrs. The lower voltage and longer running time tighten up the restriction fragment bands. Take a photograph of the gel after staining with ethidium bromide.

To blot the gel, soak the gel in 0.6 M NaCl and 0.4 M NaOH for 10 minutes at room temperature. The transfer solution is the same as the soak solution. Place the gel upside down on the filter paper wick. Cut a piece of Schleicher and Shuell membrane to size, and place the rough side down onto the gel. Place two pieces of filter paper on top of this membrane, then place a stack of terriwipes on top with a weight pressing the terriwipes down. Cover the buffer chambers with saran wrap and allow the transfer to continue overnight.

After the transfer, soak the membrane in 5X SSC buffer for five minutes; then, blot excess liquid off. Place the membrane in the Stratalinker^m light box (Stratagene) with the DNA side up, and cross link at 1200 kJ (autocrosslink mode). Check the flattened gel to make sure there is no DNA left in the gel.

Example 11

Southern Hybridization and Wash Conditions

A. Hybridization

___ Prehybridize the blot for 2 hours at 60°C using the following prehybridization solution: 5X SSC (from 20X SSC purchased from Boehringer Mannheim); 1X Denhardt's (50X contains 5 g Ficoll (Type 400, Pharmacia); 5 .g polyvinylpyrrolidone; 5 g BSA (Fraction V, Sigma); and water 500 ml; this stock to solution should distributed into 4 ml Corning tubes and stored at -20°C); 0.1% SDS; 0.05% sodium pyrophosphate (diluted from a 5% stock solution and 150 μ g/ml denatured salmon sperm. After the 2 hour prehybridization incubation, remove all of the solution, and replace with hybridization solution, which is the same as prehybridization solution with the addition of 10% dextran sulfate. Boil the probe for several minutes, quench briefly on ice, and then add this solution to the blot soaking in the hybridization

solution. Double bag the blot, and allow the hybridization to continue at 60°C overnight.

B. Wash Conditions

After the hybridization, the probe can be collected and saved in a 15 ml Corning tube at room temperature. Count 5 μ l of this probe solution in the scintillation counter. Rinse the blot twice at room temperature in the following solution: 3X SSC, 0.1% SDS, and .05% sodium pyrophosphate. Conduct four 15' washes at 60°C with this solution, which should be pre-heated to 60°C.

C. Analysis

After the Southern blot is probed and washed, the blot is analyzed to determine where the probe is bound on the membrane by PhosphoImager[™] analysis (Molecular Dynamics) or autoradiography. Probes with suitable specificity for use in Northern analyses and in situ hybridizations are identified by a hybridization pattern in which the probe has hybridized to only a few (less than 5 and preferably less than 3) bands.

Example 12

Denaturing RNA Agarose Gels for Northern Analysis

___ When making a formaldehyde gel, always use baked glassware, protective gloves, and DEPC-water. To prepare the gel, add 2 g of agarose to 140 ml of DEPC-H,O (which in turn is prepared by adding 1 ml of DEPC to 1000 ml of deionized water, resulting in a 0.1% solution of DEPC, heat the resulting solution with stirring in a hood for several hours, and autoclave); boil the resulting solution in a 500 ml flask for approximately three minutes, and then place the solution in a 60°C water bath De-ionize 65 ml of 37% formaldehyde by to equilibrate. adding several grams of mix bed resin (from Bio Rad) to

the formaldehyde in a baked ehrlenmyer flask in a hood, mixing with a stir bar for approximately five minutes, and then filtering away the resin by pouring this mixture through a funnel coated with 3MM circular WhatmannTM paper, after which the formaldehyde is ready to be used.

Next, add 44 ml of 5X running buffer to the agarose/water solution being incubated at 60°C. 5X running buffer is composed of 0.1 M MOPS (pH=7), which is prepared from a 0.5 M stock solution in which the pH has been adjusted using either NaOH or Acetic Acid and stored at room temperature in a bottle wrapped in aluminum foil; 40 mM sodium acetate; and 5 mM EDTA. Then, add 40 ml of the 37% deionized formaldehyde, mix the solution well, and pour the gel in a hood. Pre-run the gel for 5' at 150 volts in the hood; always wear gloves when handling the gel box.

Samples are prepared as follows. One typically loads 20 μ g (4.5 μ l) of RNA per lane, which is mixed prior to loading with 2 μ l of 5X running buffer; 3.5 μ l of deionized formaldehyde; and 10 µl of formamide. This mixture is heated for 15' at 65°C and then centrifuged briefly to collect the mixture at the bottom of the tube. Add 2 µl of a 10X loading buffer (which consists of: glycerol; 1 mM EDTA, pH=8; 0.25% bromophenol blue; 0.25% xylene cyanol FF) to each sample. Prepare the lambda-HindIII standards in the same way as above (denatured) using 1 µg of DNA per lane. Load the samples on the gel, and run the gel at 80 volts for 3 to 4 hours in the hood. Stop the gel half-way through the run to shake the gel lightly, circulate buffer, and add new buffer. After the gel has finished, wash the gel 3X in DEPC-H,O with each lasting 10 minutes to remove most formaldehyde. Rinse the gel in the gel-casting tray.

Cut off the lambda-<u>HindIII</u> lane with an RNA lane next to it to stain in ethidium bromide for visualization. To perform the ethidium stain, place the gel slice in a gel-

casting tray filled with DEPC- H_2O ; add a few drops of 5 mg/ml ethidium bromide to the tray; stain for 20 minutes at room temperature; perform two 10 minute washes in DEPC- H_2O ; leave overnight in water in the refrigerator to destain completely; and place the gel next to a ruler on the light box to visualize. Take a photo with the aperture set at 8 and the time set at 0.25 seconds. Take the rest of the gel and blot the gel onto a positively charged nytran membrane, 0.45 μ m pore size, "Nytran + Maximum Strength" (Schleicher and Schuell).

Next, one sets up the transfer apparatus; there is no need to pre-wet the membrane. Place the membrane onto the gel; then, place two pieces of 3MM whatman paper on top of the membrane, place a stack of terri towels on top of paper, and place a heavy weight on top of the towels. The transfer buffer is 20X SSC; the transfer should be allowed to continue overnight. Label the membrane with a pencil or ball point pen; wash the membrane in 5X SSC for several minutes; and lightly blot the membrane with filter paper. To link the RNA to the membrane, use the Stratalinker (autolink mode, 120 mJ/cm²) light box. Place the membrane on a piece of filter paper with the RNA side up. under the light.

Example 13

Northern Hybridization and Wash Conditions

Prehybridize the blot for at least 30' at 42°C using hybridization solution composed of 50% formamide; 5% Denhardt's; 0.5% SDS; 5% SSPE or SSC; and 100 μ g/ml of salmon sperm DNA. Boil the salmon sperm DNA to denature the DNA before adding it to the hybridization solution. If one is reusing hybridization solution, boil the solution for 2' and quick cool on ice before adding the solution to the blot. After boiling the probe (DD bands), one adds the probe to the blot, double-bags the

blot plus hybridization solution, and incubates the blot at 42°C overnight.

After hybridization, the probe is collected (and can be reused; store at room temperature), and the blot is rinsed twice at room temperature in the a solution composed of 1X SSC and 0.5% SDS. Then, two 30' washes at 65°C are performed using the same (but fresh) wash solution to wash away non-specifically bound probe, leaving only the specifically bound probe.

As noted above, the Northern analysis procedure can used to confirm that a genetag or probe specifically identify a senescence-related Alternatively, one can use RT-PCR (Reverse Transcriptasemediated Polymerase Chain Reaction) for this purpose. After a gene or genetag has been confirmed to be a senescence-related gene for cells in tissue culture, one can confirm that the gene is senescence-related in vivo by an in situ hybridization procedure, such as procedure described in the Novagen Suresite™ II System manual, supra, using tissues from young and old donors.

Once a probe has been confirmed as identifying specifically a senescence-related gene, the Northern analysis procedure or RT-PCR can be used to identify whether a compound can reverse, partially reverse, or modulate the pattern of expression of senescence-related genes in senescent cells. Once such a compound is identified, one can determine whether the compound has activity in vivo by analyzing tissues from treated animals with an in situ hybridization procedure. The in situ procedure can also be used to identify senescent or young cells in tissues using senescent-related gene probes of the invention that have been confirmed to identify senescence-related genes in vivo.

Example 14

B-Galactosidase Screen

A. Primary Screen.

Senescent cells are seeded in 96-well plates at 10,000 to 20,000 cells/well in DMEM medium plus 10% Bovine Calf Serum (BCS). In two preferred embodiments, senescent human embryonic lung fibroblasts (IMR90 cells) are used at Passage Doubling Level (PDL) 53, or senescent fibroblast lines derived from human foreskin (BJ cells) are used at PDL 92. Other senescent cells, in the appropriate media, can also be used. After 6 hours, the medium is removed and replaced with DMEM plus 0.5% BCS. After 3 days, the medium is replaced with fresh medium, and the sample or its vehicle is added. In a preferred embodiment, 2 μ l of sample dissolved in DMSO (1 μ M final concentration), or of DMSO alone, are added to 200 μl of medium. Other volumes. vehicles. and compound concentrations can also be used. In addition, mixtures of compounds, rather than single compounds, can be added After 4 days (or other appropriate to the cells. incubation time), the medium is again removed and the cells are fixed and stained.

To fix the cells, the medium is removed and replaced with phosphate-buffered saline (PBS). This, and all other liquid transfers can be accomplished using a Hamilton Microlab 2000™ pipeting station. Other pipeting stations, or manual pipetting, can also be used. is then removed and again replaced with PBS. The PBS is removed and replaced with freshly-prepared fixing solution (0.5% gluteraldehyde in PBS). The cells are incubated in this mixture for 2 min. at room temperature. The fixing solution is removed and replaced with PBS. The PBS is replaced with fresh PBS, and the cells are incubated for an additional 10 min. at room temperature. Other methods of fixing the cells can also be employed.

To stain the cells, 100 μl of X-gal (5-bromo-4-chloro-3-indolyl- β -D-galactosidase, at a concentration of

50 mg/ml in dimethylformamide Promega, Madison, WI), is to 10 ml of staining buffer (40 mM citric acid/Na,HPO, buffer pH=6, 5 mM potassium ferrocyanide, 5 mM potassium ferricyanide, 150 mM NaCl, and 2 mM MgCl, in distilled water). The X-gal can be replaced with other substrates for ß-galactosidase in the appropriate buffer. The PBS is removed and replaced with a sufficient volume of staining solution (50 μ l are routinely used) to cover the cells. The microtiter plate is then covered, sealed in a humidified container, and placed in an incubator at 37°C overnight. Movement of the plates between the pipetting station and the incubator can be done either manually or with a robot (e.g. Zymark XP (TM)).

Following the overnight incubation, the intensity of staining is measured by, e.g., using a plate reader at a wavelength of 540 nM. Quantitation can also be performed microscopically or with the aid of an image-analysis system. Decreased intensity of the staining in the presence of the test compound indicates reversal of the senescent phenotype. Samples that produce this affect are tested in the secondary assay. The plates can be sealed and stored at 4°C indefinitely.

B. Secondary Screen

Samples that decrease staining in senescent cells are then tested for ability to decrease staining in young cells. In two preferred embodiments, young IMR90 cells are used at a PDL lower than 35, or young BJ cells are used at a PDL lower than 55. Other young cells, in the appropriate media, can also be used. The secondary screen is carried out in the same manner as the screen in senescent cells with two modifications. First, young cells, rather than senescent cells, are used. Second, the staining buffer is adjusted to pH 4 rather than pH 6. A decrease in ß-galactosidase staining in senescent but

not young cells is interpreted as a reversal of the senescent phenotype.

Other high-throughput screens for compounds alter the expression of specific senescence-related genes include screens for fibronectin, collagen 1 (alpha 1 and 3), and elastin (see Ahmed et al., 1992, A colorimetric for glycated microassay collagen based the thiobarbituric acid method, Clinica Chimica Acta. <u> 212</u>: 133-139; Anderson and Elliot, 1991, A dye-binding assay for soluble elastin, Biochem. Soc. Trans. 19:388S; Clark al.. 1992. Monoclonal antibodies against human fibroblast collagenase and the design of an enzyme-linked immunosorbent assay to measure total collagenase, Matrix 12: 475-480; Walsh et al., 1992, Microplate reader-based quantitation of collagens, Analyt. Biochem. 203: 187-190; and Scutt et al., 1992, A semiautomated, 96-well plate assay for collagen synthesis, Analyt. Biochem. 203: 290-294). As with the beta-galactosidase screen described in this Example, one first determines whether a compound can modulate the activity or expression level (protein levels can be determined, for example, by gel analysis or by antibody-based methods) of a senescence-related gene and then determines whether the compound has the modulatory effect on a panel of senescence related genes, preferably using Northern analysis, RT-PCR or in situ hybridization with probes from known senescence-related genes, known genes that have been determined to be senescencerelated by the method of the present invention, and from senescent gene-related genetags from previously unknown genes provided by the present invention.

Example 15

Senescence-related Gene Expression Screen

___ To determine the effect of a compound on the mRNA levels of known and novel senescence-related genes

the method of the present according to invention. senescent cells are grown in a 10 or 15 cm plate using the same protocol as for the high throughput screens. One plate is incubated with the test compound, another plate is incubated with the compound test vehicle After four to 20 days of incubation, the cells are lysed in GITC (see the RNA isolation protocol. above), and RNA is prepared. The RNA is analyzed with the senescence-related gene probes of the invention by Northern analysis or by other suitable methods, such as RT-PCR. The results of this analysis will indicate the efficacy of the compound in altering the mRNA expression level on senescence-related genes. The expression levels of at least two, and preferably 3 to 5 to 10 to 20 or more, senescence-related genes will be determined.

Thus, compounds are tested to determine whether the compounds alter the expression of the young- and old-specific senescence-related genes identified by EDD and in the scientific literature. If a compound has the effect of complete reversal to a young pattern of gene expression, then the compound impacts a single common mechanism driving cell senescence. Reversal of defined groups of genes indicates that several mechanisms are operating in senescence and that different mechanisms can affect different panels of genes. A compound may also act to modulate the activity of an individual gene, suggesting the absence of a common mechanism.

An alternative screen for compounds that alter the expression of senescence-related gene involves the use of a genetic construct comprising a promoter of a senescence-related gene positioned for expression of a coding sequence from a reporter gene, such as an alkaline phosphatase gene, the expression of which can be efficiently and readily monitored. Such a construct would be used to generate stable transfectants in very young cells, such as dermal fibroblasts, and then the

cells could be used at any stage up to and including senescence to identify agents that up or down-regulate the expression of the reporter gene.

The foregoing examples describe various aspects of the invention and how the methods of the invention can be The examples are not intended to provide an exhaustive description of the many different embodiments invention. All publications and applications cited above are hereby incorporated herein by reference as if each individual publication or patent application were specifically and individually indicated to be incorporated by reference. Thus, although the foregoing invention has been described in some detail by way of illustration and example for purposes of clarity of understanding, it will be readily apparent to those of ordinary skill in the art in light of the teachings of this invention that certain changes and modifications may be made thereto without departing from the spirit or scope of the appended claims.

PCT/US95/11230 WO 96/13610

SEQUENCE LISTING

(1) GENERAL INFORMATION:

(i) APPLICANT:

Maarten H.K. Linskens, et al.

(ii) TITLE OF INVENTION:

METHODS AND REAGENTS FOR THE IDENTIFICATION AND REGULATION OF

SENESCENCE-RELATED GENES

(iii) NUMBER OF SEQUENCES:

(iv) CORRESPONDENCE ADDRESS:

(A) ADDRESSEE:

Lyon & Lyon

(B) STREET:

633 West Fifth Street

(C) CITY:

Los Angeles

(D) STATE:

California USA

(E) COUNTRY: (F) ZIP:

90017

130

(v) COMPUTER READABLE FORM:

(A) MEDIUM TYPE:

(B) COMPUTER: 3.5" Diskette, 1.44 Mb storage

(C) OPERATING SYSTEM: IBM P.C. DOS (Version 5.0)

(D) SOFTWARE: WordPerfect (Version 5.1)

(vi) CURRENT APPLICATION DATA:

(A) APPLICATION NUMBER: 08/332,420

(B) FILING DATE:

October 31, 1994

(C) CLASSIFICATION:

(vii) PRIOR APPLICATION DATA:

(A) APPLICATION NUMBER:

08/235,180

(B) FILING DATE:

April 29, 1994

(viii) ATTORNEY/AGENT INFORMATION:

(A) NAME:

Warburg, Richard J.

(B) REGISTRATION NUMBER:

32,327

(C) REFERENCE/DOCKET NUMBER: 210/041

(ix) TELECOMMUNICATION INFORMATION:

(A) TELEPHONE:

(213) 489-1600 (213) 955-0440

(B) TELEFAX: (C) TELEX:

67-3510

(2) INFORMATION FOR SEQ ID NO:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH:

22

(B) TYPE:

nucleic acid

(C) STRANDEDNESS:

single

(D) TOPOLOGY:

linear

(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 1:

93

WO 96/13610	PCT/US95/11230
GCGCAAGCTT TTTTTTTTT CT	22
(2) INFORMATION FOR SEQ ID NO: 2:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 2:	
GCGCAAGCTT TTTTTTTTT CC	22
(2) INFORMATION FOR SEQ ID NO: 3:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 3:	
GCGCAAGCTT TTTTTTTTT CG	22
(2) INFORMATION FOR SEQ ID NO: 4:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 4:	
GCGCAAGCTT TTTTTTTTT GT	22
(2) INFORMATION FOR SEQ ID NO: 5:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: :	
GCGCAAGCTT TTTTTTTTT GG	22
(2) INFORMATION FOR SEQ ID NO: 6:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	

PCT/US95/11230 WO 96/13610 (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 6: 22 GCGCAAGCTT TTTTTTTTT GA (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 22 (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 7: GCGCAAGCTT TTTTTTTTT AT 22 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 8: 22 GCGCAAGCTT TTTTTTTTT AC (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 22 (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 9: GCGCAAGCTT TTTTTTTTT AG 22 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 10: GCGCAAGCTT TTTTTTTTT AA 22 (2) INFORMATION FOR SEQ ID NO: (i) SEOUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: 22 nucleic acid (C) STRANDEDNESS: single

linear

(D) TOPOLOGY:

(ii) SEQUENCE DESCRIPTION: SEQ	ID NO: 11:
GCGCAAGCTT TTTTTTTTT CA	22
(2) INFORMATION FOR SEQ ID NO: 12:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	22 nucleic acid single linear
(ii) SEQUENCE DESCRIPTION: SEQ I	D NO: 12:
GCGCAAGCTT TTTTTTTTT GC	22
(2) INFORMATION FOR SEQ ID NO: 13:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	22 nucleic acid single linear
(ii) SEQUENCE DESCRIPTION: SEQ I	D NO: 13:
CGGGAAGCTT ATCGACTCCA AG	22
(2) INFORMATION FOR SEQ ID NO: 14:	
(i) SEQUENCE CHARACTERISTICS:	•
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	22 nucleic acid single linear
(ii) SEQUENCE DESCRIPTION: SEQ I	D NO: 14:
CGGGAAGCTT TAGCTAGCAT GG	22
(2) INFORMATION FOR SEQ ID NO: 15:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	22 nucleic acid single linear
(ii) SEQUENCE DESCRIPTION: SEQ II	O NO: 15:
CGGGAAGCTT GCTAAGACTA GC	22
(2) INFORMATION FOR SEQ ID NO: 16:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH:(B) TYPE:(C) STRANDEDNESS:(D) TOPOLOGY:	22 nucleic acid single linear

PCT/US95/11230 WO 96/13610

(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 16:	
CGGGAAGCTT TGCAGTGTGT GA	22
(2) INFORMATION FOR SEQ ID NO: 17:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 17:	
CGGGAAGCTT GTGACCATTG CA	22
(2) INFORMATION FOR SEQ ID NO: 18:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 18:	
CGGGAAGCTT GTCTGCTAGG TA	22
(2) INFORMATION FOR SEQ ID NO: 19:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 19:	
CGGGAAGCTT GCATGGTAGT CT	22
(2) INFORMATION FOR SEQ ID NO: 20:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 20:	
CGGGAAGCTT GTGTTGCACC AT	22
(2) INFORMATION FOR SEQ ID NO: 21:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 22 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	

(ii) SEQUENCE	E DESCRIPTION: SE	Q ID NO: 21:	
CGGGAAGCTT AGACGCTAC	FT GT		22
(2) INFORMATION FOR	SEQ ID NO: 22:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	22 nucleic acid single linear	1
(ii) SEQUENCE	DESCRIPTION: SE	Q ID NO: 22:	
CGGGAAGCTT TAGCTAGCA	G AC		22
(2) INFORMATION FOR	SEQ ID NO: 23:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	22 nucleic acid single linear	ı
(ii) SEQUENCE	DESCRIPTION: SEC	Q ID NO: 23:	
CGGGAAGCTT CATGATGCT	A CC		22
(2) INFORMATION FOR	SEQ ID NO: 24:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	22 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEG	Q ID NO: 24:	
CGGGAAGCTT ACTCCATGA	C TC	·	22
(2) INFORMATION FOR	SEQ ID NO: 25:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	22 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEC	Q ID NO: 25:	
CGGGAAGCTT ATTACAACG	A GG		22
(2) INFORMATION FOR	SEQ ID NO: 26:		
(i) SEQUENCE	CHARACTERISTICS:		
(B)	LENGTH: TYPE: STRANDEDNESS:	22 nucleic acid single	

(D) TOE	POLOGY:		linear		
(ii) SEQUENCE DES	CRIPTION: S	EQ ID	NO: 2	6:	
CGGGAAGCTT ATTGGATTGG TO					22
(2) INFORMATION FOR SEQ	ID NO: 27	:			
(i) SEQUENCE CHA	RACTERISTICS	:			
			22 nuclei single linear	c acid	
(ii) SEQUENCE DES	CRIPTION: SI	EQ ID	NO: 2	7:	
CGGGAAGCTT ATCTTTCTAC CC					22
(2) INFORMATION FOR SEQ	ID NO: 28	:			
(i) SEQUENCE CHA	RACTERISTICS	:			
			22 nucleic single linear	c acid	
. (ii) SEQUENCE DES	CRIPTION: SE	EQ ID	NO: 28	B:	
CGGGAAGCTT ATTTTTGGCT CC				:	22
(2) INFORMATION FOR SEQ	ID NO: 29:	:			
(i) SEQUENCE CHA	RACTERISTICS:	:			
(A) LEN (B) TYP (C) STR (D) TOP			22 nucleid single linear	c acid	
(ii) SEQUENCE DES	CRIPTION: SE	EQ ID :	NO: 29	9:	
CGGGAAGCTT TATCGATACA GG				:	22
(2) INFORMATION FOR SEQ	ID NO: 30:	:			
(i) SEQUENCE CHA	RACTERISTICS:	:			
(A) LEN (B) TYP (C) STR (D) TOP	E: ANDEDNESS:	1	22 nucleic single linear	c acid	
(ii) SEQUENCE DES	CRIPTION: SE	EQ ID	NO: 30):	
CGGGAAGCTT TATGGTAAAG GG				2	22
(2) INFORMATION FOR SEQ	ID NO: 31:	:			
(i) SEQUENCE CHA	RACTERISTICS:	:			
(A) LEN (B) TYP			22 nucleio	c acid	

	STRANDEDNESS: TOPOLOGY:	single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ II	NO: 31:	
CGGGAAGCTT TATCGGTCA	T AG		22
(2) INFORMATION FOR	SEQ ID NO: 32:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	22 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 32:	
CGGGAAGCTT TAGGTACTA	A GG		22
(2) INFORMATION FOR S	SEQ ID NO: 33:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	TYPE:	27 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 33:	
CATTTATTCA TTCATTGAGA	A CACTCAA		27
(2) INFORMATION FOR S	SEQ ID NO: 34:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)		112 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 34:	
ACAGAAAGGC CACTCAGGAT CAATAATCCT TAAACTAGCT		TG TCATTCAGCA GTCAGTCCCC	C 60 112
(2) INFORMATION FOR S	SEQ ID NO: 35:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	110 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 1:	
TTGAGTAGTT ACTGGAACCT		AG GTACTTCCAA AAAAAGGAC TA CATAATATA	C 60 110
(2) INFORMATION FOR S	SEQ ID NO: 36:		

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 36:	
CAAAGATAAG AAACCAAGGA AGAAAGCAA		30
(2) INFORMATION FOR SEQ ID NO: 37:		
(i) SEQUENCE CHARACTERISTICS:		
	30 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 37:	
CTGACGCCAN CCGCATACNC CGCANCCACA		30
(2) INFORMATION FOR SEQ ID NO: 38:		
(i) SEQUENCE CHARACTERISTICS:		
(A) LENGTH:(B) TYPE:(C) STRANDEDNESS:(D) TOPOLOGY:	60 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 38:	
AGATAAAGCA ATTAGAAGAT GCATTAAAAG ATGTGCAG	AA GAGGATGTAT GAGTCAGAAG	60
(2) INFORMATION FOR SEQ ID NO: 39:		
(i) SEQUENCE CHARACTERISTICS:		
(B) TYPE: (C) STRANDEDNESS:	86 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 39:	
ATAATAAAAC TCTTCATTTT GCGAATTATA GAAGCTAC TAGGGAAACT AAGGAGTGAC ATAGAA	TT TTTATAAAGC CATATTTTTT	60 86
(2) INFORMATION FOR SEQ ID NO: 40:		
(i) SEQUENCE CHARACTERISTICS:		
(A) LENGTH:(B) TYPE:(C) STRANDEDNESS:(D) TOPOLOGY:	43 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 40:	
AACTGCATTT TGATGTTATC GCTTATGTTT AATAGTTA	AT TCC	43
(2) INFORMATION FOR SEC ID NO. 41.		

. . .

(i) SEQUE	NCE CHARACTERISTICS:		
	(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	29 nucleic acid single linear	
(ii) SEQUE	NCE DESCRIPTION: SEC	Q ID NO: 41:	
CTATTGCCTC TCCTCC	IGCA GAGACCATG		29
(2) INFORMATION F	OR SEQ ID NO: 42:		
(i) SEQUE	NCE CHARACTERISTICS:		
	(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	40 nucleic acid single linear	
(ii) SEQUE	NCE DESCRIPTION: SEC	Q ID NO: 42:	
GAGAAGAAGA	AAGG NCACAGAGAT GGAAG	GCCA	40
(2) INFORMATION FO	OR SEQ ID NO: 43:		
(i) SEQUE	NCE CHARACTERISTICS:		
	(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	60 nucleic acid single linear	
(ii) SEQUEN	ICE DESCRIPTION: SEQ) ID NO: 43:	
GTTTCTGAAT TACATGA	ATT GTTGCAGAGC AAAGA	AACTT ATGGAAATCT TTCCATTTAT	60
(2) INFORMATION FO	OR SEQ ID NO: 44:		
(i) SEQUEN	CE CHARACTERISTICS:		
(A) LENGTH: B) TYPE: C) STRANDEDNESS: D) TOPOLOGY:	27 nucleic acid single linear	
(ii) SEQUEN	CE DESCRIPTION: SEQ) ID NO: 44:	
GTAGGCTTCT ATATTGC	ATT TAACTTG		27
(2) INFORMATION FO	R SEQ ID NO: 45:		
(i) SEQUEN	CE CHARACTERISTICS:		
(A) LENGTH: B) TYPE: C) STRANDEDNESS: D) TOPOLOGY:	72 nucleic acid single linear	
(ii) SEQUEN	CE DESCRIPTION: SEQ	ID NO: 45:	
AATGAGGTAG AAGTAGA ATAATGTGTC CC	AAG GAAGAAAAC TCAAA	GAATT CTAAAAGGAT TCATAGCAAC	60 72
		102	12

(2) INFORMATION FOR SEQ ID NO: 46: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 44 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 46: TCTCACATTC AGTCATACCC TAATGATCCC AGAAAGATAA TCAT 44 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 60 (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 47: AGAAGCCCCA GCAAGATTTA TTCCTTTTTG CTTCTTCATA ACCATGAAGC CATTGAAC 60 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 27 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 48: CTACCTCCCA CATTAATTTT CATATGT 27 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 31 (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 49: AGGGCACAGC ACCAGATGAA TTGTTGTATA T 31 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: nucleic acid (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 50: AAATTAGCTT TCATCACAGA TTTAGGAAACT TGTCT 35

(2) INFORMATION FOR SEQ ID NO: 51: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 47 (B) TYPE: nucleic acid (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 51: AAACTACTGA ACNGTTACCT AGGTTAACAAC CCTGGTTGAG TATTTGC 47 (2) INFORMATION FOR SEQ ID NO: 52: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 52: TTGNATATTG NATTTGTAGT AATATTCCAAA AGAATGTAAA TAGG 44 (2) INFORMATION FOR SEQ ID NO: 53: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 53: AAATTGTATA TTGTATTTGT AGTAATATTCC AAAAGAATGT 40 (2) INFORMATION FOR SEO ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 33 (B) TYPE: nucleic acid (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 54: TATGAATNTC ACATTTGAAT TCTTCGATCTC TAA 33 (2) INFORMATION FOR SEO ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: nucleic acid 37 (C) STRANDEDNESS: (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 55: TATGTATAAA AGCATATGTG CTACTCATCTT TGCTCAC 37

(2) INFORMATION FOR SEQ ID NO: 56: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 33 (B) TYPE: nucleic acid (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 56: AATGTCTAAT TTTCTTTCCG ACACATTTACC AAA 33 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 57: ACAACAGCAA ACAAAAAGGT GAAGTCTAAAT GAAGTGCACA 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 29 (B) TYPE: nucleic acid (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 58: AAAAGAATTG GCAGTTACAT TCATACTTT 29 (2) INFORMATION FOR SEO ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH:
(B) TYPE: 38 nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 59: AAGAATGTGC ATTCCAGTGC CATAGATAGT ATATTGAA 38 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 37 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 60: TTGCTACGGA CTTACGAAAG GACAAAGCGA AGAGCTG 37

105

(2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: 51 nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 61: AAATAATTTA TTCATTGCAG ATACTTTTTA GGTTGCATTT TATTCATTTC C 51 (2) INFORMATION FOR SEQ ID NO: 62: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 39 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 62: AGATGATGAT GTTAACCCAT TCCAGTACAG TATTCTTTT 39 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: 33 nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 63: AGTATAGTGA ATGANTATGC CTTCCTACTG CAG 33 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 33 (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY: nucleic acid single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 64: AGAAATATAA AGATTTTNAT ACCTGCCACA TGG 33 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 34 nucleic acid (B) TYPE: (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 65: GAAGANATTA TGTTGTGANC NGGAGTNACA CAAA 34

(2) INFORMATION FOR SEQ ID NO: 66 . (i) SEOUENCE CHARACTERISTICS: (A) LENGTH: 39 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEO ID NO: 66: AGGGGCACAA GAGTTTGCGG TTATTGAATC CTGAGANAA 39 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 28 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 67: GTTGAAGAGA CAGAGACAAG TAATTTGC 28 (2) INFORMATION FOR SEQ ID NO: 68: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 28 (B) TYPE: nucleic acid (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 68: CCGTGAATAC CCNTTTCTCG ACCAAAGA 28 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 31 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 69: ATGGAGTTGT GGATGAAAGC CATGTTAGNT G 31 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: 37 nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 70: GATCATATAA ACANNNCCGA GTTCTACCTC AGAGTCG 37

(2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 71: CCGTTTACCT TCACGTGGCC ATTCTCCTGT CCGTTCGCTT 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: nucleic acid (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 72: AGGAACAATA TTCACTACTC CAGGAGGAAA CCCTGCCTCT 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 73: CCGAGAGATG CTGTAGCGAC CATTTTTCTC CACGTGGTAT 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 74: CCAGGTATTT CTGGACTAAG TCCACATCTT GCTCTTGTGT 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 75: ATCAGCACCT TTGGGACCAG CATCACCTCT GTCACCCTTA 40

(2) INFORMATION FOR SEQ ID NO: 76: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 nucleic acid (B) TYPE: (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 76: AAGGTTACTG CAAGCAGCAA CAAAGTCCGC GTATCCACAA 40 (2) INFORMATION FOR SEO ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 nucleic acid (B) TYPE: (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 77: CGAGAAGTAG CCAGCTCCCC TTTTGCACAA AGCTCATCAT 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: nucleic acid (C) STRANDEDNESS: (D) TOPOLOGY: single linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 78: CCTGGATAAA AGACTCCTCC AGGAACTCCA CCAGGAATGG 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH:
(B) TYPE: 40 nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 79: AAGAAAGGAT CCTCCTCCTC CACCAGCGCC CCTGTGCTGT 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 80: GAGGAGTATT CAGAACTGGT TTCACACCGA AGGACTAGTT 40

(2) INFORMATION FOR SEQ ID NO: 81: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: nucleic acid single (C) STRANDEDNESS: (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 81: CTCGTCGGTG AGGTCACACT CTGTGTCTGT TGTGTAAAAC 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 82: CTGGTGCTCC GGTCTCGAAT TTTGGCGAAG TGCTTCTGCA 40 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 83: CTCCATATTG ATAGGCGTGC TCTATTGCTC TAGGGCTGTT 40 (2) INFORMATION FOR SEO ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH:
(B) TYPE: 40 nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 84: TCTTGAATCC CATAGCTGCT TGAATCTGCT GCTGGGTTTC (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 40 nucleic acid (B) TYPE: (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 85:

40

ATTGAGGGCA AAGAGTGTGT TTGCCACACA AAGATCCTCC

```
(2) INFORMATION FOR SEQ ID NO:
          (i) SEQUENCE CHARACTERISTICS:
                    (A) LENGTH:
                    (B) TYPE:
                                               nucleic acid
                    (C) STRANDEDNESS:
                                               single
                    (D) TOPOLOGY:
                                               linear
         (ii) SEQUENCE DESCRIPTION: SEO ID NO: 86:
TTGCTTCTGA GCACAGGGCG CAGCCATGAC TGATGTTGCT G
                                                                                 41
 (2) INFORMATION FOR SEQ ID NO:
          (i) SEQUENCE CHARACTERISTICS:
                    (A) LENGTH:
                    (B) TYPE:
                                               nucleic acid
                    (C) STRANDEDNESS:
                                               single
                    (D) TOPOLOGY:
                                               linear
         (ii) SEQUENCE DESCRIPTION: SEO ID NO: 87:
ATCTGTGGGC ATGGTACGTT TGCTGAAGGA CAGTGGCAGA
                                                                                 40
(2) INFORMATION FOR SEQ ID NO:
          (i) SEQUENCE CHARACTERISTICS:
                    (A) LENGTH:
                    (B) TYPE:
(C) STRANDEDNESS:
                                              nucleic acid
                                              single
                    (D) TOPOLOGY:
                                               linear
        (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 88:
GCGGCGGCCA TGGCGGGACA GGAGGATCCG GTGCAGCGGG AGATTCACCA
                                                                                 50
GGACTGGGCT AACCGGGAGT ACATTGAGAT AATCACCAGC AGCATCAAGA
                                                                                100
AAATCGCAGA CTTTCTCAAC TCGTTCGATA TGTCTTGTCG TTCAAGACTT
                                                                                150
GCAACACTAA ACGAGAAATT GACAGCCCTT GAACGGAGAA TAGAGTACAT
                                                                                 200
TGAAGCTCGG GTGAC
                                                                                 215
(2) INFORMATION FOR SEQ ID NO:
                                      89.
         (i) SEQUENCE CHARACTERISTICS:
                   (A) LENGTH:
                                              1434
                   (B) TYPE:
                                              nucleic acid
                   (C) STRANDEDNESS:
                                              single
                   (D) TOPOLOGY:
                                               linear
        (ii) SEQUENCE DESCRIPTION: SEO ID NO: 89:
CCAGCAATCT ATCATGGATC CTAATCAGAA CGTGAAATGC AAGATAGTTG
                                                                                50
TGGTGGGAGA CAGTCAGTGT GGAAAAACTG CGCTGCTCCA TGTCTTCGCC
                                                                                100
AAGGACTGCT TCCCCGAGAA TTACGTTCCT ACAGTGTTTG AGAATTACAC GGCCAGTTTN GAAATCGACA CACAAAGAAT AGAGTTGAGC CTGTGGGACA
                                                                                150
                                                                                200
CTTCGGGTTC TCCTTACTAT GACAATGTCC GCCCCCTCTC TTACCCTGAT
                                                                                250
TCGGATGCTG TGCTGATTTG CTTTGACATC AGTAGACCAG AGACCCTGGA
                                                                                300
CAGTGTCCTC AAAAAGTGGA AAGGTGAAAT CCAGGAATTT TGTCCCAAAT ACCAAAATGC TCTTGGTCGG CTGCAAGTCT GATCTGCGGA CAGATGTTAG
                                                                                350
                                                                                400
TACATTAGTA GAGCTCTCCA ATCACAGGCA GACGCCAGTG TCCTATGACC
                                                                                450
AGGGGGCAAA TATGGCCAAA CAGATTGGAG CAGCTACTTA TATCGAATGC
TCAGCTTTAC AGTCGGAAAA TAGCGTCAGA GACATTTTTC ACGTTGCCAC
                                                                                500
```

550

CTTGGCATGT	GTAAATAAGA	CAAATAAAAA	CGTTAAGCGG	AACAAATCAC		600
AGAGAGCCAC	AAAGCGGATT	TCACACATGC	CTAGCAGACC	AGAACTCTCG		650
GCRGTTGCTA	CGGACTTACG	AAAGGACAAA	GCGAAGAGCT	GCACTGTGAT		700
GTGAATCTTT	CATTATCTTT	AATGAAGACA	AAGGAATCTA	GTGTAAAAAA		750
CAACAGCAAA	CAAAAAGGTG	AAGTCTAAAT	GAAGTGCACA	GCCAAAGTCA		800
TGTATACCAG	AGGCTTAGGA	GGCGTTTGAG	AGGRTACTCA	TCTTTTTGGG		850
AATCCTGACC	TTAGGTTCGG	CATGTAGACC	AAGTGATGAG	AAGTGAATAC		900
ATGGAAGAGT	TTTTAAGTGT	GACTTGAAAA	ATATGCCAAA	AAATGAGAGA		950
TACAAATGAG	CTAGAGGAAG	ATGAGGGGG	ATGCGAGTAC	CTCCAAGAAG	1	000
AAAAATCACA	CTCTGAATGG	TGCTTGCATT	TTGAGGTTTN	NNNNCNNNNG	1	050
GGGTATAATC	TATCATGGAT	CTCCACTTGG	ATTAATTTA	AATGTTTAAT	1	100
CTCCTTTACA	AAAAGTATAC	GTTAATATAC	CGTCCTCAAG	GGGGACTGGC	1	150
ACTGTGACCT	AGCATTAGTT	TCTAGAGGAT	GTGATCTAAT	TCTTCTAGCT	1	200
CATCATAAAA	AGGAATTGTA	TCAGGACCCA	TGGGATATAT	CCAGAGGCAA	1	250
CTTATGAGGC	TTGGAATCTG	GCTTCCTGAA	GATAGCTGAG	TAGGATGGTC	1	300
TAAGGAAAGC	CTTGGAATCT	TGCAAGATTG	GTGGACCAGC	ACTACAAAGA	1	350
TCGCATAGAT	CAAATAGGAA	AAAAATGTCG	ATTTTTATTC	AGTCTGATGG	1	400
TTCTGTTCTC	ATGGTGGATG	GTCATAAAAA	GTGG		1	434
(2) INFORMA	TION FOR SE	EQ ID NO:	90:			
(5)	SPOTTENCE (יים ד משרייי מקמעי	rt/cg •			

(i) SEQUENCE CHARACTERISTICS:

(A)	LENGTH:	265	
(B)	TYPE:	nucleic	acid
(C)	STRANDEDNESS:	single	
(D)	TOPOLOGY:	linear	

(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 90:

GGCACGAGCG	GAATTGGACT	TGGGAGGCGC	GGTGAGGAGT	CAGGCTTAAA	5	50
ACTTGTTGGA	GGGGAGTAAC	CAGCCTGCTC	CTCTCGCTCT	CCTCCTCGTC	נ	100
TGCGCCGCGT	TTCAGAGGTT	GCCCATCAGC	CTTGTGATTT	ATTTTTATAT	3	150
CTGCTTTTTA	TAAAGAGAGA	AATATATATA	TATATATATA	TATTTTTTTT	2	200
TTCTTCTTAA	GAGAAAATTC	CTGTTCCAAG	AGAAAATAAG	GCAACATCAA	2	250
TGAAGGAGAG	AAGAG				2	265

(2) INFORMATION FOR SEQ ID NO: 91:

(i) SEQUENCE CHARACTERISTICS:

(A)	LENGTH:	223	
(B)	TYPE:	nucleic aci	ď
(C)	STRANDEDNESS:	single	
(D)	TOPOLOGY:	linear	

(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 91:

GATTGATGCA	GCATTATGCT	TTGGGCAGTA	TTACAAAATA	GCTGGCGAGT	!	50
KCTTTCTGTA	TTTAAATATT	GTAAAAAGAA	AATAAGTTAT	AACTGTTATA	:	100
AAGCAGAACT	TTTGTTGCAT	TTTTTAAACT	GTTGAAGTCA	CTGTGTATGT	:	150
TTGTTTGGTC	AATGTTTCCG	CAGTATTTAT	TAAAACATAC	TTTTTTTTT		200
CTTCAAATAA	AAAAGTAACC	ATG				223

(2) INFORMATION FOR SEQ ID NO: 92:

(i) SEQUENCE CHARACTERISTICS:

(A)	LENGTH:	65
(B)	TYPE:	nucleic acid
(C)	STRANDEDNESS:	single
(D)	TOPOLOGY:	linear

(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 92:

WO 96/13610		PCT/US95/11230
TCCTCTCGCT CTCCTCCTCG TCTGCGCCGC TTTCAGAG. TTCCAAGAGA AAATA	AG AAAATTCCTG	50 65
(2) INFORMATION FOR SEQ ID NO: 93:		
(i) SEQUENCE CHARACTERISTICS:		
(B) TYPE: (C) STRANDEDNESS:	150 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 93:	
ACCATGGAAA GTATAGTGAA TGAATATGCC TTCCTACTC AAAAAAGCCC ATGACAAATG AGAAACAAAA TTCCATTT TTCTGAGTTG TCTAAAGCCC AACTCCAAGT TAATTCAAC	rg gccaacatta	50 100 150
(2) INFORMATION FOR SEQ ID NO: 94:		
(i) SEQUENCE CHARACTERISTICS:		
	108 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 94:	
GTTTTTTTTG AGTTTAACAC AGATTTTATT GCCCTATAG GTGACCAGTG GATATCAATG AAACTTCTTA ATTATTTGA CATATTTA		50 100 108
(2) INFORMATION FOR SEQ ID NO: 95:		
(i) SEQUENCE CHARACTERISTICS:		
(C) STRANDEDNESS:	107 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 95:	
ACAAACCACA GTATTTCCAC TTTAAATATA GAACTGGTA CCTTAAACTA AAATCGGTGA CTCGGTCATC AAGAAGGTT GTGGATC		50 100 107
(2) INFORMATION FOR SEQ ID NO: 96:		
(i) SEQUENCE CHARACTERISTICS:		
(B) TYPE: (C) STRANDEDNESS:	83 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 96:	
TTTTTTTTT TTTTTTTTT TTTTTTTTT GAAAGAATA	G GTTTAATTTA	50 83
(2) INFORMATION FOR SEQ ID NO: 97:		
(i) SEQUENCE CHARACTERISTICS:		

(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	87 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 97:	
GCTAAACCAA ACCAACTCC TAGTTTAGAA CTCTCTGCA	T CTGCTTTGTC CCTTAGGTC T AGGGGTGGGA ATTAATC	eg aaaagagagg	50 87
(2) INFORMATION FOR	SEQ ID NO: 98:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	STRANDEDNESS:	67 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 98:	
TGTTATTATT AGCTTTCCA AAAAATACTC TTCTCTT	T GTGGCAGGTA TTAAAATC	TT TATATTTCTG	50 67
(2) INFORMATION FOR	SEQ ID NO: 99:		
(i) SEQUENCE	CHARACTERISTICS:		
(B)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	103 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 99:	
GATGGGGACG TCCTGATTT GCACACTAGT TTTCTCAGG TTG	'A CCAGATCAAA GTATGGTA 'A TTCAATAACC GCAAACTC	AG GCTGTAGATA TT GGTGCCCCTA	50 100 103
(2) INFORMATION FOR	SEQ ID NO: 100:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	173 nucleic acid single linear	·
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 100:	
ATGACAGCCC AATTTTTT	T TTAGGAAAGA AAGAATTC AA AATGGTTATC TTAAGTCA AT AACATCAGAT ATTTCTAA CA TGT	GG CCAGTTTTAT	50 100 150 173
(2) INFORMATION FOR	SEQ ID NO: 101:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)		253 nucleic acid single linear	

114

(ii) SEQUENCE DESCRIPTION: SEQ II	D NO: 101:	
CCCACCGGGG CCATGGCGTC CGCCGACATC GTCGTGGCCGGCGCCC TACCTCCAGG ATTATTTTAC AAATGCAAAAAAGATGC TCAGCAAGAT TACCATCTAG AATATGCAACACACAAA TAATTGAATT TACCAGAGAG CTGCATACTGACAAGAGT ATAACGGATA GCACTGTGAG AGTGATCTATGCAATGAATG	GCG GGGTGGCCCA 50 AAT AGAGAGTTGA 10 CAT GGAAAATAGC 15 CAT GTGACATAAA 20	50
(2) INFORMATION FOR SEQ ID NO: 102:		
(i) SEQUENCE CHARACTERISTICS:		
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	261 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ II	D NO: 102:	
GTTTATTAAT ATGATCCAAA CAGATTTAGC TTAAGAGA ATAAAACAAG ATTTGGCTGA GCTATATAAA GAAAACCA GATTTTCAAT TGGTGACAAT TTAAAGTCAT TTATTAAAC AAATCTACAA TAGAACCATG ACTTAAACTG AGCCAGTA GTCACCCAAA GTCAGCAAGG AACCATTCAG AACTTACA GGCTCACACG G	AAA AGAGTGACAA 10 CTA ATATAAAAAT 15 ATG TGGTACATGT 20	00000
(2) INFORMATION FOR SEQ ID NO: 103:		
(i) SEQUENCE CHARACTERISTICS:		
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	120 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ II	D NO: 103:	
TTTTTTTTT TTTTAAACTT AAAAATCAGT TTATCTTAATTAGAAGAA ATTGATTTCA CATGAAAATA TAACATTC	- -	0
(2) INFORMATION FOR SEQ ID NO: 104:		
(i) SEQUENCE CHARACTERISTICS:		
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	129 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ II	O NO: 104:	
TGAAAATGCC AAAAATAGTG AACTTAACCT TGTTGGCT ACCAGCCAAC ATTATTGTGC CTTTGCATTT CACACCAA GGTATTTTT GATGATGCAA ATGTGAAAG		0
(2) INFORMATION FOR SEQ ID NO: 105:		
(i) SEQUENCE CHARACTERISTICS:		
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: //3	214 nucleic acid single	

(D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 105: 50 GTAGCATATG ACAAATGAAT TATTATATCT CCCTGATTTT CTATTTTCA CGGCACGCCA AGACTATTCA AAGGGAATTT CTAAAATATC TAAGATATTG 100 150 ATCATAAAAT TATATATGTG TATACACATA TACACATACA TATGTATAGT GATTAACAAC TTTA 200 214 (2) INFORMATION FOR SEQ ID NO: 106: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 106: CATGTTTGTC AGGGAGATGC CTGTATATTT TGTGCATTGA AGACGATATT 50 TGCACAGTTC CAACACAGTC GAGAAAAAGC ACTTCCCTCA GATAACATAA 100 GGCATGCTCT TGCAGAAAGT TTCAAAGATG AGCAGCGATT TCAACTTGGC 150 CTTATGGATG ATGCTGCGGA GTGCTTTGAA AATATGTTGG AGAGGATTCA 200 214 TTTTCACATA GTGC (2) INFORMATION FOR SEQ ID NO: 107: (i) SEQUENCE CHARACTERISTICS: 264 (A) LENGTH: (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 107: TAAAAATACT ATCAACATTT AATGTTTCTC TCATAGGTTT CCAACCTCGG TTCCGGCTTT TACTGCTGCT ATCACAGCTG TTTCCTCTAT CCCTAGAATC 50 100 TTGGCTGCTG TTCTGTGTCA TATCCAGTGC CATTGTCACT CTTAACTTTG 150 CCAGTTATTT TCTCTCCTGG AGCAAGAATC TGGGATTTAC TTGAACTTAT TATTTGTGCA GAAGCTMGAC TYTGAGGTAC AACTCGGTCR TGTTTATATG 200 250 ATCCTTTTCC TTGA 264 (2) INFORMATION FOR SEQ ID NO: 108: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: nucleic acid single (C) STRANDEDNESS: (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 108:

(2) INFORMATION FOR SEQ ID NO: 109:

(i) SEQUENCE CHARACTERISTICS:

CCCACGCCTC GGCGGCCGCC AGTGCTCCTG GAGTGCAGA

(A) LENGTH: 166

GCGCGCCGGC AGCTGCAGGC TCAGGCCATC TCCAGTCTTC CTGCCGCGCC CGCCTCTCCT CCTGGGATTC TCCTCCTCCT CCTGGACTTC CCCGCAGCTG

50 100

139

(c)	TYPE: STRANDEDNESS: TOPOLOGY:	nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 109:	
AAGGTCGGAG TCAACGGAT	G CCACATCGCT CAGACACC. I TGGTCGTATT GGGCGCCT I GGATATTGTT GCCATCAA	GG TCACCAGGCT	50 100 150 166
(2) INFORMATION FOR S	SEQ ID NO: 110:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	196 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 110:	
CTCCCTAGGC CCCTCCCCTCGGGGGAGATTC AGTGTGGTGG	TACTTTATTG ATGGTACAT TTCAAGGGGT CTACATGG GGGACTGAGT GTGGCAGGG TCTTGTGCTCT TGCTGGGGG	CA ACTGTGAGGA GA CTCCCCAGCA	50 100 150 196
(2) INFORMATION FOR S	SEQ ID NO: 111:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	148 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 111:	
GATATCCTCT GAAAAGCATC	CTAAGGCAGA TAGACTGTC AAAATCTTCC ACTGTGAAC CTTCCAATTC TGAGCGAAC	CA CATGGGTCTC	50 100 148
(2) INFORMATION FOR S	EQ ID NO: 112:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	LENGTH: TYPE: STRANDEDNESS: TOPOLOGY:	88 nucleic acid single linear	
(ii) SEQUENCE	DESCRIPTION: SEQ ID	NO: 112:	
GCAGGGGCTG AAATAACCCA AGCCCATAGT TGGCCCCCCT		TA TACTAGATGC	50 88
(2) INFORMATION FOR S	EQ ID NO: 113:		
(i) SEQUENCE	CHARACTERISTICS:		
(B) (C)	TYPE: STRANDEDNESS:	189 nucleic acid single linear	

(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 113:	
CTGAAATGCA GCTCCCTGTC CAAGTGCCTT GGAGAACTCA CAGCAGCACG CCTTAATCAA AGGTTTTACC AGCCCTTGGA CACTATGGGA GGAGGGCCAA GAGTACACCA ATTTGTTAAA AGCAAGAAAC CACAGTGTCT CTTCACTAGT CATTTAGAAC ATGTTATCAT CCAAGACTAC TCTACCCTG	50 100 150 189
(2) INFORMATION FOR SEQ ID NO: 114:	
(i) SEQUENCE CHARACTERISTICS:	·
(A) LENGTH: 151 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 114:	
CCTCCGCTTA CAGCTCGCTG CCGCCGTCCT GCCCCGCGCC CCCAGGAGAC CTGGACCAGA CCACGATGTG GAAACGCTGG CTCGCGCTCG CGCTCGCGCT GGTGGCGGTC GCCTGGTCCC GCCGAGGAAG ACTAAGACAA ATCAAGATCT T	50 100 150 151
(2) INFORMATION FOR SEQ ID NO: 115:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 146 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 115:	
TTTTTTTTT TTWWTTTYTT TATTCATCAA TAGTATCCGA AAAGGAAGAA TCAGGAGTTA CAAAAACAAG TTAAATGCAA TATAGAAGCC TACTAAATAC AAATACAAGT CACAAACACA TATGCAAGAG AAACTTGTTT AGATTG	50 100 146
(2) INFORMATION FOR SEQ ID NO: 116:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 145 (B) TYPE: nucleic acid	
(C) STRANDEDNESS: single (D) TOPOLOGY: linear	·
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 116:	
CCGGGCCGGG GAGGCGCGCT CGCTCCGCGC TCCCTTCGCT CGCTCGTTTC CTCCTCCTC GGCAGCCGCG GCGCAGCAG GAGAAGGCGG CGGCGGCGGC TAGGGATCAG ACATGGCGGC GGATCTGAAC CTGGAGTGAT CTCCT	50 100 145
(2) INFORMATION FOR SEQ ID NO: 117:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 216 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 117:	
TTTTTTTTT TTTTTTTAA GAACATCAAC ATTTATTTAA CATGATAAAA	50

WO 96/13610	PCT/US95/11230
AAAGAAATGA GATATGAACA TTTGCATTTA AACAATAGTA AGTAGCCTTT AATACATTAC ATGTGCTCAT TGTATAATAT ATACACAATG AACATAATTA CATTTGTACA CAAACTAAGT ACCGGATTTG GAAACCTGCT TATTGCTGTA CACATGTATT CCAATG	150
(2) INFORMATION FOR SEQ ID NO: 118:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 191 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 118:	
CACCACGATC AAAAGGGACA AGCATCAAGC ACGCAGCAAT GCAGCTCAAA ACGCTTAGCC TAGCCACACC CCCACGGGAA ACAGCAGTGA TTAACCTTTA GCAATAAACG AAAGTTTAAC TAAGCTATAC TAACCCCAGG TTGGTCAATT TCGTGCCAGC CACCGCGTCA CACGATTACC AAGTCATAGA G	100
(2) INFORMATION FOR SEQ ID NO: 119:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 226 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 119:	
TTTTTTTTTT TTTTTTTTK YTYWTTTYYT TGTTTTTTG ATTGTTTTGG ATCTCTGGTT TAATTAGCAC TCTATGGTTG GGAATGTTAT TGGTTTCTTT AGTTGGTGCA TTTTCAGATG TAATCTTGTC CACTCTTYTC ACAGGTTCTG TCTGTACTAG GGCAGCATCT AACATGGCTT TCATCCACAA CTCCATTTCC TTTCCTGTAT CAGTGCAGAA AAAAGG	50 100 150 200 226
(2) INFORMATION FOR SEQ ID NO: 120:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 196 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 120:	
CGGATCTGAA CCTGGAGTGG ATCTCCCTGC CCCGGTCCTG GACTTACGGG ATCACCAGGG GCGGCCGAGT CTTCTTCATC AACGAGGAGG CCAAGAGCAC CACCTGGCTG CACCCCGTCA CCGGCGAGGC GGTGGTCACC GGACACCGCG GCAGAGCACA GATTTGCCTA CTGGCTGGGA AGAAGCATAT ACTTTT	100
(2) INFORMATION FOR SEQ ID NO: 121:	•
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 205 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	

(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 121:

PCT/US95/11230 WO 96/13610 AGATGGTCCA GCTGCCAGGA CTACTTTGGC AGGCAGCGTG CTACAGGACG 50 AAAATGTAAG AGAAGTCTAT TAAGGCTGGA CAGCCCAGGG TTATTTATAC 100 TCTCTCAGCC CCAAGTCCCC CGGACTAAAG ACCTAAAGGC TGATTGACTC 150 ATTCCTGATT GATTTAATGG AAAGTCTCCC ACCCCATCAT CATTTGCCAG 200 AGTAC 205 (2) INFORMATION FOR SEQ ID NO: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 199 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 122: CAAAACTCAG CAGTGCTTCT GGTGCTGGTG ATCAGTGCTT CTGCAACCCA TGAGGCGGAG CAGAATGACT CTGTGAGCCC CAGGAAATCC CGAGTGGCGG CTCAAAACTC AGCTGAAGTG GTTCGTTGCC TCAACAGTGC TCTACAGGTC 50 100 150 GGCTGCGGGG CTTTTGCATG CCTGGAAAAC TCCACCTGTG ACACAGATG 199 (2) INFORMATION FOR SEQ ID NO: 123: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 123: TTAGATGCAA GAAGATGCAG GCTCAAAGTC TGGTTGGACA GCCAGGCTCA 50 AGCAATTTGG TAAATGTGTC GGAAAGAAAA TTAGACATTG GAGGATCAAG ACCATAAGAC ACTAGCTCAT TAGAGATCAA GAATTCAAAT GTGACATTCA 100 150 TATTCGTCC 159 (2) INFORMATION FOR SEQ ID NO: 124: (i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 54 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 124: ATTTTAAAAG AATACTGTAC TGGAATGGGT TAACATCATC ATCTTTGGCA 50 54 (2) INFORMATION FOR SEQ ID NO: (i) SEOUENCE CHARACTERISTICS: (A) LENGTH: (B) TYPE: (C) STRANDEDNESS: nucleic acid single (D) TOPOLOGY: linear (ii) SEQUENCE DESCRIPTION: SEQ ID NO: 125:

50

100

GACATTCGCC CTGATATAAA AGATGATATA TATGACCCCA CCTACAAGGA

TAAGGAAGGC CCAAGCCCCA AGGTTGAATA TGTCTGGAGA AACATCATCC

WO 96/13610		PCT/US95/11230
TTATGTCTCT GCTACACTTG GGAGCCCTGT ATGGGATCACCTGCAAGT TCTACACCTG GCTTTGGGGG GTATTCTA	CAC TTTGATTCCT ACT ATTTTG	150 196
(2) INFORMATION FOR SEQ ID NO: 126:		
(i) SEQUENCE CHARACTERISTICS:		
	184 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ II	D NO: 126:	
TCGCACTTTG CCCCTGCTTG GCAGCGGATA AAAGGGGGCCGGGACACGGT CACCCGTTGC CAGCTCTAGC CTTTAAATGACCTCCACG CACCGCGGCT AGCGCCGACA ACCAGCTACCGCGGCTCA GCGCGGCTTCG AAAC	TTC CCGGCTCGGG	50 100 150 184
(2) INFORMATION FOR SEQ ID NO: 127:		
(i) SEQUENCE CHARACTERISTICS:		
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	217 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 127:	
TTTTTTAAAT TACAACACTT TATTGCAGCA TCGGCAAA GAAGCTGGTG AAGATTGGGC AGCATTTCCA TGTGAAAT CAAGTTTTGT TTTTTATTTA AATCTACATG CAGAAACT AAGAAAAAAT GCAAAATAGC TAGAAAAAAA GATGTAAT TACAGATGTG CTCTCCG	TGT TACAACTTTA TGA AACATGGTAA	50 100 150 200 217
(2) INFORMATION FOR SEQ ID NO: 128:		
(i) SEQUENCE CHARACTERISTICS:		
(A) LENGTH: (B) TYPE: (C) STRANDEDNESS: (D) TOPOLOGY:	174 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 128:	
CTTGTCTGAC CTAGAAGCTG AGAAACTGGT GATGTTCC ACAAACCTGG CTTGCTGATG ATGTGCTTCA TCCTGCCC TGGTATTTCT GGGGTGAAAC TTTTCAAAAC AGTGTGTT CTTGCGATAT GCTGTGGTGC TTAA	AC GCTTGTGCCC	50 100 150 174
(2) INFORMATION FOR SEQ ID NO: 129:		
(i) SEQUENCE CHARACTERISTICS:		
	234 nucleic acid single linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID	NO: 129:	
TTTTTTCAGA TTCACTTCAC TTTTATTATG AACAAACA AGTACAATTA GCTTCAGAGT TGATATTAAT AGAAATTA		50 100

WO 96/13610	PCT/US95/11230
TCTTGTCACA AGTAACTACT ATATCCCACA TAAAAAGGGA AAAAATCCCA CCCAATCACA GAAAAGGCAT CCTCTGTATG TTTCCGTGGC AATGCGTTGT TTATGTATTC TCAAATTTTG TCTGGCTAGT TATC	150 200 234
(2) INFORMATION FOR SEQ ID NO: 130:	
(i) SEQUENCE CHARACTERISTICS:	
(A) LENGTH: 261 (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
(ii) SEQUENCE DESCRIPTION: SEQ ID NO: 130:	
TTGGTGAGGG TGGGGAGGAG TGCTTCGAAG GGAGAAGCCC CAGCAAGATT TATTCCTTTT TGCTTCTTCT TCTCCCTGTC CCTGCCATAA CCATGAAGCC TTGAACAAAC CACCCAAATC TCAGGATCTT AGTGTTTTCT CTGTAAATTG TAATATGAAC TTATAAAGAT CCTCCATTGC TGATAGTCTC AGGTTCTGTG AGTAACAGCA AAAAAACTTT GTATCTAACT TCAACCAGAG CAGGCTGTAC CCTTAAGCTC T	50 100 150 200 250 261

We Claim:

1. A method for identifying and isolating a senescence-related gene, which method comprises the steps of:

- (a) isolating mRNA from a senescent cell and a young quiescent cell;
- (b) amplifying in separate reaction mixtures said mRNA from said senescent cell and said young quiescent cell to produce amplified gene sequences;
- (c) separating said amplified gene sequences in each reaction mixture by size and/or charge; and
- (d) analyzing said amplified gene sequences separated in step (c) to identify an amplified gene sequence from young quiescent and young dividing cells that is present at a level different from that observed from said senescent cells.
- 2. The method of claim 1 wherein said amplifying is by a polymerase chain reaction.
- 3. The method of claim 1 wherein said separating is by gel electrophoresis.
- 4. The method of Claim 1, wherein said amplifying step comprises the steps of:
- (i) aliquoting each different mRNA preparation into at least four and up to twelve different aliquots and adding to each aliquot reverse transcriptase and a primer selected from a group of primers consisting of primers having a 3'-end defined by a sequence $5'-T_{6-12}R_1R_2-3'$, R_1 is A, G, or C; and R_2 is A, G, C, or T, under conditions such that reverse transcriptase-mediated primer extension can occur;
- (ii) further subdividing the aliquots prepared in step (1) into at least four and up to twenty different reaction vessels and adding to each vessel DNA polymerase

and a primer selected from a group of primers consisting of primers each having a different 3'-end from other primers in said group, said 3'-end defined by a random sequence of 10 to 14 nucleotides, under conditions such that DNA polymerase-mediated primer extension can occur.

- 5. The method of Claim 1, wherein said analysis comprises physically removing an amplified gene sequence from the gel and determining its nucleotide sequence, either directly or after cloning into a recombinant DNA vector.
- 6. The method of Claim 4, wherein said primers are selected from the group consisting of primers in Table 1, above.
- The method of Claim 6, 7. wherein said DNA polymerase is a thermostable DNA polymerase; step (ii) further comprises extending said primers in a polymerase chain reaction of at least ten cycles of primer extending, annealing, and denaturing; said cycles comprise a first set of two to four cycles of lowtemperature primer annealing and a second set of at least six to eight cycles of high-temperature primer annealing; and said low-temperature primer annealing is conducted at a temperature at least 10°C below said high-temperature annealing.
- 8. A purified oligonucleotide probe comprising a sequence of at least 12 contiguous nucleotides identical or complementary to a contiguous sequence from a nucleic acid in Table 4 or Table 6, above.
- 9. A method for detecting senescent cells and for distinguishing senescent cells from non-senescent cells, said method comprising the steps of:

(a) contacting mRNA present in a cell or tissue with a labelled nucleic acid probe that comprises a sequence of a senescence-related gene under conditions such that complementary nucleic acids hybridize to one another;

- (b) determining whether specific hybridization has occurred; and
- (c) correlating the presence of senescent and non-senescent cells with the occurrence of hybridization.
- 10. The method of Claim 9, wherein said labelled nucleic acid probe comprises a sequence of at least 12 contiguous nucleotides identical or complementary to a contiguous sequence from a nucleic acid in Table 4 or Table 6, above.
- 11. The method of Claim 10, wherein said mRNA is in a tissue section, and said contacting step is conducted in situ.
- 12. A method for screening compounds to identify compounds that can alter gene expression in senescent cells, which method comprises:
 - (a) contacting senescent cells with a compound;
- (b) determining mRNA expression patterns in said senescent cells by determining expression levels of mRNA of two or more senescence-related genes; and
- (c) correlating an alteration in mRNA expression of a senescence-related gene with a compound that can alter gene expression in senescent cells.
- 13. The method of Claim 12, wherein said senescencerelated genes are selected from the group consisting of genes in Table 2, above.

14. The method of Claim 12, wherein said compounds have, prior to step (a), been tested for ability to modulate activity or expression levels of a first senescence-related gene product and determined to have said ability.

- 15. The method of Claim 12, wherein step (b) comprises determining expression levels of mRNA of more than two but less than five senescence-related genes.
- 16. The method of Claim 12, wherein step (b) comprises determining expression levels of mRNA of more than five but less than ten senescence-related genes.
- 17. The method of Claim 12, wherein step (b) comprises determining expression levels of mRNA of more than ten but less than twenty senescence-related genes.
- 18. The method of Claim 14, wherein said first senescence-related gene product is beta-galactosidase.
- 19. The method of Claim 14, wherein said first senescence-related gene product is collagenase.
- 20. The method of Claim 14, wherein said first senescence-related gene product is IFN gamma.
- 21. Method for extending the proliferative capability of a population of cells from a donor comprising the steps of:
 - a) isolating cells from said donor;
- b) separating senescent cells from young cells;and,
 - c) reintroducing said young cells into said donor.

22. The method of claim 21 wherein said method further comprises, after step b), expanding said young cells in the presence of an agent which slows the replicative senescence of said cells.

- 23. A method for destroying a cell expressing a senescence-related gene product, wherein said method comprises the step of providing a toxic substance which specifically destroys said cell..
- 24. The method of claim 23 wherein said toxic substance is a toxic substrate, wherein said toxic substrate is activated by a senescence-related gene product.
- 25. The method of claim 24 wherein said toxic substate is activated by beta-galactosidase.
- 25. A method for destroying a cell expressing a senescence-related gene product wherein said method comprises the step of providing a antibody to said senecence-related gene product, wherein said antibody is linked with a toxic substance.

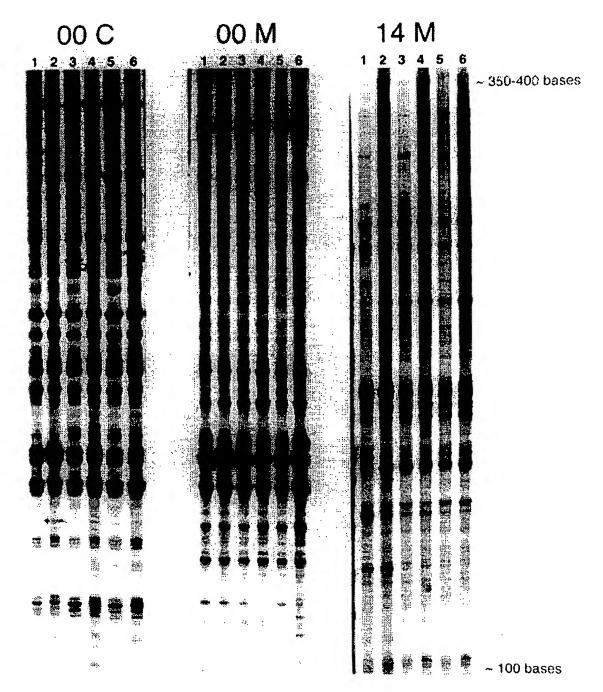
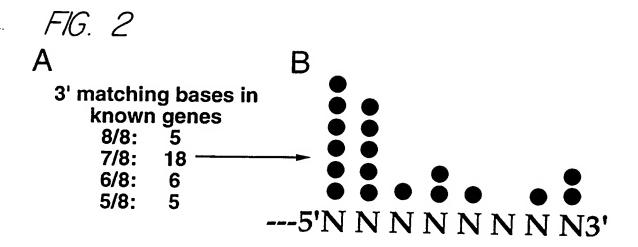
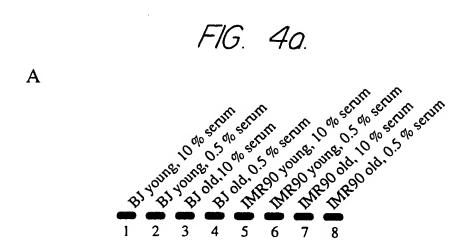


FIG. 1.

SUBSTITUTE SHEET (RULE 26)







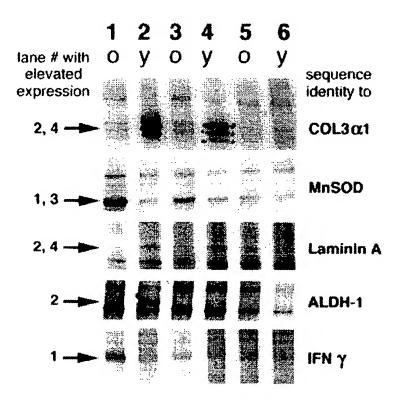


FIG. 3.

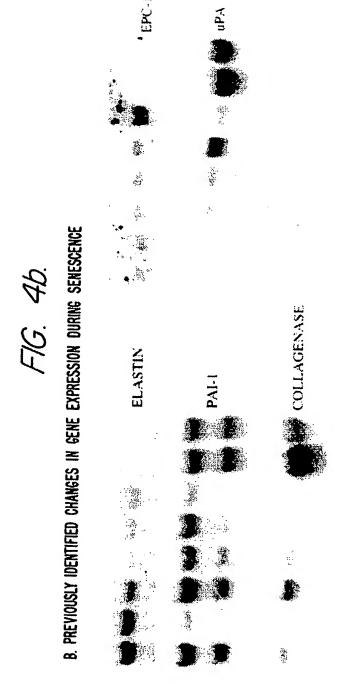

SUBSTITUTE SHEET (RULE 26)

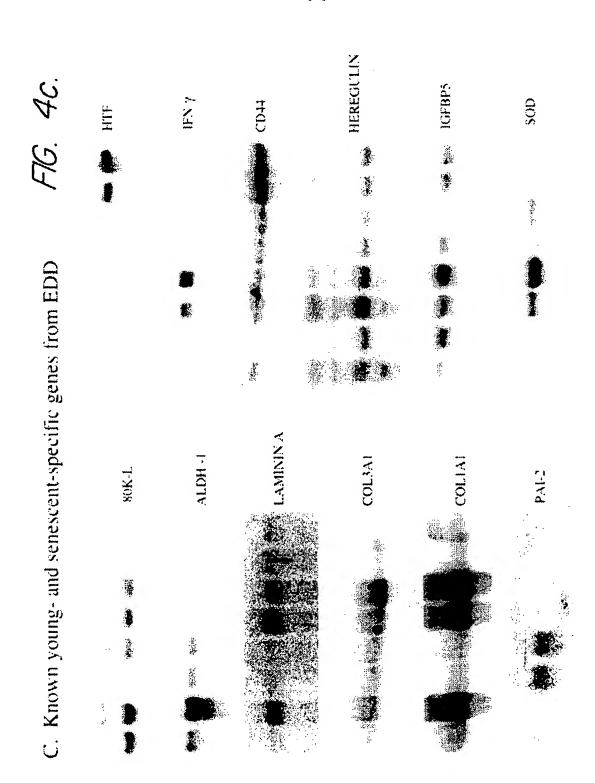

SUBSTITUTE SHEET (RULE 26)

